




Instruction For Use

Rotors and Tubes

For Beckman Coulter
Preparative Ultracentrifuges

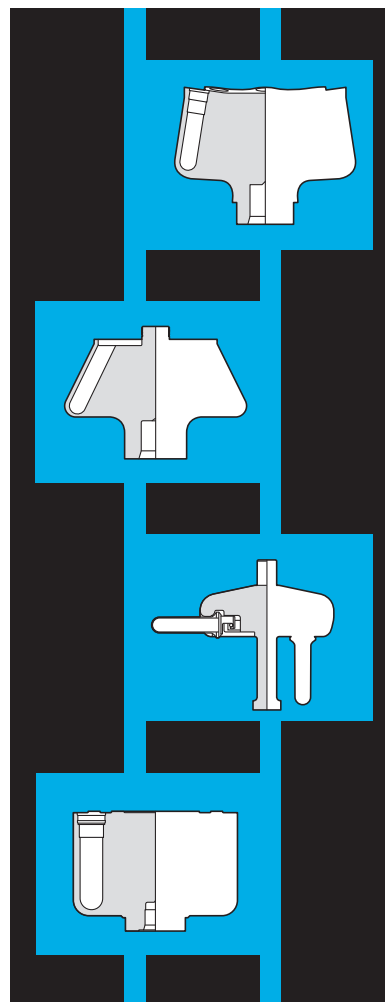

LR-IM-24AK
August 2022



Beckman Coulter, Inc.
250 S. Kraemer Blvd.
Brea, CA 92821 U.S.A.



**Rotors and Tubes for Beckman Coulter Preparative
Ultracentrifuges Instructions for Use**
LR-IM-24AK (August 2022)

© 2022 Beckman Coulter, Inc.
All rights reserved.

Trademarks

Beckman Coulter, the stylized logo, and the Beckman Coulter product and service marks mentioned herein are trademarks or registered trademarks of Beckman Coulter, Inc. in the United States and other countries.

All other trademarks, service marks, products, or services are trademarks or registered trademarks of their respective holders.

Contact Us

If you have any questions, contact our Customer Support Center.

- Worldwide, find us via our website at www.beckman.com/support/technical
- In the USA and Canada, call us at 1-800-369-0333.
- In Austria, call us at 0810 300484
- In Germany, call us at 02151 333999
- In Sweden, call us at +46 (0)8 564 859 14
- In Netherlands, call us at +31 348 799 815
- In France, call us at 0825838306 6
- In the UK, call us at +44 845 600 1345
- In Ireland, call us at +353 (01) 4073082
- In Italy, call us at +39 0295392 456
- In other locales, contact your local Beckman Coulter Representative.

Find us on the World Wide Web at:
www.beckman.com



Beckman Coulter Eurocenter S.A.
22, rue Juste-Olivier
Case Postale 1044
CH - 1260 Nyon 1, Switzerland
Tel: +41 (0) 22 365 36 11

Glossary of Symbols is available at
beckman.com/techdocs (PN C24689).

*May be covered by one or more pat. - see
www.beckman.com/patents*

Original Instructions

Revision History

This document applies to the latest version listed and higher versions. When a subsequent version changes the information in this document, a new issue will be released to the Beckman Coulter website. For updates, go to www.beckman.com/techdocs and download the latest version of the manual or system help for your instrument.

Revision AD, July 2017

- Scope
- CHAPTER 2, Introduction
- APPENDIX C, The Use of Cesium Chloride Curves

Revision AE, October 2018

- Table 2.1, Characteristics and Chemical Resistances of Tube and Bottle Materials

Revision AF, November 2018

- CHAPTER 2, Labware Selection Criteria
- CHAPTER 2, Certified Free Tubes
- CHAPTER 2, Sterile Tubes
- CHAPTER 7, Sterilization and Disinfection

Revision AG, October 2019

- Table 3.5, Available Bottles, Assembly and Operation
- Table 4.2, Maximum Run Speeds and Tube Volumes for Uncapped Tubes in Fixed-Angle Rotors

Revision AH, February 2022

- Removed references to discontinued instruments and obsolete rotors throughout the manual.
- Added references / information for the VTi 50.1 rotor throughout the manual.
- Table 3.2, OptiSeal Tubes and Accessories
- Table 3.3, Tube Cap Assemblies for Open-Top Tubes in Fixed-Angle Rotors
- Removed Appendix A, *Chemical Resistances for Beckman Coulter Centrifugation Products*. Referenced publication IN-175 instead.
- Removed outdated references in APPENDIX D, References.

Revision AJ, May 2022

- Table 1.1, Beckman Coulter Preparative Rotors by Use
- Table 4.1, General Specifications for Beckman Coulter Preparative Fixed-Angle Rotors

Revision AK, August 2022

- [Table 3.2, OptiSeal Tubes and Accessories](#)
- [Table 4.2, Maximum Run Speeds and Tube Volumes for Uncapped Tubes in Fixed-Angle Rotors](#)

Note: *Changes that are part of the most recent revision are indicated in text by a bar in the margin of the amended page.*

Safety Notice

Read all product manuals and consult with Beckman Coulter-trained personnel before attempting to operate instrument. Do not attempt to perform any procedure before carefully reading all instructions. Always follow product labeling and manufacturer's recommendations. If in doubt as to how to proceed in any situation, contact your Beckman Coulter Representative.

Alerts for Danger, Warning, Caution, Important, and Note



DANGER indicates an imminently hazardous situation which, if not avoided, will result in death or serious injury.



WARNING indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.



CAUTION indicates a potentially hazardous situation, which, if not avoided, may result in minor or moderate injury.

IMPORTANT IMPORTANT is used for comments that add value to the step or procedure being performed. Following the advice in the Important adds benefit to the performance of a piece of equipment or to a process.

NOTE NOTE is used to call attention to notable information that should be followed during installation, use, or servicing of this equipment.

This safety notice summarizes information basic to the safe operation of the rotors and accessories described in this manual. The international symbol displayed above is a reminder that all safety instructions should be read and understood before use or maintenance of rotors or accessories. When you see the symbol on other pages, pay special attention to the safety information presented. Also observe any safety information contained in applicable rotor and centrifuge manuals. Observance of safety precautions will help to avoid actions that could cause personal injury, as well as damage or adversely affect the performance of the centrifuge/rotor/tube system.

Chemical and Biological Safety

Normal operation may involve the use of solutions and test samples that are pathogenic, toxic, or radioactive. Such materials should not be used in these rotors, however, *unless all necessary safety precautions are taken.*

- Observe all cautionary information printed on the original solution containers prior to their use.
- Handle body fluids with care because they can transmit disease. No known test offers complete assurance that they are free of micro-organisms. Some of the most virulent—Hepatitis (B and C) and HIV (I–V) viruses, atypical mycobacteria, and certain systemic fungi—further emphasize the need for aerosol protection. Handle other infectious samples according to good laboratory procedures and methods to prevent spread of disease. Because spills may generate aerosols, observe proper safety precautions for aerosol containment. Do not run toxic, pathogenic, or radioactive materials in the rotor without taking appropriate safety precautions. Biosafe containment should be used when Risk Group II materials (as identified in the World Health Organization *Laboratory Biosafety Manual*) are handled; materials of a higher group require more than one level of protection.
- Dispose of all waste solutions according to appropriate environmental health and safety guidelines.
- If disassembly reveals evidence of leakage, you should assume that some fluid escaped the container or rotor. Apply appropriate decontamination procedures to the centrifuge, rotor, and accessories.

Mechanical Safety

- Use only the rotors, components, and accessories designed for use in the rotor and ultracentrifuge being used (refer to the applicable rotor manual). *The safety of rotor components and accessories made by other manufacturers cannot be ascertained by Beckman Coulter. Use of other manufacturers' components or accessories in Beckman Coulter rotors may void the rotor warranty and should be prohibited by your laboratory safety officer.*
- Do not use rotors in ultracentrifuges with any classification except those indicated in the rotor manual or engraved on the rotor.
- Rotors are designed for use at the speeds indicated; however, speed reductions may be required because of weight considerations of tubes, adapters, and/or the density of the solution being centrifuged. Be sure to observe the instructions in the applicable rotor manual.
- NEVER attempt to slow or stop a rotor by hand.
- The strength of containers can vary between lots, and will depend on handling and usage. We highly recommend that you pretest them in the rotor (using buffer or gradient of equivalent density to the intended sample solution) to determine optimal operating conditions. Scratches (even microscopic ones) significantly weaken glass and polycarbonate containers.

To help prevent premature failures or hazards by detecting stress corrosion, metal fatigue, wear or damage to anodized coatings, and to instruct laboratory personnel in the proper care of rotors, Beckman Coulter offers the Field Rotor Inspection Program (FRIP). This program involves a visit to your laboratory by a specially trained Beckman Coulter representative, who will inspect all of your rotors for corrosion or damage. The representative will recommend repair or replacement of at-risk rotors to prevent potential rotor failures. Contact your local Beckman Coulter office to request this service.

It is your responsibility to decontaminate the rotors and accessories before requesting service by Beckman Coulter Field Service.

Contents

Revision History, iii

Safety Notice, v

Scope, xv

Classification Program, xvii

CHAPTER 1: Rotors, 1-1

Introduction, 1-1

General Description, 1-1

Rotor Designations, 1-1

Materials, 1-2

Drive Pins, 1-3

Rotor Selection, 1-3

Pelleting (Differential Separation), 1-7

Isopycnic Separations, 1-10

Rate Zonal Separations, 1-12

General Operating Information, 1-13

Rotor Balance, 1-13

Overspeed Protection, 1-14

Allowable Run Speeds, 1-15

CHAPTER 2: Tubes, Bottles, and Accessories, 2-1

Introduction, 2-1

Labware Selection Criteria, 2-1

Labware Material Compatibility with Solvents and Sample, 2-2

Gradient Formation and Fractionation, 2-3

Certified Free Tubes, 2-3

Sterile Tubes, 2-3

Labware Types, 2-3

Polypropylene Tubes, 2-3

Open-Top Polypropylene Tubes, 2-4

OptiSeal Tubes, 2-4

Quick-Seal Polypropylene Tubes, 2-5

Polycarbonate Tubes, 2-5

Polyethylene Tubes, 2-5

Ultra-Clear Tubes, 2-6

Cellulose Propionate Tubes, 2-6

Stainless Steel Tubes, 2-6
 konical Tubes, 2-6
 Bottles, 2-7

Temperature Limits, 2-7

Spacers and Floating Spacers, 2-8

Adapters, 2-9

CHAPTER 3: Using Tubes, Bottles, and Accessories, 3-1

Introduction, 3-1

Gradient Preparation, 3-1

Cesium Chloride Gradients, 3-2

General Filling and Sealing or Capping Requirements, 3-3

Filling and Plugging OptiSeal Tubes, 3-4

Filling the Tubes, 3-5

Seating the Tube Plugs, 3-7

Filling the Tubes, 3-8

Layering Sample, 3-9

Layer sample onto the tube contents as follows. , 3-9

Sealing the Tubes, 3-10

Handy Tips, 3-15

Loading the Rotor, 3-16

Filling Open-Top Tubes, 3-17

Open-Top Polypropylene Tubes, 3-17

Swinging-Bucket Rotors, 3-17

Fixed-Angle Rotors, 3-17

Other Open-Top Tubes, 3-17

Polycarbonate, 3-17

UltraClear, 3-17

Polypropylene, 3-18

Polyethylene, 3-18

Stainless Steel, 3-18

Capping Tubes, 3-18

Tube Cap Assemblies, 3-19

Titanium Caps, 3-19

Aluminum Caps, 3-21

Inspecting and Lubricating Tube Caps, 3-22

Assembling Tube Caps, 3-23

Filling and Capping Tubes, 3-26

Sample Recovery, 3-27

Capped Tubes, 3-28

OptiSeal Tubes, 3-29

Removing Plugs from Tubes, 3-31

Quick-Seal Tubes, 3-32

Making Ultra-Clear Tubes Wettable, 3-33

CHAPTER 4: Using Fixed-Angle Rotors, 4-1

Introduction, 4-1

Description, 4-1

Tubes and Bottles, 4-3

Rotor Preparation and Loading, 4-3

Prerun Safety Checks, 4-4

Rotor Preparation and Loading, 4-5

Operation, 4-6

Installing the Rotor, 4-7

Removal and Sample Recovery, 4-7

CHAPTER 5: Using Swinging-Bucket Rotors, 5-1

Introduction, 5-1

Description, 5-1

Tubes and Bottles, 5-3

Rotor Preparation and Loading, 5-3

Prerun Safety Checks, 5-3

Rotor Preparation and Loading, 5-4

Operation, 5-7

Removal and Sample Recovery, 5-9

CHAPTER 6: Using Vertical-Tube and Near-Vertical Tube Rotors, 6-1

Introduction, 6-1

Description, 6-1

Vertical-Tube Rotors, 6-1

Near-Vertical Tube Rotors, 6-2

Tubes and Bottles, 6-4

Rotor Preparation and Loading, 6-4

Prerun Safety Checks, 6-4

Rotor Preparation and Loading, 6-5

Operation, 6-8

Removal and Sample Recovery, 6-9

CHAPTER 7: Care and Maintenance, 7-1

Introduction, 7-1

Rotor Care, 7-1

Decontamination, 7-3

Sterilization and Disinfection, 7-3

Inspection, 7-3

	Field Rotor Inspection Plan, 7-5
	Lubrication, 7-5
	Overspeed Disk Replacement, 7-6
	Tube, Bottle, and Accessory Care, 7-7
	Cleaning, 7-7
	Decontamination, 7-7
	Sterilization and Disinfection, 7-8
	Tube and Bottle Storage, 7-10
	Removing Jammed or Collapsed Tubes, 7-10
	Tube Cap Care, 7-11
	Cleaning, 7-11
	Decontamination, 7-12
	Sterilization and Disinfection, 7-12
	Lubrication, 7-13
	Inspection, 7-13
	Nylon Insert Replacement, 7-14
	Returning a Rotor or Accessory to the Factory, 7-14
	Diagnostic Hints, 7-15
APPENDIX A:	Use of the w2t Integrator, A-1
	w2t Integrator, A-1
	Reproducing Band Positions (Refer to Figure B-1), A-2
	Calculating Sedimentation Coefficients, A-4
APPENDIX B:	The Use of Cesium Chloride Curves, B-1
	Cesium Chloride Curves, B-1
	Typical Examples for Determining CsCl Run Parameters, B-4
APPENDIX C:	Gradient Materials, C-1
	Description, C-1
APPENDIX D:	References, D-1
	List of References, D-1
	Glossary
	Ultracentrifuge Rotor Warranty

Illustrations

- 1.1 Fixed-Angle, Swinging-Bucket, Vertical-Tube, and Near-Vertical Tube Rotors, 1-2
- 1.2 Particle Separation in Fixed-Angle, Swinging-Bucket, Vertical-Tube, and Near-Vertical Tube Rotors, 1-4
- 1.3 Sedimentation Coefficients (in Svedberg Units) for Some Common Biological Materials, 1-9
- 1.4 Nomogram, 1-11
- 1.5 Arranging Tubes Symmetrically in a Rotor, 1-14
- 3.1 Filling OptiSeal Tubes, 3-7
- 3.2 Layering Sample, 3-9
- 3.3 North American Markets, 3-11
- 3.4 EU/UK Markets, 3-11
- 3.5 Seal Former, 3-11
- 3.6 Seal Guide, 3-12
- 3.7 North American Market, 3-12
- 3.8 EU/UK Market, 3-12
- 3.9 Heat Sink, 3-13
- 3.10 Remove Tool, 3-13
- 3.11 North American Market, 3-14
- 3.12 EU/UK Market, 3-14
- 3.13 Heat Sink, 3-14
- 3.14 Removal Tool, 3-14
- 3.15 Test for Seal Leaks, 3-15
- 3.16 Spacers, 3-16
- 3.17 Tools Used to Assemble Tube Caps, 3-23
- 3.18 Tube Cap Installation, 3-25
- 3.19 Tube Cap Vise, 3-26
- 4.1 Fixed-Angle Rotors, 4-2
- 4.2 Removal Tools Used in Fixed-Angle Rotors, 4-9
- 5.1 Swinging-Bucket Rotors, 5-2
- 5.2 Install Gaskets Inside the Buckets, 5-5
- 5.3 Install a Spacer Over the Plugged Tube, 5-5
- 5.4 Install Spacers Over Sealed Tubes, 5-6
- 5.5 Place Filled Tubes in at Least Two Opposing Buckets, 5-7

- 5.6 Rotor Adapter Indentations Indicate Drive Pin Locations, 5-7
- 5.7 Lift Rotor With Two Hands and Lower Onto the Drive Hub, 5-8
- 5.8 Checking Hook-on Bucket Positions After the Rotor is Installed, 5-8
- 5.9 Lift Rotor Straight Up, 5-9
- 5.10 Remove Adapters With the Removal Tool, 5-9
- 6.1 Vertical-Tube Rotors, 6-3
- 6.2 Near-Vertical Tube Rotors, 6-3
- 6.3 Preparing a Vertical-Tube or Near-Vertical Tube Rotor, 6-7
- A.1 The sw2t Charts for the SW 60 Ti Rotor, A-3
- B.1 Precipitation Curves for the Type 90 Ti Rotor, B-2
- B.2 CsCl Gradients at Equilibrium, B-3

Tables

1.1	Beckman Coulter Preparative Rotors by Use, 1-6
2.1	Characteristics and Chemical Resistances of Tube and Bottle Materials, 2-2
2.2	Quick-Seal Tube Spacers, 2-9
2.3	Dimensions of Delrin Adapters, 2-12
3.1	Filling and Capping Requirements for Tubes and Bottles, 3-3
3.2	OptiSeal Tubes and Accessories, 3-5
3.3	Tube Cap Assemblies for Open-Top Tubes in Fixed-Angle Rotors, 3-20
3.4	Required Tools and Torque Values, 3-24
3.5	Available Bottles, Assembly and Operation, 3-27
4.1	General Specifications for Beckman Coulter Preparative Fixed-Angle Rotors, 4-3
4.2	Maximum Run Speeds and Tube Volumes for Uncapped Tubes in Fixed-Angle Rotors, 4-4
5.1	General Specifications for Beckman Coulter Preparative Swinging-Bucket Rotors, 5-3
6.1	General Specifications for Beckman Coulter Preparative Vertical-Tube and Near-Vertical Tube Rotors, 6-2
6.2	Rotor Plugs and Tools Used for Vertical-Tube and Near-Vertical Tube Rotors, 6-8
7.1	Tube and Bottle Sterilization and Disinfection, 7-9
7.2	Troubleshooting Chart, 7-15
C.1	Commonly Used Gradient Materials with Their Solvents, C-2
C.2	Density, Refractive Index, and Concentration Data—Cesium Chloride at 25°C, Molecular Weight = 168.37, C-3
C.3	Density, Refractive Index, and Concentration Data—Sucrose at 20°C, Molecular Weight = 342.3, C-4
C.4	Density Conversion for Cesium and Rubidium Salts at 20°C, C-5

Scope of this Manual

This manual contains general information for properly preparing a rotor for centrifugation in a Beckman Coulter preparative ultracentrifuge. This manual should be used with the individual rotor instruction manual. The rotor manuals provide specific information for each rotor, including special operating procedures and precautions; tube, bottle, and adapter part numbers; and equations to calculate maximum allowable rotor speeds. Each manual has a code number in the upper right-hand corner of the cover page that can be used for reordering. To reorder, contact customer service at 1-800-742-2345 in the United States; outside the U.S., contact your local Beckman Coulter representative.

A lot of information is compiled in this manual, and we urge you to read it carefully—especially if this is your first experience with Beckman Coulter products.

- [CHAPTER 1](#) describes, by usage, Beckman Coulter's currently produced preparative ultracentrifuge rotors; this should help you determine the appropriate rotor to use for a particular application. Also included in this section is a discussion of rotor materials, components, and centrifugation techniques.
- [CHAPTER 2](#) describes various tubes, bottles, adapters, and spacers to help you choose a particular tube or bottle for your application.
- [CHAPTER 3](#) provides instructions for using tubes and related accessories.
- [CHAPTER 4](#) contains step-by-step procedures for preparing a fixed angle rotor for a centrifuge run. Similar information for swinging bucket rotors is in [CHAPTER 5](#), and [CHAPTER 6](#) contains the same type of information for vertical tube and near-vertical tube rotors. (Analytical, continuous flow, and zonal rotors are not covered in this manual.)
- [CHAPTER 7](#) provides rotor, tube, bottle, and accessory care and maintenance information, as well as some diagnostic hints. Please read it. Proper rotor care results in longer rotor life.
- Several appendixes contain information that may be of special interest:
 - [APPENDIX A](#) contains information about the use of the w2t integrator.
 - [APPENDIX B](#) describes the use of cesium chloride curves.
 - [APPENDIX C](#) contains reference information on some commonly used gradient materials.
 - [APPENDIX D](#) lists references for further reading.
 - [Glossary](#) provides a glossary of terms.

Scope

Scope of this Manual

Classification Program

Introduction

All Beckman Coulter preparative ultracentrifuges are classified according to the size and protective barrier of the rotor chamber, the type of overspeed detection system, and the degree of updating the instruments have, if any. Preparative ultracentrifuges should have a decal above the rotor chamber opening on top of the instrument or on the chamber door, indicating their classification letter. Beckman Coulter rotors are then specified for use in particular instrument classes.

In June, 1984, a major reclassification program was established to ensure continued safety to users of older ultracentrifuges and/or rotors. This reclassification of instruments and rotors is outlined below. It is essential that you use this program to determine which rotors may be safely run in which instruments. (Rotors in parentheses are no longer manufactured.)



Rotors without mechanical overspeed devices should not be used in ultracentrifuges classified other than G, H, R, or S.

Instrument Classification	Rotors that may be Used in this Instrument^a
All Model L's, classified "A"	(Type 40), (Type 40.2), (Type 40.3), (SW 50.1), (SW 25.1), and (A1-15).
All Model L's, classified "B"	(Type 50 Ti), (Type 50.3 Ti), (Type 50), (Type 40), (Type 40.2), (Type 40.3), (SW 50.1), (SW 30), (SW 30.1), (SW 25.1), and zonals.
All Model L2-50's, classified "C"	(Type 50 Ti), (Type 50.3 Ti), (Type 50), (Type 40), (Type 40.2), (Type 40.3), Type 25, (Type 15), SW 41 Ti, (SW 30), (SW 30.1), SW 28, SW 28.1, (SW 27), (SW 27.1), (SW 25.1), (SW 25.2), and zonals.
All Model L2-50's, classified "D"	(Type 50 Ti), (Type 50.3 Ti), (Type 50), (Type 40), (Type 40.2), (Type 40.3), Type 25, (Type 15), (SW 50.1), SW 41 Ti, (SW 30), (SW 30.1), SW 28, SW 28.1, (SW 27), (SW 27.1), (SW 25.1), (SW 25.2), and zonals.
All Model L2-50's, classified "F"	(Type 50 Ti), Type 50.2 Ti, (Type 50.3 Ti), Type 50.4 Ti, (Type 50), Type 45 Ti, (Type 40), (Type 40.2), (Type 40.3), Type 25, (Type 15), (SW 50.1), SW 41 Ti, (SW 30), (SW 30.1), SW 28, SW 28.1, (SW 27), (SW 27.1), (SW 25.1), (SW 25.2), and zonals.
All Model L2-65's, classified "D"	(Type 50 Ti), (Type 50.3 Ti), (Type 50), (Type 40), (Type 40.2), (Type 40.3), Type 25, (Type 15), (SW 50.1), SW 41 Ti, (SW 30), (SW 30.1), SW 28, SW 28.1, (SW 27), (SW 27.1), (SW 25.1), (SW 25.2), and zonals.

Instrument Classification	Rotors that may be Used in this Instrument^a
All Model L2-65's, classified "F"	(Type 50 Ti), Type 50.2 Ti, (Type 50.3 Ti), (Type 50), Type 45 Ti, (Type 40), (Type 40.2), (Type 40.3), Type 25, (Type 15), (SW 50.1), SW 41 Ti, SW 28, SW 28.1, (SW 27), (SW 27.1), (SW 25.1), (SW 25.2), and zonals.
All Model L2-65B's and Model L2-75B's, classified "G"	(Type 50 Ti), Type 50.2 Ti, (Type 50.3 Ti), (Type 50), Type 45 Ti, (Type 40), (Type 40.2), (Type 40.3), Type 25, (SW 50.1), SW 41 Ti, SW 40 Ti, (SW 30), (SW 30.1), SW 28, SW 28.1, (SW 27), (SW 27.1), (SW 25.1), (SW 25.2), and zonals.
All Model L3-40's and Model L3-50's, classified "F"	(Type 50 Ti), Type 50.2 Ti, (Type 50.3 Ti), (Type 50), Type 45 Ti, (Type 40), (Type 40.2), (Type 40.3), Type 25, (SW 50.1), SW 41 Ti, (SW 30), (SW 30.1), SW 28, SW 28.1, (SW 27), (SW 27.1), (SW 25.1), (SW 25.2), and zonals.
All Model L3-40's and Model L3-50's, classified "G"	(Type 50 Ti), Type 50.2 Ti, (Type 50.3 Ti), (Type 50), Type 45 Ti, Type 42.2 Ti, (Type 40), (Type 40.2), (Type 40.3), Type 25, (SW 50.1), SW 41 Ti, (SW 30), (SW 30.1), SW 28, SW 28.1, (SW 27), (SW 27.1), (SW 25.1), (SW 25.2), and zonals.
Model L4's, classified "Q"	(Type 50 Ti), (Type 50.3 Ti), (Type 50), Type 45 Ti, (Type 40), (Type 40.2), (Type 40.3), (SW 50.1), (SW 30), (SW 30.1), (SW 25.1), and zonals.
Model L5's, L5B's, L8's, and L8M's, all classified "H"	Any Beckman Coulter preparative rotor (including zonal and continuous flow rotors) EXCEPT the following: (a) all (Type 15) rotors and (b) all (Type 35) and (Type 42.1) rotors with serial numbers 1299 or lower (see Special Action below). (Type 16) and (Type 28) rotors in Model L8's and L8M's only.
Model L7's and Optima L's, LE's all classified "R"	Any Beckman Coulter preparative rotor EXCEPT the (Type 15) rotor and zonal and continuous flow rotors.
Optima XL's, L-XP's, and Optima XE, XPN classified "S"	Any Beckman Coulter preparative rotor, including zonal and continuous flow rotors.

a. To the maximum speed of the ultracentrifuge as applicable.

Special Action on Older Type 35 and Type 42.1 Rotors

We have found that there is a high risk associated with Type 35 rotor and Type 42.1 rotors having serial numbers 1299 and lower. These rotors were originally stamped "Type 42" or "Type 50.2" and were derated over 15 years ago. THESE ROTORS ARE NOW OVER 20 YEARS OLD AND MUST BE RETIRED IMMEDIATELY, REGARDLESS OF THE INSTRUMENTS IN WHICH THEY ARE USED.

Introduction

This chapter is an introduction to the Beckman Coulter family of preparative ultracentrifuge rotors, providing general information on rotor design, selection, and operation. Rotor designs described are fixed angle, swinging bucket, vertical tube, and near vertical tube type. Specific instructions for using each type of rotor are contained in [CHAPTER 4](#), [CHAPTER 5](#) and [CHAPTER 6](#). Care and maintenance information for all of these rotors is contained in [CHAPTER 7](#). Analytical, continuous flow, and zonal rotors are not covered in this manual; they are described in detail in their respective rotor instruction manuals.

General Description

Rotor Designations

Beckman Coulter preparative rotors are named according to the type of rotor, the material composition, and the rotor's maximum allowable revolutions per minute (in thousands), referred to as rated speed. For example, the SW 28 is a swinging bucket rotor with a maximum speed of 28,000 RPM. Decimal units that are sometimes part of the rotor name, as in the Type 50.2 Ti and the Type 50.4 Ti, make it possible to distinguish between different rotors that have the same maximum allowable speed. An example of each rotor type is shown in [Figure 1.1](#).

Tubes in *fixed-angle rotors* (designated **T***ype*) are held at an angle to the axis of rotation in numbered tube cavities. The bodies of some large, heavy rotors are fluted to eliminate unnecessary weight and minimize stresses.

In *swinging-bucket rotors* (designated **SW**), containers are held in rotor buckets or attached to the rotor body by hinge pins or a crossbar. The buckets swing out to a horizontal position as the rotor accelerates, then seat against the rotor body for support.

In *vertical-tube rotors* (designated **V**), tubes are held parallel to the axis of rotation. These rotors (and the near-vertical tube rotors) have plugs, screwed into the rotor cavities over sealed tubes, that restrain the tubes in the cavities and provide support for the hydrostatic forces generated by centrifugation.

Figure 1.1 Fixed-Angle, Swinging-Bucket, Vertical-Tube, and Near-Vertical Tube Rotors



Tubes in *near-vertical tube rotors* (designated **NV**), are also held at an angle to the axis of rotation in numbered tube cavities. However, the reduced tube angle of these rotors (typically 7 to 10 degrees) reduces run times from fixed-angle rotors (with tube angles of 20 to 45 degrees) while allowing components that do not band under separation conditions to either pellet to the bottom or float to the top of the tube. As in vertical-tube rotors, rotor plugs are used in these rotors to restrain the tubes in the cavities and provide support for the hydrostatic forces generated by centrifugation.

Materials

Beckman Coulter rotors are made from either aluminum or titanium, or from fiber-reinforced composites. A titanium rotor is designated by T or Ti, as in the Type 100 Ti, the SW 55 Ti, or the NVT 90 rotor. Rotors without the T or Ti designation (such as the Type 25) are fabricated from an aluminum alloy. Titanium rotors are stronger and more chemical resistant than the aluminum rotors.


Exterior surfaces of titanium are finished with black polyurethane paint. Titanium buckets and lids of high-performance rotors are usually painted red for identification.

On some swinging-bucket rotors a solid film lubricant coating is added to the bucket flange where the bucket contacts the rotor body. The purpose of the coating, which is a dull gray in color, is to minimize friction and enable the bucket to swing into the rotor bucket pocket more smoothly. With use and handling, all or part of this coating may wear off; this should not affect the rotor performance, as the bucket swing-up will wear in with use.

Aluminum rotors are anodized to protect the metal from corrosion. The anodized coating is a thin, tough layer of aluminum oxide formed electrochemically in the final stages of rotor fabrication. A colored dye may be applied over the oxide for rotor identification.

The O-rings or gaskets in fixed-angle rotor assemblies or lids, and in swinging-bucket caps, are usually made of Buna N elastomer and maintain atmospheric pressure in the rotor if they are kept clean and lightly coated with silicone vacuum grease. Plug gaskets in vertical tube or near-vertical tube rotors are made of Hytrel and do not require coating.


Drive Pins

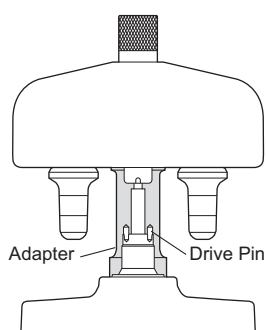

Relatively light rotors have drive pins in the drive hole that mesh with pins on the ultracentrifuge drive hub when the rotor is installed to ensure that the rotor does not slip on the hub during initial acceleration. (Heavier rotors do not require the use of drive pins.) For swinging-bucket rotors, an indentation on the rotor adapter or the position of the mechanical overspeed cartridges (see [Overspeed Protection](#) below) indicates the location of the drive pins. In this way, the pins can be properly aligned without lifting the rotor and dislocating the buckets.


Rotor Selection

Selection of a rotor depends on a variety of conditions, such as sample volume, number of sample components to be separated, particle size, desired run time, desired quality of separation, type of separation, and the centrifuge in use. Fixed-angle, swinging-bucket, vertical-tube, and near-vertical tube rotors are designed to provide optimal separations for a variety of sample types. (For especially large sample volumes, continuous flow and zonal rotors are available.)

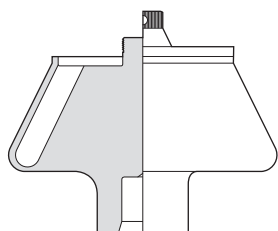


Fixed-angle rotors are general-purpose rotors that are especially useful for pelleting subcellular particles and in short-column banding of viruses and subcellular organelles. Tubes are held at an angle (usually 20 to 45 degrees) to the axis of rotation in numbered tube cavities. The tube angle shortens the particle pathlength (see [Figure 1.2](#)), compared to swinging-bucket rotors, resulting in reduced run times. Refer to [CHAPTER 4](#) for specific information about the use of fixed-angle rotors.


Figure 1.2 Particle Separation in Fixed-Angle, Swinging-Bucket, Vertical-Tube, and Near-Vertical Tube Rotors*

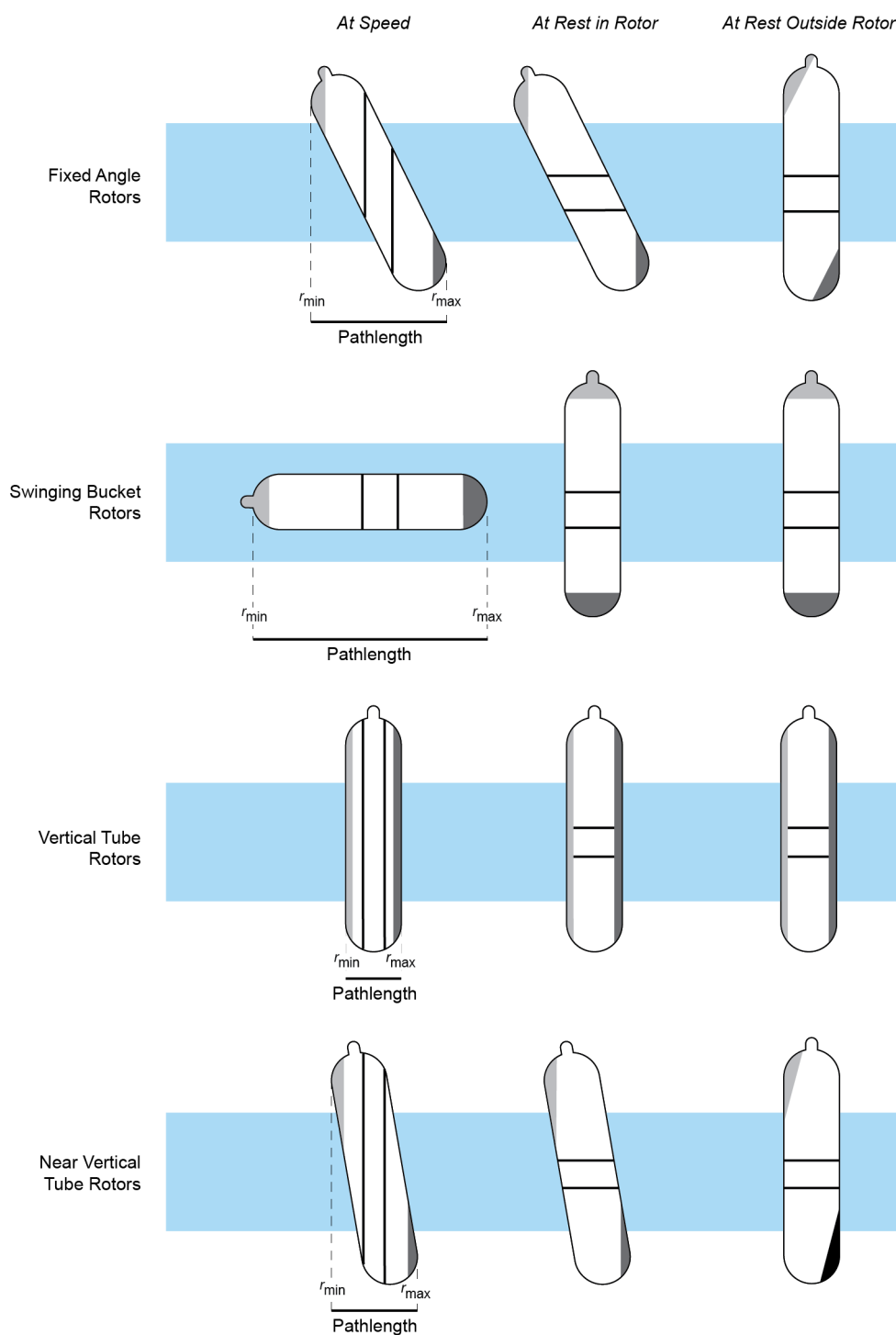


* Dark gray represents pelleted material, light gray is floating components, and bands are indicated by black lines.

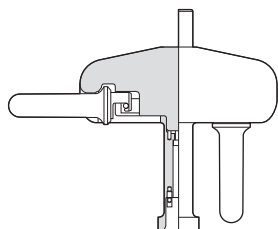

Swinging-bucket rotors are used for pelleting, isopycnic studies (separation as a function of density), and rate zonal studies (separation as a function of sedimentation coefficient). Swinging-bucket rotors are best applied for rate zonal studies in which maximum resolution of sample zones are needed, or pelleting runs where it is desirable for the pellet to be in the exact center of the tube bottom. Gradients of all shapes and steepness can be used. Refer to [CHAPTER 5](#) for specific information about the use of swinging-bucket rotors.

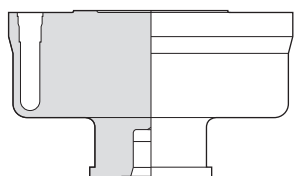

Vertical-tube rotors hold tubes parallel to the axis of rotation; therefore, bands separate across the diameter of the tube rather than down the length of the tube (see [Figure 1.2](#)). Vertical-tube rotors are useful for isopycnic and, in some cases, rate zonal separations when run time reduction is important. Only Quick-Seal and OptiSeal tubes are used in vertical-tube rotors, making tube caps unnecessary. Refer to [CHAPTER 6](#) for specific information about the use of vertical-tube rotors.

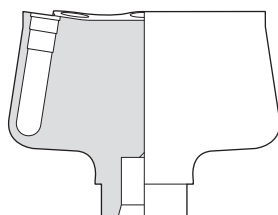

Near-vertical tube rotors are designed for gradient centrifugation when there are components in a sample mixture that do not participate in the gradient. The reduced tube angle of these rotors significantly reduces run times from the more conventional fixed-angle rotors, while allowing components that do not band under separation conditions to either pellet to the bottom or float to the top of the tube. Like the vertical-tube rotors, near-vertical tube rotors use only Quick-Seal and OptiSeal tubes. Refer to [CHAPTER 6](#) for specific information about the use of near-vertical tube rotors.

[Table 1.1](#) lists Beckman Coulter preparative rotors by use.

Table 1.1 Beckman Coulter Preparative Rotors by Use

Rotor	Maximum Speed ^a (rpm)	Relative Centrifugal Field ^b ($\times g$) at $r_{\max}$	k Factor	Number of Tubes $\times$ Nominal Capacity (mL) of Largest Tube	Nominal Rotor Capacity (mL)	For Use in Instruments Classified
Rotors for Centrifuging Extremely Small Particles						
NVT 100	100,000	750,000	8	8 $\times$ 5.1	40.8	R, S
Type 100 Ti	100,000	801,920	15	8 $\times$ 6.8	54	R, S
NVT 90	90,000	645,000	10	8 $\times$ 5.1	40.8	H,R,S
Type 90 Ti	90,000	694,000	25	8 $\times$ 13.5	108	H,R,S
VTi 90	90,000	645,000	6	8 $\times$ 5.1	40.8	H,R,S
NVT 65.2	65,000	416,000	16	16 $\times$ 5.1	81.6	H,R,S
NVT 65	65,000	402,000	21	8 $\times$ 13.5	108	H,R,S
VTi 65.2	65,000	416,000	10	16 $\times$ 5.1	81.6	H,R,S
VTi 65.1	65,000	402,000	13	8 $\times$ 13.5	108	H,R,S
Rotors for Centrifuging Small Particles in Volume						
Type 70 Ti	70,000	504,000	44	8 $\times$ 38.5	308	G ^c ,H,R,S
Type 50.2 Ti	50,000	302,000	69	12 $\times$ 39	468	F,G ^c ,H,R,S
VTi 50.1	50,000	251,000	34	12 $\times$ 39	468	R, S
VTi 50	50,000	242,000	36	8 $\times$ 39	312	H,R,S
Type 45 Ti	45,000	235,000	133	6 $\times$ 94	564	F,G ^c ,H,Q,R,S
Rotors for Differential Flotation						
Type 50.4 Ti	50,000	312,000 ^d	33	44 $\times$ 6.5	286	G ^c ,H,R,S
Type 42.2 Ti	42,000	223,000	12	72 $\times$ 230 μ L	16.5	G ^c ,H,R,S
Type 25	25,000	92,500 ^e	62	100 $\times$ 1	100	C,D,F,G,H,R,S
Rotors for Centrifuging Large Particles						
Type 70.1 Ti	70,000	450,000	36	12 $\times$ 13.5	162	G ^c ,H,R,S
Rotors for Centrifuging Large Particles in Volume						
Type 19	19,000	53,900	951	6 $\times$ 250	1500	H,R,S
Rotors for Isopycnic and Rate-Zonal Gradients						
SW 60 Ti	60,000	485,000	45	6 $\times$ 4	24	G ^c ,H,R,S

Table 1.1 Beckman Coulter Preparative Rotors by Use (Continued)

Rotor	Maximum Speed ^a (rpm)	Relative Centrifugal Field ^b ($\times g$) at $r_{\max}$	k Factor	Number of Tubes $\times$ Nominal Capacity (mL) of Largest Tube	Nominal Rotor Capacity (mL)	For Use in Instruments Classified
SW 55 Ti	55,000	368,000	48	6 $\times$ 5	30	G ^c ,H,R,S
Rotors with Long, Slender Tubes for Rate-Zonal Gradients						
SW 41 Ti	41,000	288,000	124	6 $\times$ 13.2	79.2	C,D,F,G,H,R,S
SW 40 Ti	40,000	285,000	137	6 $\times$ 14	84	G ^c ,H,R,S
SW 32 Ti	32,000	175,000	204	6 $\times$ 38.5	231	H, R, S
SW 28.1	28,000	150,000	276	6 $\times$ 17	102	C,D,F,G,H,R,S
Rotors for Larger-Volume Density Gradients						
SW 32.1	32,000	187,000	228	6 $\times$ 17	102	H, R, S
Rotors for Larger-Volume Density Gradients (continued)						
SW 28	28,000	141,000	245	6 $\times$ 38.5	231	C,D,F,G,H,R,S

- a. Maximum speeds are based on a solution density of 1.2 g/mL in all rotors except for the near-vertical tube and vertical-tube rotors, which are rated for a density of 1.7 g/mL.
- b. Relative Centrifugal Field (RCF) is the ratio of the centrifugal acceleration at a specified radius and speed ($r\omega^2$) to the standard acceleration of gravity (g) according to the following formula: $RCF = r\omega^2/g$ where r is the radius in millimeters, ω is the angular velocity in radians per second ($2\pi\text{RPM}/60$), and g is the standard acceleration of gravity (9807 mm/s^2). After substitution: $RCF = 1.12r (\text{RPM}/1000)^2$.
- c. Class G, Model L3 only.
- d. Maximum RCF measured at outer row.
- e. Maximum RCF measured at the third (outermost) row. Radial distances are those of the third row.

Pelleting (Differential Separation)

Pelleting separates particles of different sedimentation coefficients, the largest particles in the sample traveling to the bottom of the tube first. Differential centrifugation is the successive pelleting of particles of decreasing sedimentation velocities, using increasingly higher forces and/or long run times. The relative pelleting efficiency of each rotor is measured by its k factor (clearing factor):

EQ 1

$$k = \frac{\ln(r_{\max}/r_{\min})}{\omega^2} \times \frac{10^{13}}{3600}$$

where ω is the angular velocity of the rotor in radians per second ($2\pi\text{RPM}/60$, or $\omega = 0.10472 \times \text{rpm}$), $r_{\max}$ is the maximum radius, and $r_{\min}$ is the minimum radius.

After substitution,

EQ 2

$$k = \frac{(2.533 \times 10^{11}) \ln(r_{\max}/r_{\min})}{\text{rpm}^2}$$

This factor can be used in the following equation to estimate the time t (in hours) required for pelleting:

EQ 3

$$t = \frac{k}{s}$$

where s is the sedimentation coefficient* of the particle of interest in Svedberg units. (Because s values in seconds are such small numbers, they are generally expressed in Svedberg units (S), where 1 S is equal to 10^{-13} seconds). It is usual practice to use the standard sedimentation coefficient $s_{20,w}$ based on sedimentation in water at 20°C. Clearing factors can be calculated at speeds other than maximum rated speed by use of the following formula:

EQ 4

$$k_{\text{adj}} = k \left(\frac{\text{rated speed of rotor}}{\text{actual run speed}} \right)^2$$

Run times can also be calculated from data established in prior experiments when the k factor of the previous rotor is known. For any two rotors, a and b :

EQ 5

$$\frac{t_a}{t_b} = \frac{k_a}{k_b}$$

where the k factors have been adjusted for the actual run speed used.


Figure 1.3 lists sedimentation coefficients for some common biological materials. The k factors at maximum speeds for Beckman Coulter preparative rotors are provided in the table of general specifications in each rotor use section.

* $1 s = dr/dt \times 1/w2r$, where dr/dt is the sedimentation velocity.

The centrifugal force exerted at a given radius in a rotor is a function of the rotor speed. The nomogram in [Figure 1.4](#) allows you to determine relative centrifugal field (RCF) for a given radius and rotor speed.

Run times can be shortened (in some rotors) by using the *g*-Max system. The short pathlength means less distance for particles to travel in the portion of the tube experiencing greatest centrifugal force, and hence shortened run times. Run times can also be shortened (in some rotors) by using partially filled thickwall polypropylene and polycarbonate tubes. The *k* factors for half-filled tubes can be calculated by using an approximate $r_{\max}$ and r_{av} in *k*-factor equation (1).

Figure 1.3 Sedimentation Coefficients (in Svedberg Units) for Some Common Biological Materials

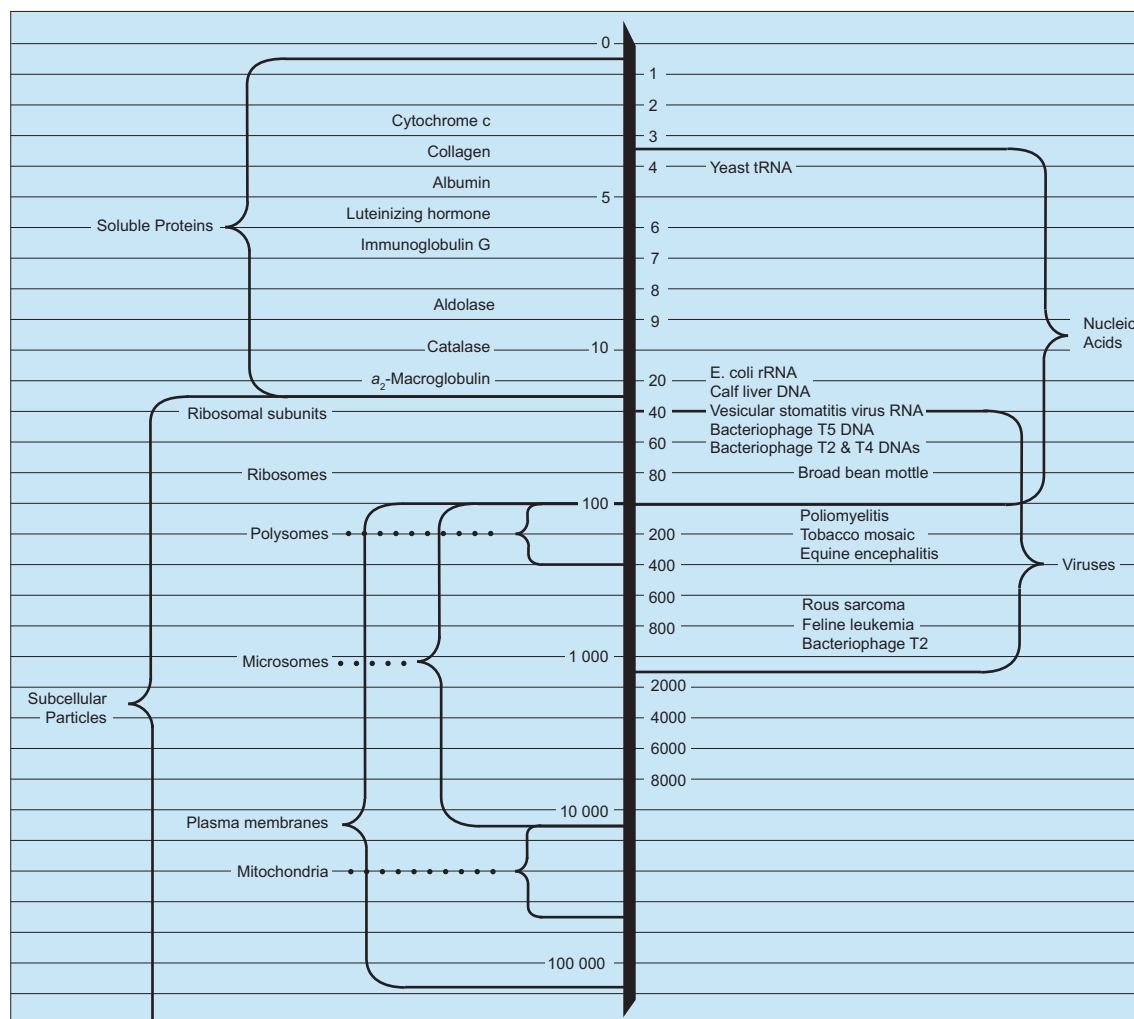


Isopycnic Separations

A sedimentation-equilibrium, or isopycnic, method separates particles on the basis of particle buoyant density. Each component in the sample travels through the gradient until it reaches an equilibrium position. Particle velocity due to differences in density is given in the following expression:

EQ 6

$$v = \left[\frac{d^2(\rho_p - \rho_c)}{18\mu} \right] \times g$$

Figure 1.4 Nomogram*

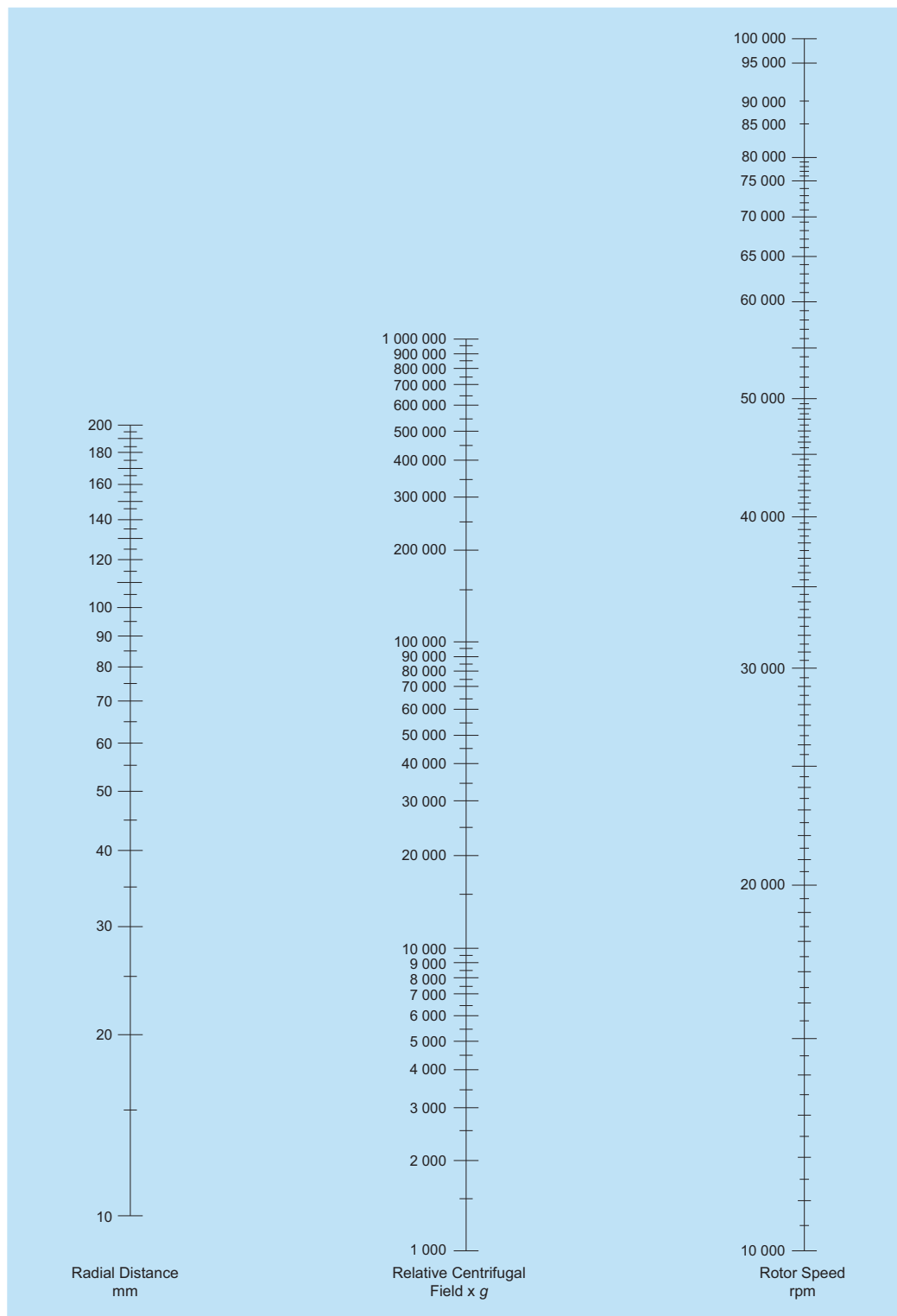

where

* Align a straightedge through known values in two columns; read the figure where the straightedge intersects the third column.

- v = sedimentation velocity (dr/dt)
- d = particle diameter
- ρ_p = particle density
- ρ_c = solution density
- μ = viscosity of liquid media
- g = standard acceleration of gravity

At equilibrium, $\rho_p - \rho_c$ is zero, and particle velocity is therefore zero.

The gradient may be preformed before the run or generated during centrifugation. For gradients formed by centrifugation, the time it takes to form a gradient depends on the sedimentation and diffusion coefficients of the gradient material, the pathlength, and the rotor speed. For a given gradient material, the shorter the pathlength and the higher the rotor speed, the faster the gradient will form. In general, the time required for gradients to reach equilibrium in swinging bucket rotors will be longer than in fixed angle rotors. One way to reduce run times is to use partially filled tubes. Refer to the appropriate rotor instruction manual to determine the maximum allowable speed and solution density when using partially filled tubes.

Rate Zonal Separations

Particle separation achieved with rate zonal separation is a function of the particles' sedimentation coefficient (density, size, and shape) and viscosity of the gradient material. Sucrose is especially useful as a gradient material for rate zonal separation because its physical characteristics are well known and it is readily available. Samples are layered on top of the gradient. Under centrifugal force, particles migrate as zones. Rate zonal separation is time dependent; if the particles are more dense than the most dense portion of the gradient, some or all of the particles will pellet unless the run is stopped at the appropriate time.

A separation is sometimes a combination of rate zonal and isopycnic. Depending on particle buoyant densities and sedimentation coefficients, some particles may be separated by their differential rates of sedimentation, while others may reach their isopycnic point in the gradient.

Clearing factors of swinging-bucket rotors at maximum speeds and various particle densities have been calculated for 5 to 20% (wt/wt) linear sucrose gradients at 5°C. These are called k' factor, and are given in [Table 5.1](#) in [CHAPTER 5](#). These constants can be used to estimate the time, t (in hours), required to move a zone of particles of known sedimentation coefficient and density to the bottom of a 5 to 20% gradient:

EQ 7

$$t = \frac{k'}{s}$$

where s is the sedimentation coefficient in Svedberg units, S .

In most cases, when banding two or three components by rate zonal separation, run times can be considerably reduced by using reduced fill levels. Tubes are partially filled with gradient, but the sample volume is not changed (however, gradient capacity will be reduced). Thickwall tubes should be used when this technique is employed, since thinwall tubes will collapse if not full.

If swinging bucket rotors are used with preformed shallow gradients (< 5 to 20%), or if fixed angle, vertical tube, or near-vertical tube rotors are used with any preformed gradient, use the slow acceleration control on your ultracentrifuge. Slow acceleration will protect the sample-to-gradient interface, and slow deceleration will maintain the integrity of the separation during the reorientation process.

General Operating Information


Careful centrifugation technique is essential, because forces generated in high-speed centrifugation can be enormous. For example, 1 gram at the bottom of an SW 60 Ti rotor bucket, rotating at 60,000 rpm, exerts the gravitational equivalent of 0.5 ton of centrifugal mass at the bottom of the bucket.

Note the classification letter of the ultracentrifuge to be used, and be sure the rotor is appropriate for the instrument (see the [Classification Program](#) chart at the beginning of this manual and [Table 1.1](#)).

Acceptable classification letters are engraved on rotor lids, handles, stands, or bodies.

NOTE Specific information about filling, sealing, and capping containers, loading rotors, etc., can be found in later sections.

Rotor Balance

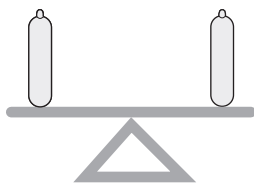


The mass of a properly loaded rotor will be evenly distributed on the ultracentrifuge drive hub, causing the rotor to turn smoothly with the drive. An improperly loaded rotor will be unbalanced; consistent running of unbalanced rotors will reduce ultracentrifuge drive life. To balance the rotor load, fill all opposing tubes to the same level with liquid of the same density. Weight of opposing tubes must be distributed equally. Place tubes in the rotor symmetrically, as illustrated in [Figure 1.5](#).



CAUTION

For swinging bucket rotors, attach ALL buckets, whether loaded or empty. For vertical tube and near-vertical tube rotors, insert spacers and rotor plugs ONLY in holes containing loaded tubes.

Figure 1.5 Arranging Tubes Symmetrically in a Rotor*

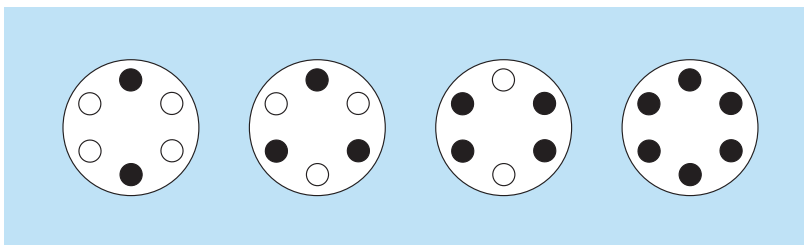

If sample quantity is limited and the rotor is not balanced, do one of the following to balance the rotor, depending on the rotor in use:

- Load the opposite rotor cavities or buckets with tubes containing a liquid of the same density as opposing tubes.
- Use smaller tubes with adapters or smaller Quick-Seal tubes with floating spacers to distribute the sample symmetrically.
- Use thickwall tubes partially filled to distribute sample to additional tubes.
- Layer a low-density, immiscible liquid, such as mineral oil, on top of the sample to fill opposing tubes to the same level. (Do not use an oil overlay in Ultra-Clear tubes.)


Overspeed Protection

Rotors are specifically designed to withstand a maximum load (that is, volume and density of the rotor contents) at maximum rated speed. At greater speeds, or at rated speeds with heavier loads, rotors are subject to failure. It is the operator's responsibility to limit rotor speed when centrifuging dense solutions or when using heavy tubes; refer to [Allowable Run Speeds](#), below.


Rotors are protected from exceeding their maximum rated speed to help prevent failure and damage to the rotor and the instrument. Two overspeed protection systems are used in Beckman Coulter preparative ultracentrifuges.

- Optima L and LE (classified R), Optima XL and L-XP (classified S), Optima XE and XPN, and Model L7 (classified R), have a photoelectric overspeed system. This system includes a photoelectric device in the rotor chamber next to the drive hub and an overspeed disk on the rotor bottom.
- Earlier model ultracentrifuges — classified other than G, H, R, or S (and some F) — have a mechanical overspeed system.

* For example, two, three, four, or six tubes can be arranged symmetrically in a six-place rotor.

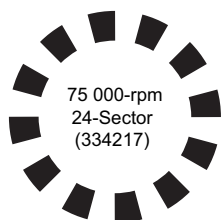

All Beckman Coulter preparative rotors are shipped with an overspeed disk attached, and are therefore protected from overspeeding in instruments with the photoelectric system. These instruments will not operate unless an overspeed disk is attached to the installed rotor. The disk has alternating sectors of reflecting and nonreflecting material. The number of sectors on the disk is a function of the rotor's maximum allowable speed. During centrifugation, if the reflective segments pass over the photoelectric pickup faster than the indicated set speed, the drive will automatically decelerate to the allowed speed.

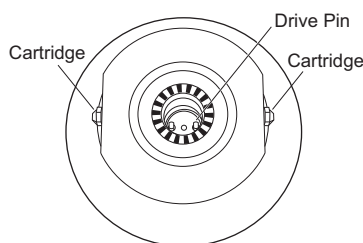

The earlier model ultracentrifuges—classified A, B, C, D, N, O, P, Q, and some F—with the mechanical overspeed system have a knockout pin in the rotor chamber. Rotors that are equipped for the mechanical system have overspeed cartridges installed in the sides of the rotor base. If overspeeding occurs, a small pin is forced out of the cartridge and knocks out the overspeed pin in the chamber, causing the instrument to shut down.

CAUTION

Rotors without mechanical overspeed cartridges should not be used in ultracentrifuges classified other than G, H, R, or S.

The overspeed device should be replaced if a rotor is regularly being used at speeds below its rated speed due to the use of adapters, stainless steel tubes, CsCl gradients, etc. Instructions for replacing overspeed disks are provided in Section 7 of this manual.

Allowable Run Speeds

Under some conditions, the maximum allowable speed of the rotor (indicated by the rotor name) must be reduced to ensure that neither the rotor nor the labware are overstressed during centrifugation. Check the recommended run speed for your rotor before centrifuging dense solutions, CsCl gradients, stainless steel tubes, polycarbonate bottles, uncapped plastic tubes in fixed angle rotors, and sleeve-type adapters.

- **Dense Solutions.** To protect the rotor from excessive stresses due to the added load, reduce run speed when centrifuging a solution with a density greater than the allowable density rating of the rotor (specified in the rotor instruction manual). When using dense solutions in plastic labware, determine maximum run speed using the following square-root reduction formula:

$$\text{reduced run speed} = \text{maximum rated speed} \sqrt{\frac{A}{B}}$$

- where A is the maximum permissible density of the tube contents for a particular rotor (from the rotor instruction manual), and B is the actual density of the tube contents to be centrifuged.
- When using dense solutions in *stainless steel tubes*, refer to the individual rotor instruction manual or *Run Speeds for Stainless Steel Tubes* (publication L5-TB-072) for allowable speeds.
- *Cesium Chloride Gradients*. Run speed often must be reduced to avoid the precipitation of CsCl during centrifugation of concentrated CsCl solutions. Use the CsCl curves provided in the individual rotor instruction manual to determine run speeds. An example of the use of CsCl curves is in [APPENDIX B](#) of this manual.
- *Uncapped Thickwall Plastic Tubes in Fixed-Angle Rotors*. Speed limitations are required to prevent tube collapse when thickwall plastic tubes are centrifuged without the support of tube caps in fixed-angle rotors (refer to [CHAPTER 4](#)).
- *Polycarbonate and Polypropylene Bottles*. Speed limitations are required to prevent the bottle material from overstressing and deforming (refer to [CHAPTER 2](#)).
- *Adapters*. When small tubes are used with Delrin adapters, run speed often must be reduced due to the increased density of Delrin (1.4 g/mL). The formula for speed reduction is described in [CHAPTER 2](#). Consult individual rotor manuals for allowable run speeds.
- *Stainless Steel Tubes*. Reduce run speed when centrifuging stainless steel tubes to prevent the rotor from overstressing due to the added weight. The criteria for speed reduction percentage depends on the tube-cap material and the strength of the rotor in use; consult the individual rotor manual or publication L5-TB-072.

CHAPTER 2

Tubes, Bottles, and Accessories

Introduction

This chapter describes various labware used in Beckman Coulter preparative rotors. General instructions for using containers follow in [CHAPTER 3](#). Care and maintenance instructions are in [CHAPTER 7](#). General rotor use instructions are in [CHAPTER 4](#), [CHAPTER 5](#), and [CHAPTER 6](#). The individual rotor manual provides specific instructions on the tubes, bottles, and accessories that can be used in a particular rotor. A table of chemical resistances can be found in IN-175.*

Labware Selection Criteria

No single tube or bottle design or material meets all application requirements. Labware choice is usually based on a number of factors.

- The centrifugation technique to be used, including the rotor in use, volume of sample to be centrifuged, importance of band visibility, and so forth
- Chemical resistance—the nature of the sample and any solvent or gradient media
- Sterility
- Certified Free of Contaminates
- Temperature and speed considerations
- Whether tubes or bottles are to be reused

[Table 2.1](#) contains an overview of some of the characteristics of tube and bottle materials.

* A complete list of tubes, bottles, and accessories is provided in the latest edition of the Beckman Coulter Ultracentrifuge Rotors, Tubes & Accessories catalog (BR-8101), available at www.beckman.com.

Table 2.1 Characteristics and Chemical Resistances of Tube and Bottle Materials^a

Tube or Bottle Type	Optical Property	Puncturable	Sliceable	Reusable ^b	Acids (dilute or weak)	Acids (strong)	Alcohols (aliphatic)	Aldehydes	Bases	Esters	Hydrocarbons (aliphatic)	Hydrocarbons (aromatic and halogenated)	Ketones	Oxidizing Agents (strong)	Salts
thinwall polypropylene	transparent	yes	yes	no	S	U	U	M	S	U	U	U	U	U	S
thickwall polypropylene	translucent	no	no ^c	yes	S	S	S	M	S	M	M	U	M	U	S
Ultra-Clear	transparent	yes	yes	no	S	U	U	S	U	U	U	U	U	U	M
polycarbonate	transparent	no	no	yes	M	U	U	M	U	U	U	U	U	M	M
polypropylene	translucent/transparent	no	no ^b	yes	S	S	S	M	S	M	S	M	M	M	S
polyethylene	transparent/translucent	yes	no	yes	S	S	S	S	S	S	U	M	M	M	S
cellulose propionate	transparent	no	no ^b	no	S	U	U	U	U	M	S	S	U	M	S
stainless steel	opaque	no	no	yes	S	U	S	S	M	S	S	S	M	S	M

S = satisfactory resistance M = marginal resistance U = unsatisfactory resistance

a. Refer to IN-175 for information about specific solutions.

b. Inspect all tubes for cracking, crazing, stretching, or other damage. Do not reuse tubes that have these defects.

c. Polypropylene, polypropylene, and cellulose propionate tubes with diameters of 5 to 13 mm may be sliced using the Centritube Slicer (part number 347960) and appropriate adapter plate.

NOTE This information has been consolidated from a number of sources and is provided only as a guide to the selection of tube or bottle materials. Soak tests at 1 g (at 20°C) established the data for most of the materials; reactions may vary under the stress of centrifugation, or with extended contact or temperature variations. To prevent failure and loss of valuable sample, ALWAYS TEST SOLUTIONS UNDER OPERATING CONDITIONS BEFORE USE.



WARNING

Do not use flammable substances in or near operating centrifuges.

Labware Material Compatibility with Solvents and Sample

The chemical compatibility of tube or bottle materials with the gradient-forming medium or other chemicals in the solution is an important consideration. Although neutral sucrose and salt solutions cause no problems, alkaline solutions cannot be used in Ultra-Clear tubes or in polycarbonate tubes and bottles. Polycarbonate and Ultra-Clear tubes are incompatible with DMSO, sometimes used in the preparation of sucrose gradients for sedimentation of denatured DNA.

Gradient Formation and Fractionation

Consideration should be given to gradient formation and fractionation when choosing a tube for a density gradient run. If the bands or zones formed during centrifugation are indistinct, they may not be visible through a translucent material such as polypropylene. If optimum band visualization is important, Ultra-Clear, polycarbonate, or cellulose propionate tubes should be used. Whenever collection of bands or zones must be done by slicing or puncturing the tube, a thin, flexible tube wall is required. Ultra-Clear or polypropylene tubes should be used in these cases, depending on the need for transparency.


Certified Free Tubes



Based on Sample
Results Below
Detectable Limit

Certified free tubes are lot traceable to testing that confirms the absence of endotoxin, DNase, RNase, and human & mouse DNA below a detectable limit.

Sterile Tubes

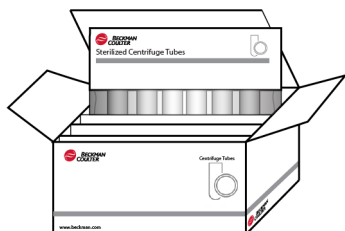

Sterile tubes are sterilized via ethylene oxide in compliance with ISO 11135. Cartons include several peel packages, each containing a typical run quantity of tubes. Packaging meets requirements of ISO 11607.

Labware Types

Polypropylene Tubes

Polypropylene tubes are translucent or transparent in appearance, depending on wall thickness, and are non-wettable (although some polypropylene tubes can be chemically treated to make them wettable). Polypropylene tubes are reusable unless deformed during centrifugation or autoclaving. Polypropylene tubes have good tolerance to gradient media, including alkalines. They are

satisfactory for many acids, bases, alcohols, DMSO, and some organic solvents. They can be used with or without caps in fixed-angle rotors. Speed reduction is sometimes required with these tubes if run with less than full volume (refer to your rotor manual). Several types of polypropylene tubes are available.

Open-Top Polypropylene Tubes




Thinwall open-top tubes are used in swinging bucket and fixed-angle rotors. In swinging-bucket rotors, thinwall tubes should be filled to within 2 or 3 mm of the tube top for proper tube support. Caps are usually required in fixed-angle rotors. Thinwall tubes are designed for one-time use and should be discarded after use.




Thickwall open-top tubes offer the convenience of centrifuging partially filled tubes without tube caps in fixed-angle and swinging-bucket rotors. Because the solution re-orientes during centrifugation, the maximum partial fill volume depends on the tube angle. For greater fill volumes, use tubes with caps. Refer to the applicable rotor manual for fill volumes and speed reduction requirements. Thickwall polypropylene tubes are typically reusable unless deformed during centrifugation or autoclaving.

OptiSeal Tubes

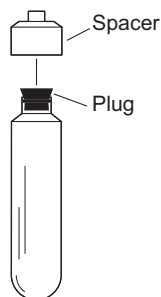

OptiSeal tubes, single-use tubes designed for use in certain rotors, are available in dome-top and bell-top styles. These tubes, which come with plastic sealing plugs, can be quickly and easily prepared for use without tools or heat. Spacers are used to seal the tubes and to support the tops of the tubes during centrifugation. With the tube plug and spacer (and rotor plug, if required) in place, the g forces during centrifugation ensure a tight, reliable seal that protects your samples. For a detailed discussion on the use of OptiSeal tubes, refer to Using OptiSeal Tubes (publication IN-189), included with each box of tubes.

Quick-Seal Polypropylene Tubes

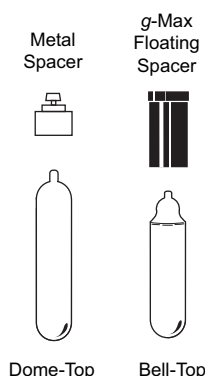

Heat-sealed Quick-Seal tubes are used in swinging bucket, vertical tube, near vertical tube, and in most fixed-angle rotors. Single-use Quick-Seal tubes are a convenient form of sealable tube; they are especially useful for the containment of radioactive or pathogenic samples. There are two Quick-Seal tube designs, dome-top and bell-top.

- The bell-top simplifies removal of materials that float during centrifugation.
- Dome-top tubes hold more volume than their bell-top equivalents.

Detailed information about Quick-Seal tubes is contained in publication IN-181.

Polycarbonate Tubes



Polycarbonate is tough, rigid, nonwetable, and glass-like in appearance. Polycarbonate tubes are used with or without caps in fixed-angle rotors, and at least half full in swinging-bucket rotors. Speed reduction may be required in some rotors if the tubes are not completely filled.

Although polycarbonate tubes may be autoclaved, doing so greatly reduces the usable life of these tubes. Cold sterilization methods are recommended. Washing with alkaline detergents can cause failure. Crazing—the appearance of fine cracks in the tube—is the result of stress “relaxation” and can affect tube performance. These cracks will gradually increase in size and depth, becoming more visible. Tubes should be discarded before cracks become large enough for fluid to escape. These tubes have good tolerance to all gradient media except alkalines (pH greater than 8). They are satisfactory for some weak acids, but are unsatisfactory for all bases, alcohol, and other organic solvents.

Polyethylene Tubes

Polyethylene tubes are translucent or transparent and have a good tolerance for use with strong acids and bases. They are reusable but cannot be autoclaved. In swinging-bucket rotors, they are used without caps, and with or without caps in fixed-angle rotors.

Ultra-Clear Tubes



Ultra-Clear tubes, made of a tough thermoplastic, are thinwall and not wettable (but can be made wettable; see [CHAPTER 3](#)). Ultra-Clear tubes are available in two types—open-top and Quick-Seal. They are transparent centrifuge tubes, offering easy location of visible banded samples. Standard straight-wall Ultra-Clear tubes must be filled completely and capped for use in fixed-angle rotors.

Ultra-Clear tubes are designed to be used one time only. These tubes have good resistance to most weak acids and some weak bases, but are unsatisfactory for DMSO and most organic solvents, including all alcohols. Ultra-Clear tubes should not be autoclaved.

Cellulose Propionate Tubes


Cellulose propionate tubes, used in some fixed-angle rotors, are transparent and designed for onetime use. They are used without caps and should be full for centrifuging. They should not be autoclaved or sterilized with alcohol. These tubes have good tolerance to all gradient media including alkalines. They are unsatisfactory for most acids and alcohols.

Stainless Steel Tubes



Stainless steel tubes offer excellent resistance to organic solvents and heat, but should not be used with most acids or bases. They offer only marginal resistance to most gradient-forming materials other than sucrose and glycerol. Stainless steel tubes are very strong and can be centrifuged when filled to any level. Because of their weight, however, run speeds must often be reduced (see publication L5-TB-072). Stainless steel tubes can be used indefinitely if they are undamaged and not allowed to corrode. They may be autoclaved after use as long as they are thoroughly dried before storage.

Konical Tubes

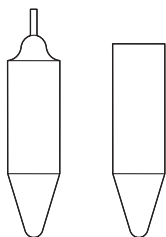

konical tubes, used with conical adapters in swinging-bucket rotors to optimize pelleting separations, have a conical tip that concentrates the pellet in the narrow end of the tube. The narrow bottom also reduces the tube's nominal volume and minimizes the amount of gradient material needed when pelleting through a dense cushion. They are available in polypropylene and Ultra-Clear. The konical tubes come in both open-top and Quick-Seal tube designs. The Quick-Seal type have bell-shaped tops to fit the floating spacers in the g-Max system for smaller volume runs with faster pelleting.

Bottles




Bottles are available in polycarbonate (hard and clear) and polypropylene (translucent).

- Threaded-top polycarbonate bottles are available for many fixed-angle rotors. They have a liquid-tight cap assembly and are easy to use. Caps (and plugs, if applicable) should always be removed before autoclaving.
- The Type 19 rotor uses a polypropylene bottle with a three-piece cap assembly consisting of a Noryl plug, a neoprene O-ring, and a Delrin cap.

Information about these bottles can be found in the individual rotor manuals.


Temperature Limits

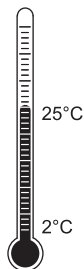

Each labware material has a specified temperature range. Although some high-speed centrifuges can achieve temperatures as high as 45°C, only certain tube or bottle materials can be run under these conditions. Most containers are made of thermoplastic materials that soften at elevated temperatures. This temperature-induced softening, together with such factors as the centrifugal force, the run duration, the type of rotor, previous run history, and the tube angle, can cause labware to collapse. Therefore, if high-temperature runs—above 25°C—are required, it is best to pretest labware under the actual experimental conditions, using buffer or gradient of similar density rather than a valuable sample. (Stainless steel tubes can be centrifuged at any temperature.)

- Plastic labware has been centrifuge tested for use at temperatures between 2 and 25°C. For centrifugation at other temperatures, pretest tubes under anticipated run conditions.
- If plastic containers are frozen before use, make sure that they are thawed to at least 2°C prior to centrifugation.

Spacers and Floating Spacers

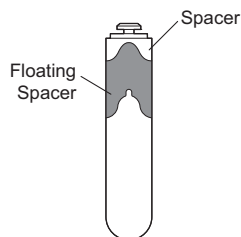


- OptiSeal tubes must be used with the appropriate spacer to seal properly. (OptiSeal spacers are listed in [Table 3.2](#).)
- Quick-Seal tubes use a spacer ([Table 2.2](#)), one or more floating spacers, or a combination of both (depending on the size of the tube) to support the top of the tube during centrifugation. The particular combination depends on the type of rotor being used. In swinging bucket and fixed-angle rotors, the top of the tube must be supported. In near-vertical tube and vertical-tube rotors, the entire tube cavity must be filled

The g-Max system uses a combination of short bell-top Quick-Seal tubes and floating spacers (also referred to as g-Max spacers). The floating spacers sit on top of the Quick-Seal tubes so there is no reduction of maximum radial distance, and therefore, no reduction of g force. The shorter pathlength of the tubes also permits shorter run times. For more information on the *g-Max Optima* Application Note at www.beckman.com/resources/reading-material/application-notes/g-max-system-for-ultracentrifugation.

Plastic spacers have been tested for centrifugation between 2 and 25°C. If spacers are centrifuged at temperatures significantly greater than 25°C, deformation of the spacer and tube may occur.

Adapters*

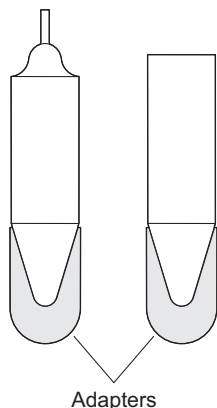



Many rotors can accommodate a variety of tube sizes by using adapters that line the tube cavity or bucket.

- Small, open-top tubes use Delrin* adapters, which line the tube cavity or bucket.
- Adapters with conical cavities must be used to support both open-top and Quick-Seal konical tubes.

Tubes used with adapters can be filled (and capped) according to the type of tube and the design of the rotor being used. Many of the small, straightwall tubes, when used with adapters, require speed reductions due to the added density of Delrin (1.4 g/mL). Additional speed reductions for heavy tube loads may also be required (refer to [Allowable Run Speeds](#) in [CHAPTER 1](#)). In vertical-tube rotors, $r_{\min}$ is unchanged (see the illustration in [Figure 1.2](#)).


However, in fixed angle and near-vertical tube rotors, $r'_{\min}$ must be calculated:

Table 2.2 Quick-Seal Tube Spacers

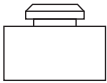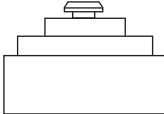Part Number	Spacer Description
342883	black-anodized aluminum 
342418	clear-anodized aluminum 
342696	clear-anodized aluminum 
342695	red-anodized aluminum 
342699	red-anodized aluminum 
342417	clear-anodized aluminum 
342697	titanium 

* Delrin is a registered trademark of E. I. Du Pont de Nemours & Company.

Table 2.2 Quick-Seal Tube Spacers

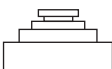Part Number	Spacer Description	
344389	white Delrin	
344634 344635	white Delrin	
344676	black Noryl	
345828	black Noryl	
349289	blue-anodized aluminum	
358164	black Delrin	

EQ 9

$$r'_{\min} = r_{\max} - \frac{d}{2} (1 - \sin \theta + \cos \theta) - L \sin \theta$$


where

$r_{\max}$ = the distance in millimeters from the axis of rotation to the farthest part of the tube cavity,

d = diameter of the tube in millimeters,

L = length of the tube in millimeters, and

θ = tube angle of the rotor being used

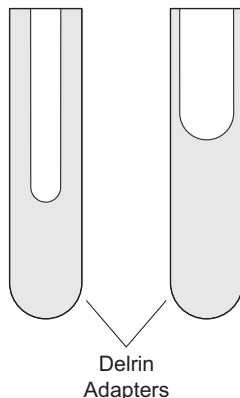

A Delrin adapter in a rotor cavity or bucket will significantly change the radial distances measured in the tube. The equations below can be used to determine $r'_{\max}$ and $r'_{\min}$ for a given rotor with a Delrin adapter. [Table 2.3](#) lists adapter dimensions used in the equations

EQ 10

$$r'_{\max} = r_{\max} - \left(\frac{d_1 - d_2}{2} \right) - \left(t - \frac{d_1 - d_2}{2} \right) \sin \theta$$

EQ 11

$$r'_{\min} = r_{\max} - \frac{d_1}{2} - \left(t - \frac{d_1}{2} + L \right) \sin \theta - \frac{d_2}{2} \cos \theta$$

where

- $r_{\max}$ = the distance in millimeters from the axis of rotation to the farthest part of the tube cavity,
- d_1 = outside diameter of the adapter,
- d_2 = inside diameter of the adapter,
- L = adapter cavity length,
- t = thickness of the adapter bottom, and
- θ = tube angle of the rotor being used

The values of $r'_{\max}$ and $r'_{\min}$ can be used to calculate the k factor and the relative centrifugal field when adapters are used (see the equations in the [Glossary](#)).

Table 2.3 Dimensions of Delrin Adapters^a

Delrin Adapter		Dimensions (mm)			
Tube Size (mL)	Part Number	d_1	d_2	L	t
0.8	305527	13.08	5.49	42.09	3.99
	356860	18.08	5.36	43.51	3.99
2	303376	16.23	8.66	46.25	6.93
	303823	13.08	8.66	46.25	6.93
	303699	13.08	8.66	46.25	33.91
3	303401	16.23	13.34	26.97	44.73
	303956	16.23	13.34	26.97	31.50
3.5	350781	38.25	11.10	71.12	14.30
4	303402	16.23	13.34	36.50	35.20
	303957	16.23	13.34	35.50	22.23
6.5	303313	16.23	13.34	58.72	12.98
	303392	25.65	13.34	58.72	25.40
	303449	38.23	13.34	58.72	37.31
	303687	25.65	13.34	69.85	11.13
10.5	303459	38.23	13.34	84.12	11.91
13.5	303307	25.65	16.51	71.42	12.70
	303448	38.23	16.51	71.42	24.61

a. Use these values to calculate radial distances for tubes in Delrin adapters


Using Tubes, Bottles, and Accessories

Introduction

*This chapter contains general instructions for filling and capping the labware used in Beckman Coulter preparative rotors, for selecting and using the appropriate accessories, and for recovering samples after a run. Individual rotor manuals provide specific instructions on tubes, bottles, and accessories that can be used in a particular rotor.**

Rotor use instructions are in [CHAPTER 4](#) for fixed-angle rotors, in [CHAPTER 5](#) for swinging-bucket rotors, and in [CHAPTER 6](#) for vertical-tube and near-vertical tube rotors. A table of chemical resistances is in IN-175. Reference information on some commonly used gradient materials is in [APPENDIX C](#).

Gradient Preparation†

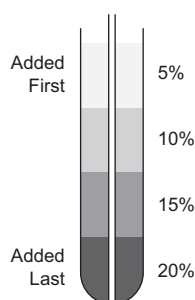

Many commercial gradient formers are available. These devices usually load a tube by allowing the gradient solutions to run down the side of the tube. The heaviest concentration is loaded first, followed by successively lighter concentrations. This method is acceptable for wettable tubes; however, loading a nonwetable tube (such as Ultra-Clear, polypropylene,[†] and polycarbonate) by allowing solutions to run down the side of the tube can cause mixing.


Gradients in nonwetable tubes can be prepared using a gradient former by placing a long syringe needle or tubing to the tube bottom and reversing the gradient chambers. In that way the lightest gradient concentration is loaded first, underlayed by increasingly heavier concentrations.

* A complete list of tubes, bottles, and adapters is provided in the latest edition of the Beckman Coulter Ultracentrifuge Rotors, Tubes & Accessories catalog (BR-8101), available at www.beckman.com.

† It has been reported, however, that polypropylene tubes have been made wettable by soaking them in a chromic acid bath for about 30 minutes (see Preparation of Polypropylene Centrifuge Tubes for Density Gradients, Anal. Biochem. 32:334-339. H. Wallace, 1969). Also, a method of making Ultra-Clear tubes wettable that has proven successful for some users is described at the end of this chapter.




You can also prepare preformed step gradients by hand, using a pipette. Carefully layer solutions of decreasing concentration by placing the tip of the pipette at the angle formed by the tube wall and the meniscus, or float the lighter gradient concentrations up by adding increased density solutions to the tube bottom using a hypodermic syringe with a long needle such as a pipetting needle.

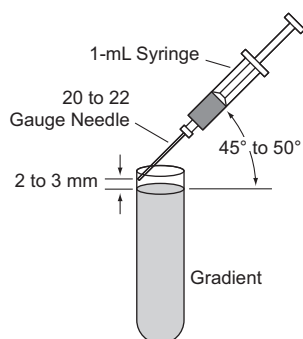

Another way to form a linear gradient is to allow a step gradient to diffuse to linearity. Depending on the concentration differential between steps and the cross-sectional area, allow 3 to 6 hours for diffusion at room temperature, and about 16 hours at 0 to 4°C. For diffusion of step gradient in Quick-Seal and capped straightwall tubes, slowly lay the tube on its side (tube contents will not spill, but make sure the tube does not roll). After two hours at room temperature, slowly set the tube upright.

Once the gradient is prepared, layer the sample on top of the gradient.

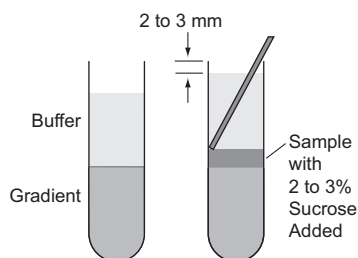

For thinwall tubes only partially filled with gradient, add a buffer solution to fill the tube to provide tube wall support. Although the gradient volume is reduced, sample volume is not changed.

NOTE If a partially filled thickwall tube is centrifuged, the tube does not require liquid support, and therefore, the buffer solution is not required.

Cesium Chloride Gradients

Cesium chloride gradients can be made by filling the tube with a homogeneous solution of CsCl and sample. Select a homogeneous CsCl solution density so that when it is distributed, its density range will encompass the density of the particle(s) of interest. Refer to [APPENDIX B](#) for an explanation of the use of the CsCl curves.

General Filling and Sealing or Capping Requirements

See [Table 3.1](#) for general filling and sealing or capping requirements for tubes and bottles used in preparative rotors. Maximum fill volume includes sample and gradient. Refer to individual rotor manuals for specific filling and capping requirements.

Table 3.1 Filling and Capping Requirements for Tubes and Bottles

Tubes or Bottles	Filling Level Requirements		
	Swinging-Bucket Rotors	Fixed-Angle Rotors	Vertical and Near-Vertical Tube Rotors
Polypropylene			
thinwall tubes	within 2–3 mm of top	full with cap	—
thickwall tubes	at least 1/2 full	1/2 full to max capless level or full with cap (Table 3.3)	—
OptiSeal tubes	full or plugged	full and plugged	full and plugged
Quick-Seal tubes	full and heat sealed	full and heat sealed	full and heat sealed
konical Quick-Seal tubes	full and heat sealed	—	—
konical open-top tubes	within 2–3 mm of top	—	—
bottles	—	min to max with screw-on cap or cap assembly (Table 3.3)	—
Ultra-Clear			
open-top tubes	within 2–3 mm of top	full with cap	—
Quick-Seal tubes	—	full and heat sealed	full and heat sealed
Polycarbonate			
thickwall tubes	at least 1/2 full	1/2 full to max capless level or full with cap or cap assembly (Table 3.3)	—
thickwall bottles	—	min to max with screw-on cap or cap assembly (Table 3.3)	—
Stainless Steel			
tubes	any level	any level with cap or cap assembly (Table 3.3)	—
Cellulose Propionate			
tubes	full	1/2 to max capless level; no cap	—

Table 3.1 Filling and Capping Requirements for Tubes and Bottles (*Continued*)

Tubes or Bottles	Filling Level Requirements		
	Swinging-Bucket Rotors	Fixed-Angle Rotors	Vertical and Near-Vertical Tube Rotors
Polyethylene tubes	at least $\frac{1}{2}$ full	$\frac{1}{2}$ to max capless level or full with cap	—
Corex/Pyrex tubes and bottles	at least $\frac{1}{2}$ full	$\frac{1}{2}$ to max capless	—

 WARNING

Handle body fluids with care because they can transmit disease. No known test offers complete assurance that they are free of micro-organisms. Some of the most virulent —Hepatitis (B and C) and HIV (I–V) viruses, atypical mycobacteria, and certain systemic fungi—further emphasize the need for aerosol protection. Handle other infectious samples according to good laboratory procedures and methods to prevent spread of disease. Because spills may generate aerosols, observe proper safety precautions for aerosol containment. Do not run toxic, pathogenic, or radioactive materials in these rotors without taking appropriate safety precautions. Biosafe containment should be used when Risk Group II materials (as identified in the World Health Organization *Laboratory Biosafety Manual*) are handled; materials of a higher group require more than one level of protection.

Filling and Plugging OptiSeal Tubes

OptiSeal tubes are not sealed prior to centrifugation; a Noryl plug, furnished with each tube, is inserted into the stem of filled tubes. When the tubes are loaded into the rotor with tube spacers (and rotor plugs, in vertical-tube and near-vertical tube rotors) in place, the g-force during centrifugation ensures a tight, reliable seal that protects your samples. For a detailed discussion on the use of OptiSeal tubes, refer to *Using OptiSeal Tubes* (publication IN-189).

Filling the Tubes

For filling convenience, use the appropriate eight-tube rack listed in [Table 3.2](#).

- 1 Use a pipette or syringe to fill each tube, leaving no fluid in the stem (see [Figure 3.1](#)).
 - Overfilling the tube can cause overflow when the plug is inserted; however, too much air can cause the tube to deform and disrupt gradients and sample bands, as well as increasing the force required to remove the tube from the cavity after centrifugation.

NOTE If air bubbles occur in the tube shoulder area, tilt and rotate the tube before it is completely filled to wet the tube.

- a. Homogeneous solutions of gradients and sample may be loaded into the tubes and centrifuged immediately. See [Gradient Preparation](#) above.
 - b. If the sample is to be layered on top, be sure to allow enough room for the sample so that there is *no fluid in the tube stem*.

- 2 After filling the tube, make sure that there is no fluid in the stem.



- a. (Draw off excess fluid with a syringe or pipette.
 - If necessary, wipe the inside of the stem with a lintless tissue.)

- 3 Fill the remaining tubes in the same manner.

Table 3.2 OptiSeal Tubes and Accessories^a


Size (mm)	Volume (mL)	Part Number ^b (pkg/56)	Spacer	Rack Assembly	Rotor
13 × 33	3.3	361627 	361678 (pkg/2) amber Ultem ^c 	361650	SW 55 Ti
13 × 48	4.7	361621 Bell-top 	361676 (pkg/2) amber Ultem 	361638	Type 50.4 Ti

Table 3.2 OptiSeal Tubes and Accessories^a (Continued)

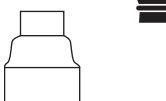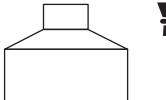Size (mm)	Volume (mL)	Part Number ^b (pkg/56)	Spacer	Rack Assembly	Rotor
13 × 51	4.9	362185 	362198 gold aluminum 	361638	VTi 90, VTi 65.2, NVT 90, NVT 65.2
			362199 black Noryl 	361638	VTi 65
16 × 60	8.9	361623 Bell-top 	361670 (pkg/2) amber Ultem 	361642	Type 90 Ti, Type 70.1 Ti
16 × 70	11.2	362181 	362202 gold aluminum 	360538	NVT 65, VTi 65.1
25 × 77	29.9	361625 Bell-top 	361669 (pkg/2) amber Ultem 	361646	Type 70 Ti, Type 50.2 Ti
			392883 (pkg/2) amber Ultem 	361646	SW 32 Ti, SW 28
25 × 89	36.2	362183 	362204 gold aluminum 	360542	VTi 50 VTi 50.1

- a. Spacers are shown in the correct orientation for placement onto tubes.
b. Disposable plastic plugs included.
c. Ultem is a registered trademark of SABIC GLOBAL TECHNOLOGIES B.V.


Figure 3.1 Filling OptiSeal Tubes*

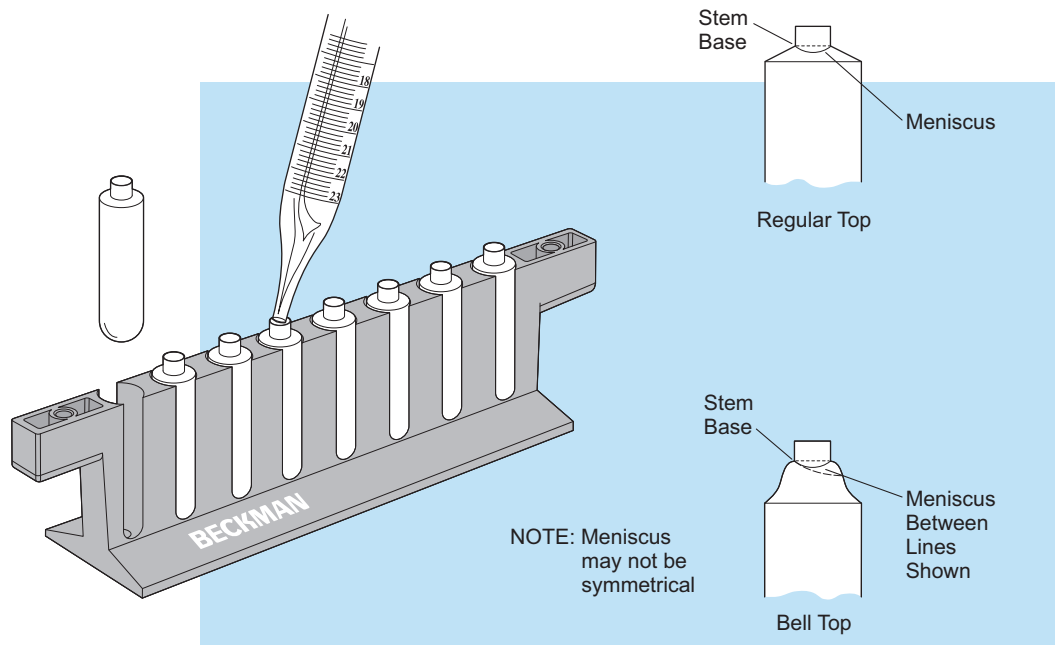

Seating the Tube Plugs

Eight tubes can be prepared for use at once in the specially designed racks listed in [Table 3.2](#).


NOTE The Ultem spacers (361678) snap onto the 3.3-mL tubes (361627). To avoid disturbing the sample or splashing out liquid, put the spacers on these tubes before inserting the plugs.

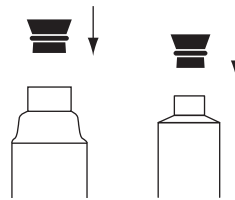

- 1 Make sure that no fluid is in the tube stem and that the stem is clean and dry.
- 2 Insert a Noryl plug assembly (plug and O-ring—shipped assembled) in each tube stem.
- 3 Set the plug seating bar on the rack, ensuring that the pegs at each end fit into the rack openings.

* Stems are large enough to accept standard pipettes.

- 4 Press firmly straight down all along the top of the bar.
 - When you remove the bar, the plugs should be straight and seated into the stems.

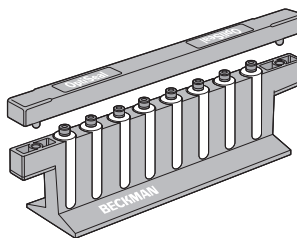

- 5 Check the tubes to be sure all plugs are seated. If any plugs are not seated, seat them individually.



Filling the Tubes

Fill each tube to the base of the stem, using a syringe with a 14-gauge or smaller needle. Do not leave a large air space — too much air can cause excessive tube deformation, disrupting gradients or sample.



Homogeneous solutions of gradients and sample may be sealed into the tubes and centrifuged immediately. Step gradients should be loaded with a long needle inserted to the bottom of the tube. Because the tubes are not wettable, load the light end of the preformed gradient first and float up with more dense solution. If the sample is to be layered on top, be sure to allow enough room for the sample so that the tube stem is not filled.

Form a continuous gradient quickly from a step gradient as follows.

- 1 Load a step gradient into the tube.
- 2 Before layering on the sample, slowly tip the tube sideways and lay it on its side (tube contents will not spill, but be sure that the tube doesn't roll). The increased surface area between gradient steps allows a linear gradient to form by diffusion more quickly.


- 3 After 2 hours at room temperature, slowly set the tube upright.
 - a. Alternately, the tube containing the step gradient can stand upright in a refrigerator overnight and a continuous gradient will be formed.)
- 4 Layer the sample onto the gradient as described below.

Layering Sample

Layer sample onto the tube contents as follows.*

- 1 Slip the 5-mL polyethylene funnel (342415) down over the tube stem.
- 2 Compress the tube until the gradient is forced up into the funnel (Figure 3.2). (A sample application block [342694] is available for holding and compressing the tubes.) Force only enough gradient into the funnel to provide a surface for sample application.

Figure 3.2 Layering Sample

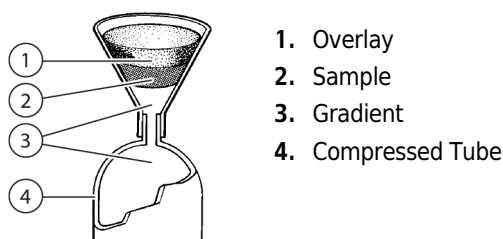

- 3 Layer the sample onto the meniscus of the gradient in the funnel. If an overlay is being used, carefully layer it on top of the sample. (Do not use an oil overlay in Ultra-Clear tubes.)
- 4 Gently release the pressure compressing the tube and allow the sample and overlay to drain back into the tube.
- 5 After filling the tube, remove the funnel and wipe off the exterior of the filled tube; make sure that no fluid is trapped in the stem. (The tube stem should be clean and dry before sealing.)

* Method and apparatus for layering sample onto tube contents patented: U.S. Pat. No. 4,167,955; Swiss Pat. No. 631089; British Pat. No. 2,022,455; Italian Pat. No. 1,121,265

- 6 Fill the remaining tubes in the same manner and place the tubes in the appropriate tube rack. The tube racks are color-coded to indicate the tube diameter that fits in each.

Sealing the Tubes

The following procedures provide the two methods for heat-sealing Quick-Seal tubes using the Tube Topper.

NOTE Separate instructions are provided for the North American and EU/UK markets.

WARNING

Risk of personal injury. Touching the heated tip of the Tube Topper will cause burns. When the push button is pressed the tip heats almost immediately. To avoid burns:


NORTH AMERICAN MARKET — When not in use, be sure the push button is turned to **LOCK** position unless you are actually sealing a tube.

EU/UK MARKET — When not in use, be sure the circular safety switch is turned to **OFF (O)** unless you are actually sealing a tube.

- 1 Remove the Tube Topper from the charging stand.
 - a. **NORTH AMERICAN MARKET** — Make sure the push button is turned to the **OFF (LOCK)** position. Insert the ends of the Tube Topper tip into the two openings of the copper strips at the end of the Tube Topper device (Figure 3.3).


EU/UK MARKET — Make sure the circular switch is turned to the **OFF (O)** position. Insert the ends of the Tube Topper tip through the plastic sleeve and into the two openings at the end of the Tube Topper device. Slide the plastic sleeve over the end of the Tube Topper. Tighten the screws to secure the probes in place (Figure 3.4).

Figure 3.3 North American Markets

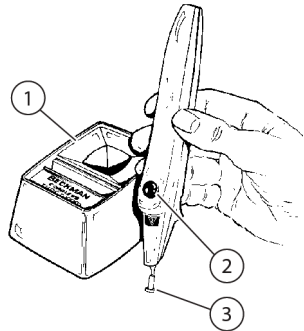

1. Charging Stand
2. Push button
3. Tip


Figure 3.4 EU/UK Markets

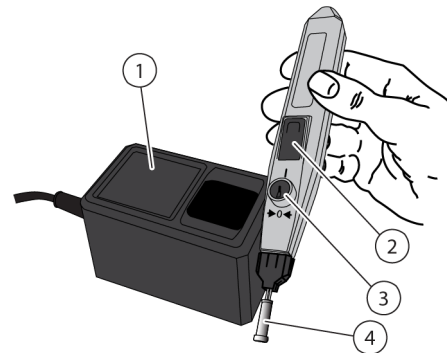

1. Charging Stand
2. Push button
3. Circular Safety Switch
4. Tip

- 2 Place a seal former on each tube stem (Figure 3.5). (The Teflon* coating on the seal formers is permanent. Do not scratch the interior of the formers, as you may damage this coating.)

Figure 3.5 Seal Former

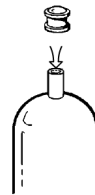

- 3 Seal each tube using [METHOD A — WITH THE SEAL GUIDE](#) or [METHOD B — WITHOUT THE SEAL GUIDE](#), below. Method A is preferable when sealing smaller tubes or when resealing a tube that leaks.

CAUTION

Always keep the Tube Topper in its charging stand when not in use. Do not lay the unit against any surface after use until the tip has cooled (approximately 5 minutes).

* Registered trademark of E.I. du Pont de Nemours & Company

METHOD A — WITH THE SEAL GUIDE

- a. Place a seal guide (with the flat side down) over the seal former (Figure 3.6).

Figure 3.6 Seal Guide



- b. **NORTH AMERICAN MARKET** — Turn the Tube Topper push button to the **ON (USE)** position. Press the push button and wait 3 to 5 seconds for the tip to heat.
- EU/UK MARKET** — Turn the Tube Topper circular safety switch to the **ON (I)** position. Press the **ON** button and wait approximately 9 seconds for the tip to heat.
- c. Apply the tip of the Tube Topper vertically to the seal former (Figure 3.7 and Figure 3.8). Press down gently for about 10 seconds. The seal guide should move down the tube stem until it rests on the tube shoulder. Using the seal guide prevents the seal former from being pressed into the tube shoulder. (For tube 344625 only, the seal former should move only to about 2 mm above the tube shoulder.)

Figure 3.7 North American Market


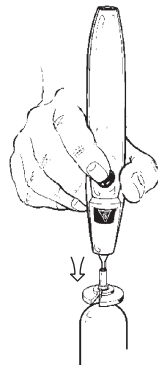


Figure 3.8 EU/UK Market

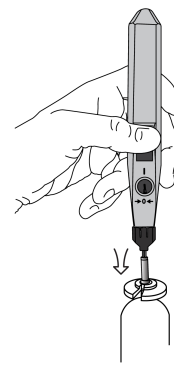


NOTE Always apply the tip of the Tube Topper vertically to the seal former. Apply gentle pressure when sealing the tube.

- d. When the seal guide has moved to the correct position, remove the Tube Topper and pinch the circular seal guide to hold the seal former in place.


- e. Place the heat sink (small end) over the cap for 2 to 3 seconds while the plastic cools — do NOT let the seal former pop up (Figure 3.9). (If the seal former does pop up, the tube may not have an adequate seal and may need to be resealed.)

Figure 3.9 Heat Sink

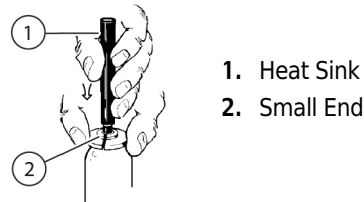

- f. Remove the heat sink and seal guide. When the seal former cools, remove it by hand or with the removal tool (361668) (Figure 3.10). Save the seal guide and former for future use.

Figure 3.10 Remove
Tool

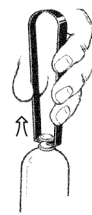

METHOD B — WITHOUT THE SEAL GUIDE

- a. **NORTH AMERICAN MARKET** — Turn the Tube Topper push button to the **ON (USE)** position. Press the push button and wait 3 to 5 seconds for the tip to heat.

EU/UK MARKET — Turn the Tube Topper circular safety switch to the **ON (I)** position. Press the **ON** button and wait 9 seconds for the tip to heat.


NOTE Always apply the tip of the Tube Topper vertically to the seal former. Apply gentle pressure when sealing the tube.

- b. Apply the tip of the Tube Topper vertically to the seal former for about 10 seconds (Figure 3.11 and Figure 3.12). The seal former should move down the tube stem until it just rests on the tube shoulder. Be careful NOT to press the seal former into the tube shoulder — it may cause the tube to leak. (For tube 344625 only, the seal former should move only to about 2 mm above the tube shoulder.)

Figure 3.11 North American Market




Figure 3.12 EU/UK Market

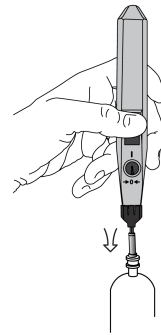

NOTE It is very important to apply the heat sink immediately. To do so, we recommend that you have it in one hand, ready to apply as soon as needed.

- c. Remove the Tube Topper. IMMEDIATELY place the large end of the heat sink over the seal former (Figure 3.13). Hold it there for a few seconds while the plastic cools — do NOT let the seal former pop up. (If the seal former does pop up, the tube may not have an adequate seal and may need to be resealed.)


Figure 3.13 Heat Sink

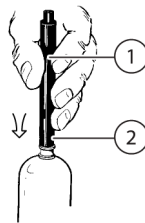

Immediately
1. Heat Sink
2. Large End


- d. Remove the heat sink. When the seal former cools, remove it by hand or with the removal tool (361668)(Figure 3.14).

Figure 3.14 Removal Tool

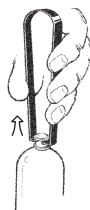

-
- 4 After completing either heat-sealing method, squeeze the tube gently (if the tube contents may be disturbed) to test the seal for leaks (Figure 3.15). If the tube leaks, try resealing it using [METHOD A — WITH THE SEAL GUIDE](#).

Figure 3.15 Test for Seal
Leaks

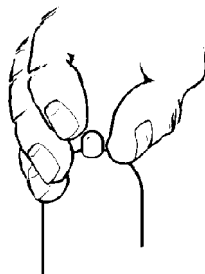

-
- 5 The tube is now ready for centrifugation. Seal the remaining tubes.
-
- 6 Return the Tube Topper to its charging stand when finished.
-


Handy Tips

- Always keep the Tube Topper in its charging stand when not in use. Do not lay the unit against any surface after use until the tip has cooled (approximately 5 minutes).
- **NORTH AMERICAN MARKET** — Always return the push button to **OFF (LOCK)** position when not in use.
EU/UK MARKET — Always return the circular safety switch to the **OFF (O)** position when not in use.
- Always apply the tip of the Tube Topper vertically to the seal former. Apply gentle pressure when sealing the tube.
- Never allow the seal former to pop up during the sealing process. If the seal former does pop up, the tube may not have an adequate seal and may need to be resealed.
- New users may find the sealed stem slightly slanted or elongated after removing the seal former. Reseal the tube only if the tube spacer won't fit on top of the tube.
- You may get an overflow of melted plastic around the seal former during the sealing procedure. Check that the tube spacer fits onto the sealed tube.
- The Teflon coating on the seal formers is permanent. Do not scratch the interior of the formers, as you may damage this coating.

Loading the Rotor

To support the top of the tube during centrifugation, each Quick-Seal tube is used with a spacer, a floating spacer, or a combination of both (refer to the applicable rotor manual for the correct combination). Insert the sealed tubes into the rotor and install the correct spacers and/or floating spacers as described in the rotor manual. The particular combination depends on the type of rotor you are using. In fixed angle and swinging bucket rotors it is particularly important that the top of each tube be supported. In a near vertical tube or vertical tube rotor, it is especially important that the entire tube cavity be filled. Therefore, in near vertical tube or vertical tube rotors, metal plugs are inserted over the spacers and then screwed in to seal the tube cavities ([Figure 3.16](#)).

Figure 3.16 Spacers

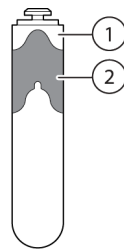


1. Spacer
2. g-Max Floating Spacer

Filling Open-Top Tubes

Open-Top Polypropylene Tubes

Open-top polypropylene tubes are used in swinging-bucket and fixed-angle rotors.


Swinging-Bucket Rotors

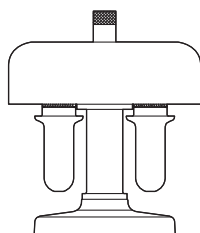

Fill all opposing tubes to the same level.

- Thinwall Tubes—Fill to within 2 or 3 mm of the top for proper tube wall support.
- Thickwall Tubes—Fill at least half full.

Fixed-Angle Rotors

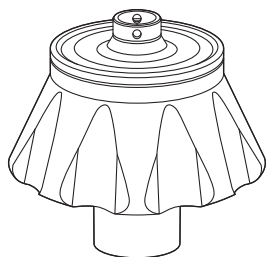

Fill all opposing tubes to the same level

- Thinwall Tubes — Must be completely filled; liquid and cap for support of the tube wall is critical.
- Thickwall Tubes—Can be partially filled and centrifuged as indicated in the applicable rotor manual. Speed reductions may be required for these partially filled tubes. For greater fill volumes and faster speeds, tube caps should be used. Refer to the applicable rotor manual for fill volumes and speed limitations.

Other Open-Top Tubes

Open-top tubes of other materials can also be used in fixed angle and swinging-bucket rotors. (Vertical-tube and near-vertical tube rotors use only OptiSeal or Quick-Seal tubes.) Fill these tubes as indicated below.

Polycarbonate

Thickwall polycarbonate tubes can be centrifuged partially filled. Observe maximum rotor speeds and fill volumes listed in the applicable rotor manual.

UltraClear

Fill all opposing tubes to the same level.

- For *swinging-bucket* rotors, fill to within 2 or 3 mm of the top of the tube.

- Fill thickwall polypropylene tubes at least half full to maximum level in *fixed-angle* rotors. Speed reduction is required. Refer to the applicable rotor manual.

Polypropylene

Fill all opposing tubes to the same level.

- For *swinging-bucket* rotors, fill to within 2 or 3 mm of the top of the tube.
- Fill thickwall polypropylene tubes at least half full to maximum level in *fixed-angle* rotors. Speed reduction is required. Refer to the applicable rotor manual.

Polyethylene

For *swinging-bucket* and *fixed-angle* rotors, fill these tubes from half full to maximum level. Refer to the applicable rotor manual.

Stainless Steel


Because of their strength, stainless steel tubes can be centrifuged while filled to any level (with all opposing tubes filled to the same level). However, run speeds must be reduced due to their weight. The criteria for speed reduction depends on the tube-cap material and the strength of the rotor being used. Refer to the applicable rotor manual or *Run Speeds for Stainless Steel Tubes* (publication L5-TB-072) for correct run speeds.

Capping Tubes

Caps must be used with thinwall polypropylene and Ultra-Clear tubes in fixed-angle rotors. To prevent spillage, thickwall polypropylene, polycarbonate, and stainless steel tubes must be capped when fill levels exceed the maximum level for uncapped tubes as listed in the applicable rotor manual.

Cap requirements depend on the tube or bottle material, diameter, and wall thickness, as well as on the rotor. The applicable rotor instruction manual specifies which cap should be used with a particular tube or bottle.

Tube Cap Assemblies

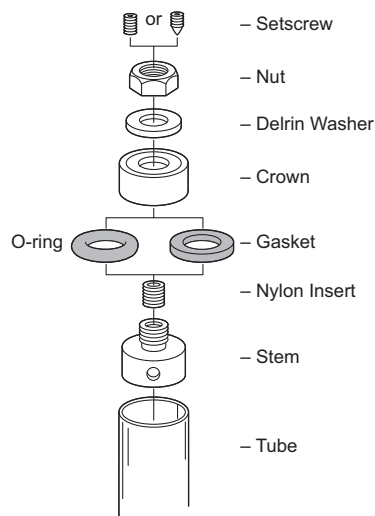

A tube-cap assembly includes

- a stem,
- a nylon insert,
- an O-ring or flat gasket,
- a crown,
- a Delrin crown washer (in red, blue, and black aluminum caps),
- a hex-shaped nut, and
- a stainless steel setscrew.

The stem supports the upper portion of the tube. To provide tube support during centrifugation, the stem is longer for thinwall tubes than for thickwall or metal tubes. Some stems have an abraded surface to increase friction between the O-ring and the stem, minimizing rotation of the stem when the cap nut is tightened. The O-ring or gasket seals the cap-to-tube interface.

The crown seats on the rotor tube cavity counterbore and supports the stem and the nut during centrifugation. In some high-performance rotors, tube caps have crown washers. The washer minimizes friction, which would reduce the effective tightening of the cap nut, and also protects the nut and the crown. After the tube has been capped, tightened, and filled, the setscrew is used to seal the filling hole in the stem by seating against the nylon insert.

Refer to [Table 3.3](#) for detailed information about tube caps.

CAUTION

Do not interchange tube caps or tube-cap components, even if they appear to be the same. Tube caps are designed specifically for a particular tube in a particular rotor. Cap stems and crowns are often machined differently for each type of rotor to ensure proper sealing and support and to withstand stresses experienced during centrifugation. The uneven weight difference between an O-ring cap and a comparable flat-gasket cap (as much as 0.7 gram) could damage the rotor. Store tube caps assembled, dry, and classified according to the tube and rotor for which they are designed.

Titanium Caps

High-strength titanium cap assemblies for thinwall Ultra-Clear and polypropylene tubes are required for maximum rotor speeds in the Type 90 Ti and 70.1 Ti rotors. Titanium caps can be identified by the darker gray, shiny metal. The cap crown is specially machined to lock onto the cap stem. To ensure proper compression of the O-ring, these caps must be tightened with a torque wrench while the capped tube is held in the tube-cap vise.

A special crimp-lock cap assembly is required to provide the reliable seal necessary for maximum rotor speed in the Type 70.1 Ti rotor. The 25 × 83-mm thinwall polypropylene tube is crimped between the titanium crown and the aluminum stem. Instructions for assembling the tube and cap are in the Type 70 Ti rotor instruction manual. A special tool kit (338841) is required.

Table 3.3 Tube Cap Assemblies for Open-Top Tubes in Fixed-Angle Rotors^a

Tube Cap Assembly ^b	Hex Nut	Crown	Set-screw	Insert	O-ring or Gasket	Stem	Tube Type	Rotor Type
8 mm (5/16 in.)								
303624	303379	303809	—	—	303730	303377	UC ^c	90 Ti, 70.1 Ti
303658	303379	303810	—	—	303370	303377	UC	
13 mm (1/2 in.)								
303113	301870	307004	—	—	344672	307005	SS	70 Ti, 50.4 Ti, 50.2 Ti, 45 Ti
305022 ^d	301870	307004	—	—	344672	302331	SS	70.1 Ti, 70 Ti, 50.4 Ti, 50.2 Ti, 45 Ti
346256	301870	307004	803543	302312	344672	346246	thinwall PP, UC, SS	90 Ti, 70.1 Ti, 70 Ti, 50.4 Ti, 50.2 Ti, 45 Ti
16 mm (5/8 in.)								
303319	301870	307006	338864	302312	301869	302266	SS	90 Ti, 70 Ti, 50.2 Ti, 45 Ti
330860	301870	330774	803543	302312	858046	330788	thinwall PP, UC	70 Ti, 50.2 Ti, 45 Ti
338907 ^e	301870	338911	338864	302312	878572	338910	thickwall PP, PC	90 Ti, 70 Ti, 50.2 Ti, 45 Ti
341968 ^f	335320	335319	338864	302312	858046 870380	341969 341969	thinwall PP UC	90 Ti, 70.1 Ti 90 Ti, 70.1 Ti
25 mm (1 in.)								
302133	301870	302169	338864	302312	301473	302138	SS	70 Ti, 50.2 Ti
331151	330791	331153 ^h	338864	302312	334280	331152	thinwall PP, UC	70 Ti, 50.2 Ti
338906 ^e	330791	338915 ^h	338864	302312	878188	338908	thickwall PP, PC	70 Ti, 50.2 Ti
337927 ^g	330791	338863 ^h	338864	302312	—	338865	thinwall PP	70 Ti

Table 3.3 Tube Cap Assemblies for Open-Top Tubes in Fixed-Angle Rotors^a (Continued)

Tube Cap Assembly ^b	Hex Nut	Crown	Set-screw	Insert	O-ring or Gasket	Stem	Tube Type	Rotor Type
38 mm (1 1/2 in.)								
326905	301870	326890	338864	302312	801761	326899	SS	45 Ti
330901	330791	330793 ^h	338864	302312	346242	330900	thinwall PP, UC	45 Ti
338905 ^e	330791	338913 ^g	338864	302312	341767	338909	thickwall PP, PC	45 Ti

a. Tube caps are not available

b. Tube caps are aluminum unless otherwise noted.

c. Abbreviations: PP = polypropylene; PC = polycarbonate; SS = stainless steel; UC = UltraClear

d. Aluminum and stainless steel

e. Tube cap is optional. Use a tube cap when centrifuging a thickwall tube at its maximum fill capacity.

f. Titanium

g. Aluminum and titanium

h. Washer, part number 330899, is also required.

Aluminum Caps

Aluminum caps are anodized for corrosion resistance, with colored crowns for identification.




Red-anodized. Aluminum caps (aluminum stem and crown) with red-anodized crowns are used with thinwall Ultra-Clear and polypropylene tubes in high-performance rotors. These extra-strength caps are designed for the greater forces generated in the high-performance rotors. The cap nut should be tightened with a torque wrench while the tube is held in the tube-cap vise.

Blue-anodized. Aluminum caps with blue-anodized crowns are used with thickwall polypropylene and polycarbonate tubes for centrifugation at their maximum fill volumes in high-performance rotors. The cap nut should be tightened with a torque wrench while the tube is held in the tube-cap vise.

Clear- and black-anodized

- Clear-anodized crown aluminum caps that use O-rings for sealing are used in many rotors with many types of tubes. Refer to [Table 3.3](#). The caps should be hand tightened with a hex driver while the tube is held in the tube-cap vise (refer to [Assembling Tube Caps](#), below).
- Aluminum caps that use flat gaskets for sealing are used with small-diameter (13-mm) thinwall Ultra-Clear and polypropylene tubes in all fixed-angle rotors except Types 25 and 19. They are also used with stainless steel tubes. Some caps for very small-diameter (less than 13-mm) tubes do not have filling holes (nylon insert or setscrew). The tube crown is made from a lighter-weight aluminum alloy than that used for other clear aluminum caps; therefore, *do not interchange cap parts or use these caps in place of O-ring caps, since the weight difference can cause rotor imbalance*. The caps should be hand tightened with a hex driver while the tube is held in the tube-cap vise.

Inspecting and Lubricating Tube Caps

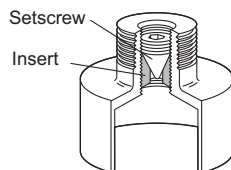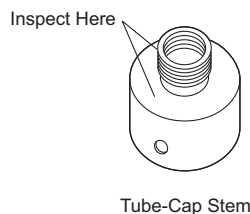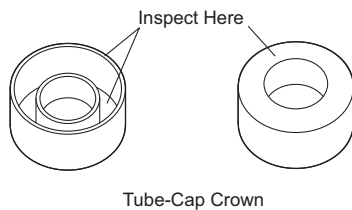

1. Inspect cap components before each use. Replace any damaged components.
 - Inspect the cap crown for stress cracking, and check the stem and nut threads for damage or signs of wear and for adequate lubrication.
 - Inspect the O-ring or gasket for cracks, nicks, or flattened areas.
 - Inspect the underside of the stem; the white nylon insert should not protrude below the filling hole.
 - If the cap assembly has a filling hole, run the setscrew in against the nylon insert, making sure the setscrew will not displace the insert.
 - Check the setscrew hex socket for damage that would prevent tightening or removal.
2. Regularly apply a thin, uniform coat of Spinkote lubricant (306812) on the stem threads.

NOTE Keep the O-ring or flat gasket dry and free from lubricant during assembly. Wet or greased O-rings or gaskets will slip when the cap nut is tightened and the cap will not seal properly.

Assembling Tube Caps

See [Figure 3.3](#), [Figure 3.4](#), and [Table 3.4](#) for required tools and torque requirements.

CAUTION

Do not use damaged wrenches or hex drivers, or tools that have burrs. A burred tool can score the crown, which could then fail and damage the rotor.

Figure 3.17 Tools Used to Assemble Tube Caps

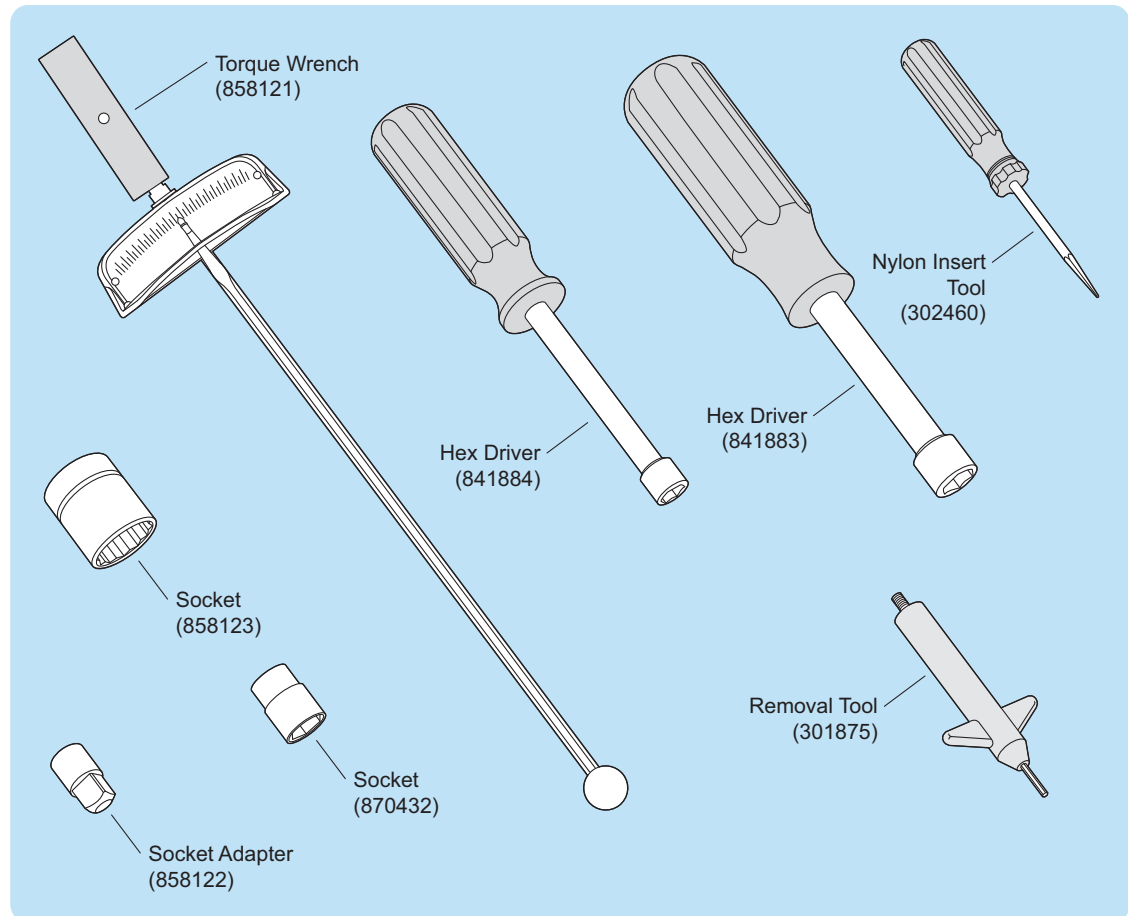


Table 3.4 Required Tools and Torque Values


Tightening Tool	Tube Caps ^a	Cap Nut ^b Size/ Part Number	Torque Value
Torque wrench (858121) Socket (870432)	titanium cap, 341968	11 mm (7/16 in.) 335320 (titanium)	10.2 to 11.3 N•m (90 to 100 in.-lb)
Torque wrench (858121) Socket (858122) Socket (858123)	331151 (red) 330901 (red) 338905 (blue) 338906 (blue)	20 mm (3/4 in.) 301870	11.3 to 13.6 N•m (100 to 120 in.-lb) for the first four runs; 11.3 N•m (100 in.-lb) starting with the fifth run
Hex driver (841884)	303624 303658	8 mm (5/16 in.) 303379	hand tighten
Hex driver (841883)	303113, 346256, 305022, 330860, 338907, 303319, 326891, 302133, 326905, 337927	11 mm (7/16 in.) 301870	hand tighten

a. Unless otherwise indicated, caps are clear-anodized aluminum.

b. Unless otherwise indicated, cap nuts are aluminum.

1 If possible, fill tubes one-half to three-quarters full before capping.

- Small-diameter tubes that use caps without filling holes (caps 303624, 303658, 303113, and 305022) must be completely filled before capping.

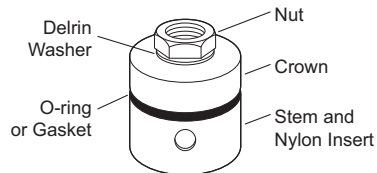

2 Loosely assemble the stem, the O-ring or gasket, the crown, the crown washer (if applicable), and the nut.


- The nylon insert should already be installed in the stem.*

a. For titanium caps, turn the crown slightly to be sure it is properly seated on the stem.

* Nylon inserts are installed in the stems of cap assemblies purchased as a unit. Stems ordered separately do not contain an insert. See Section 7 for installation.

- 3** Slide the tube up around the stem PAST the O-ring or gasket as shown in [Figure 3.18](#), slightly rotating the cap assembly.
- The tube wall should pass between the O-ring or gasket and the crown so that the top of the tube rests on the underside of the crown.
- a.** Tighten the nut by hand just enough to hold the tube cap in place.

Figure 3.18 Tube Cap Installation*

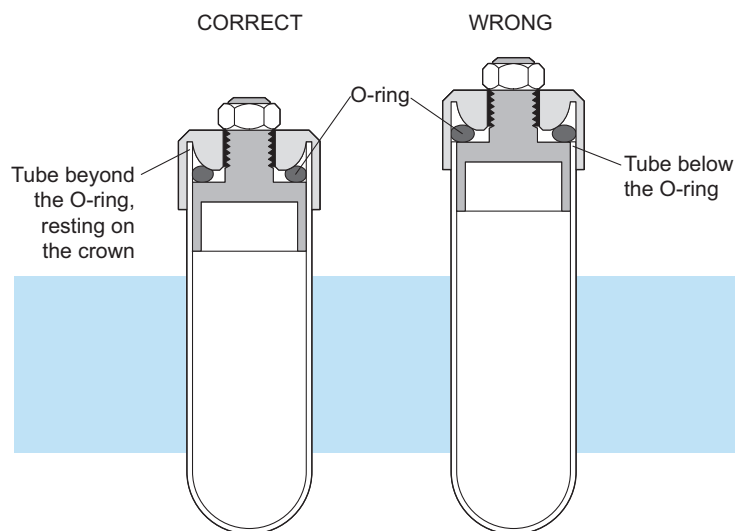


- 4** Position the capped tube in the appropriate-sized hole from the underside of the tube-cap vise (305075).
- The vise must be correctly mounted to the bench with the clamping positioned on the right (see [Figure 3.19](#)), or crimping of the crown may result.

NOTE While holding the tube with one hand, tighten the vise around the crown by using the clamping knob. Make sure that the cap and the tube are level (horizontal).

- 5** Tighten the cap nut as described in [Table 3.4](#).
- 6** Use a syringe to finish filling the tube through the filling hole in the stem.
- Thinwall tubes must be as full as possible to prevent tube collapse.
 - Thickwall tubes may be filled to within 13 mm of the top, but may still collapse if not completely full.
 - Stainless steel tubes may be filled to any level.
 - Tubes placed opposite each other in the rotor must be filled to the same level.

* The tube must be pushed up past the O-ring so that the crown will clamp the tube and NOT the O-ring.

Figure 3.19 Tube Cap Vise*

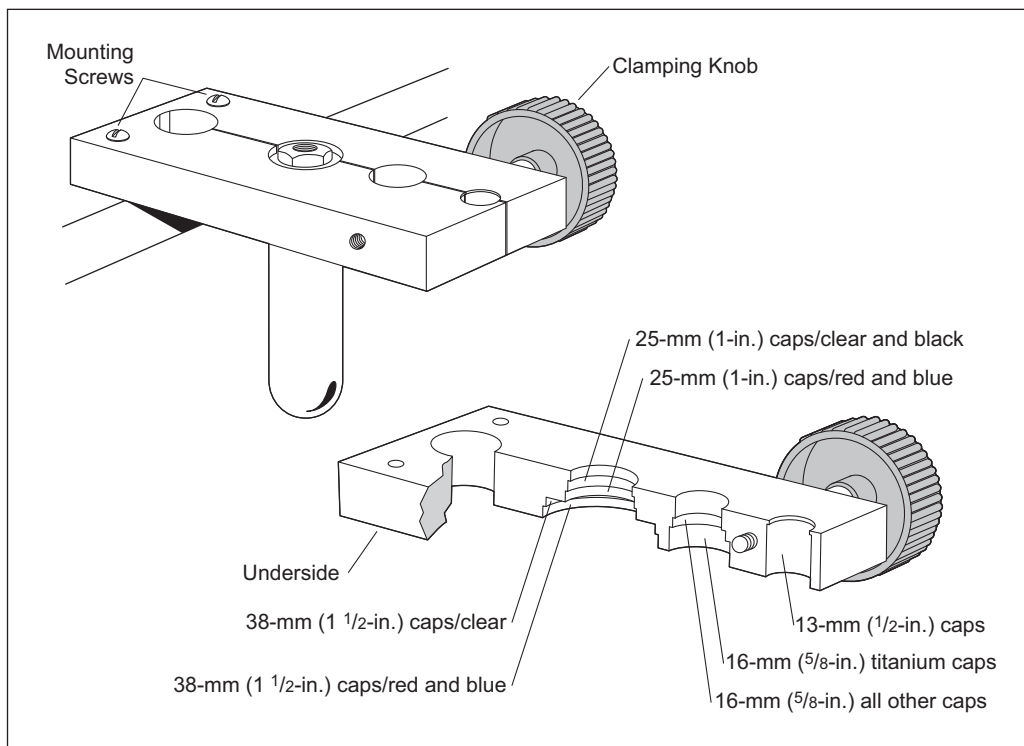

Filling and Capping Tubes

To prevent spillage and provide support, polycarbonate and polypropylene bottles used in fixed-angle rotors must be capped when fill levels exceed the maximum level allowed for uncapped bottles. Bottles should be filled to maximum fill levels when spun at maximum rated speeds. Unless specified otherwise, the minimum recommended volume for bottles is half full; this will require reduced rotor speed for optimum labware performance. Refer to [Table 3.5](#) and the applicable rotor manual for bottle fill levels and cap requirements.


* Screw the vise to a bench or table top for operation. The vise must be correctly mounted, with the clamping knob positioned on the right, or crimping of the crown may result.

Table 3.5 Available Bottles, Assembly and Operation^a

Bottle		Required Cap Assembly		Bottle and Cap Assembly	Volume (mL)		Rotor	Maximum speed (rpm) ^b	Required Adapter
Part Number	Dimensions (mm)	Material	Part Number		Max.	Min.			
355651	16 × 76	Noryl	355604	355603	10.4	5 ^b	Types 90 Ti, 70.1 Ti	65 000	
355654	25 × 89	aluminum	355619	355618	26.3	16 ^b	Types 70 Ti	60 000	
							Type 50.2 Ti	50 000	
355655	38 × 102	aluminum	355623	355622	70	35 ^b	Type 45 Ti	45 000	
355627	60 × 120	Delrin (w/Noryl plug)	362247	334025	250	250	Type 19	19 000	

a. Bottles are polycarbonate unless otherwise indicated.

b. Several rotors must be centrifuged at reduced speeds when bottles are filled below maximum fill volume: Types 90 Ti and 70.1 Ti at 60 000 rpm; Types 70 Ti and 50.2 Ti at 50 000 rpm; Type 45 Ti at 35 000 rpm.

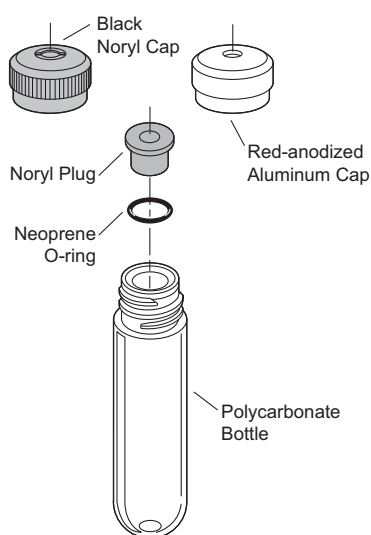

Cap bottles with three-piece cap assemblies as follows:

1. Be sure the O-ring, plug, and bottle lip are dry and free of lubrication.
2. Place the O-ring on the underside of the plug.
3. Insert the plug into the neck of the bottle, ensuring that no fluid contacts the O-ring.
4. Tighten the cap by hand.


Sample Recovery

CAUTION

If disassembly reveals evidence of leakage, you should assume that some fluid escaped the container or rotor. Apply appropriate decontamination procedures to the centrifuge, rotor, and accessories.

Sample recovery depends on the type of labware used, the component(s) isolated, and the analysis desired. The Beckman Coulter Universal Fraction Recovery System (343890) can be useful when recovering sample from tubes (see publication L5-TB-081).

Capped Tubes

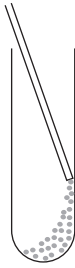


The usual methods of recovering supernatants or pellets include decanting or withdrawing the gradient and scraping pellets from the tube bottom.

- Remove tube caps carefully to avoid sample mixing.
- If tubes will be reused, scrape pellets out with a plastic or wooden tool; scratches on tube interiors caused by abrasive or sharply pointed tools can result in tube failure during subsequent runs.

OptiSeal Tubes

Centrifugation exerts high forces on plastic labware. The effect of these forces on OptiSeal labware is compression of the tube, characterized by tube deformation that, even if slight, causes a decrease in internal volume. OptiSeal labware is designed to contain the resulting slight pressure increase during separation, as well as during normal post-separation handling. However, a small volume ($\approx 50 \mu\text{L}$) of fluid may occasionally “ooze” from around the plug onto the tube stem area as a plug is removed. Therefore, we recommend using a tissue to contain escaped fluid when extracting plug assemblies from tubes.

- 1 After centrifugation, use the spacer removal tool (338765) or a hemostat to carefully remove the spacers, taking care not to scratch the rotor cavities.
 - (A tube will sometimes come out of the rotor cavity along with the spacer.)
 - a. Separate the tube from the spacer with a twisting motion.)

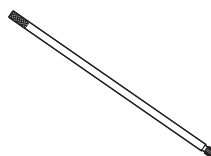


Spacer Removal Tool

NOTE SW 32 Ti and SW 28 rotors only—Use the spacer removal tool (338765) to remove the spacer and tube together from the rotor bucket. Place the tubes in the rack. Grasp the tube and use the spacer removal tool in a lifting and twisting motion to remove the spacer.

NOTE Centrifugation causes a slight vacuum to build up in the tube cavity, occasionally resulting in a suction effect when removing the tubes from the rotor. This effect is especially pronounced in a rotor that has been centrifuged at a low temperature. A brief delay (approximately 5 minutes) after the rotor comes to rest before removing the tubes will make tube removal easier. If you experience difficulties in removing the tubes from the rotor, use a gentle twisting or rocking motion, and remove the tube slowly to avoid sample mixing.

- 2 Remove the tube with the extraction tool (361668), grasping the base of the stem only—do NOT try to remove the tubes by pulling on the plugs.
 - Some tube deformation occurs during centrifugation, which causes a slight internal pressure to develop inside the tube.

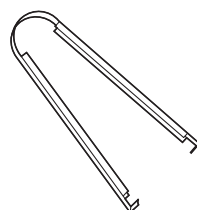

Extraction Tool
(361668)


3 Place the tubes back into the tube rack.

- Openings in the rack allow the tubes to be pierced either from the bottom or sides, permitting fractions to be easily collected regardless of the type of separation.


NOTE If you want to collect particles from the tube side or bottom, first create an air passage by removing the tube plug (see instructions below) or inserting a hollow hypodermic needle in the top of the tube.

4 Use one of the following methods to retrieve the sample:


- a.** Puncture the side of the tube just below the sample band with a needle and syringe and draw the sample off.
- Take care when piercing the tube to avoid pushing the needle out the opposite side.

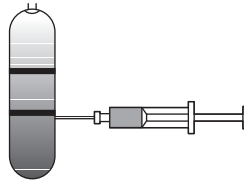

- b.** Puncture the bottom of the tube and collect the drops

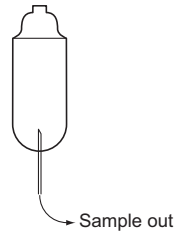


- c.** Aspirate the sample from the tube top by removing the tube plug (see instructions below), then aspirating the sample with a Pasteur pipette or needle and syringe.

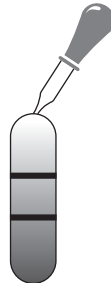

- d.** Slice the tube, using the Beckman CentriTube Slicer (303811).
- Refer to publication L-TB-010 for instructions for using the CentriTube Slicer.

- 1) Use CentriTube Slicer (347960) and CentriTube Slicer Adapter (354526) for 13-mm tubes.
 - (Tubes are pressurized after centrifugation, so pierce the tube top with a needle to relieve pressure before slicing.)

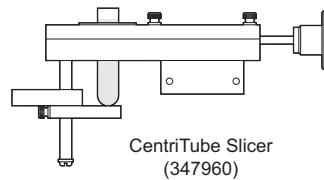


Removing Plugs from Tubes

- 1 Place the tube rack insert over the tubes in the rack.
- 2 Press down on the rack insert on each side of the tube being unplugged to hold the tube in place during plug removal.






NOTE Do not hold onto or squeeze the tubes. Tube contents will splash out when the plug is removed if pressure is applied to the tube.

- 3 While pressing down on the rack insert, use the extraction tool to firmly grasp the plug.

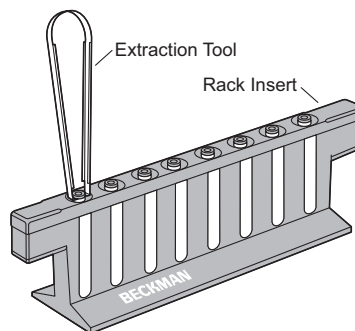

- 4 Use a slight twisting motion to slowly release any residual internal pressure when pulling the plug assembly from the tube.
 - 5 Repeat for each tube.
-

Quick-Seal Tubes

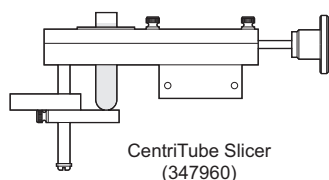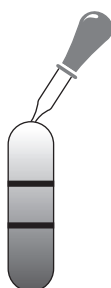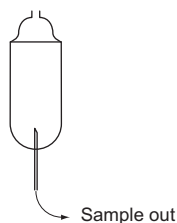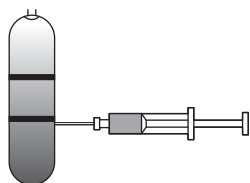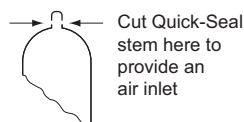

There are several methods of recovering fractions from Quick-Seal tubes. One of the following procedures may be used.

NOTE If you plan to collect particles from the tube side or bottom, first create an air passage by snipping the stem or inserting a hollow hypodermic needle in the top of the tube.

- Puncture the side of the tube just below the band with a needle and syringe and draw the sample off. Take care when piercing the tube to avoid pushing the needle out the opposite side.
- Puncture the bottom of the tube and collect the drops.
- Aspirate the sample from the tube top by snipping off the tube stem, then aspirating the sample with a Pasteur pipette or needle and syringe.
- Slice the tube, using the Beckman CentriTube Slicer (347960). Refer to publication L-TB-010 for instructions for using the CentriTube Slicer.

For additional information on fraction recovery systems available from Beckman Coulter, refer to the latest edition of Ultracentrifuge Rotors, Tubes & Accessories (publication BR-8101).

Making Ultra-Clear Tubes Wettable

The following method of making Ultra-Clear tubes wettable has proven successful for some users:

1. Polyvinyl alcohol (2 g) was dissolved in distilled water (50 mL) by stirring and heating to gentle reflux.
2. Isopropanol (50 mL) was slowly added to the hot solution and stirring and heating continued until a clear solution was obtained.
3. The solution was then allowed to cool to room temperature.
4. Ultra-Clear tubes were filled with the coating solution, then aspirated out with a water pump after 15 minutes, leaving a thin film on the tube walls. A small amount of solution that collected in the tube bottoms after standing was removed with a pipette.
5. The tubes were left open to dry at room temperature overnight, then filled with distilled water. After standing overnight at room temperature, the distilled water was poured out.
6. Finally, the tubes were briefly flushed with water, tapped to remove excess liquid, and left to dry.

CHAPTER 4


Using Fixed-Angle Rotors

Introduction

This chapter contains instructions for using fixed-angle rotors in preparative ultracentrifuges. In addition to these instructions, observe procedures and precautions provided in the applicable rotor and ultracentrifuge manuals.

Refer to [CHAPTER 2](#) for labware selection information, and [CHAPTER 3](#) for recommended filling and sealing or capping requirements and for sample recovery procedures. Refer to [CHAPTER 7](#) for information on the care of rotors and accessories.

Description

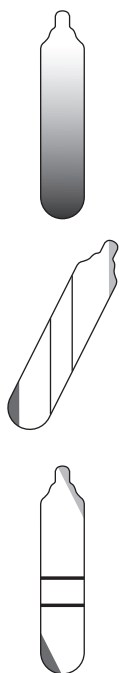

Fixed-angle rotors (see [Figure 4.1](#)) are general-purpose rotors that are especially useful for pelleting and isopycnic separations. Refer to [Table 4.1](#) for general rotor specifications. Tubes in fixed-angle rotors are held at an angle (usually 20 to 35 degrees) to the axis of rotation in numbered tube cavities. The tube angle shortens the particle pathlength compared to swinging bucket rotors, resulting in reduced run times.

Most fixed-angle rotors have a lid secured by a handle. Most handles have holes so that a screwdriver or metal rod can be used to loosen the lid after centrifugation. The lids of some high-performance rotors have either two or four small holes to provide a temporary vent, which prevents rotor damage by allowing liquid to escape in the event of tube leakage.

O-rings, made of Buna N rubber, are located in the rotor lid. The O-rings help to maintain atmospheric pressure inside the rotor during centrifugation, if they are properly lubricated.

Some rotors have fluted bodies, designed to eliminate unnecessary weight and minimize stresses.

Figure 4.1 Fixed-Angle Rotors

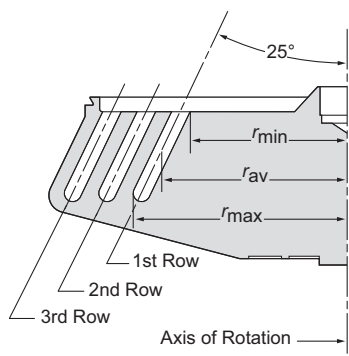
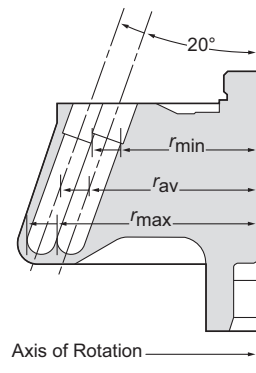
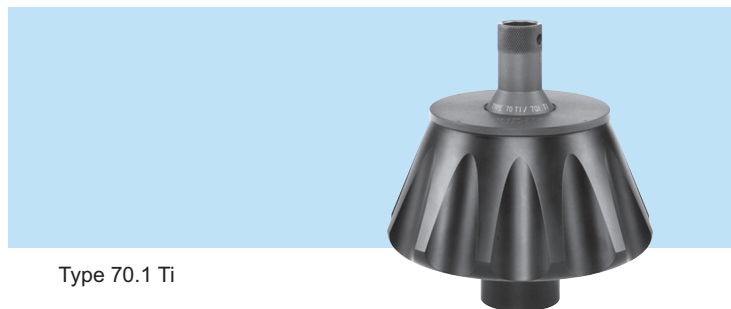
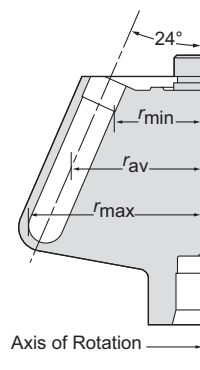
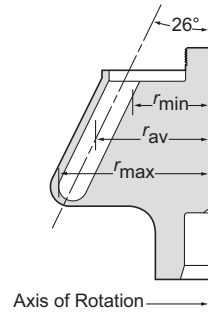


Table 4.1 General Specifications for Beckman Coulter Preparative Fixed-Angle Rotors

Rotor Type	Maximum Speed ^a (rpm)	Relative Centrifugal Field ($\times g$) at $r_{\max}$	Tube Angle (degrees)	Radial Distances (mm)			k Factor	Number of Tubes $\times$ Tube Capacity (mL)
				$r_{\max}$	r_{av}	$r_{\min}$		
100 Ti	100,000	801,920	26	71.6	55.5	39.5	15	8 $\times$ 6.8
90 Ti	90,000	694,000	25	76.5	55.4	34.2	25	8 $\times$ 13.5
70.1 Ti	70,000	450,000	24	82.0	61.2	40.5	36	12 $\times$ 13.5
70 Ti	70,000	504,000	23	91.9	65.7	39.5	44	8 $\times$ 38.5
50.4 Ti	50,000	312,000 ^b	20	111.5	96.2	80.8	33	44 $\times$ 6.5
50.2 Ti	50,000	302,000	24	107.9	81.2	54.4	69	12 $\times$ 39
45 Ti	45,000	235,000	24	103.8	69.8	35.9	133	6 $\times$ 94
42.2 Ti	42,000	223,000	30	113.0	108.5	104.0	9	72 $\times$ 230
25	25,000	92,500 ^c	25	132.1	122.8	113.4	62	100 $\times$ 1
19	19,000	53,900	25	133.4	83.9	34.4	951	6 $\times$ 250

a. Maximum speeds are based on a solution density of 1.2 g/mL in all fixed-angle rotors.

b. Maximum RCF measured at outer row.

c. Maximum RCF measured at the third row. Radial distances are those of the third row.

NOTE Although rotor components and accessories made by other manufacturers may fit in the Beckman Coulter rotor you are using, their safety in the rotor cannot be ascertained by Beckman Coulter. Use of other manufacturers' components or accessories in a Beckman Coulter rotor may void the rotor warranty, and should be prohibited by your laboratory safety officer. Only the components and accessories listed in the applicable rotor manual should be used.

Tubes and Bottles

Fixed-angle rotors can accommodate a variety of tube types, listed in the rotor manual. Refer to [CHAPTER 3](#), for tube filling and sealing or capping requirements. Observe the maximum rotor speeds and fill volumes listed in the applicable rotor instruction manual.

Fill volumes, maximum rotor speeds, and capping requirements for ultracentrifuge bottles are listed in [CHAPTER 3](#). Some rotors must be centrifuged at reduced speeds when bottles are run partially filled. Refer to the applicable rotor manual for specific minimum and maximum fill volumes and rotor speeds.

When running uncapped tubes, observe the maximum rotor speeds and fill volumes listed in [Table 4.2](#).

Rotor Preparation and Loading


For runs at other than room temperature, refrigerate or warm the rotor beforehand for fast equilibration.

Prerun Safety Checks



Read all safety information in the rotor manual before using the rotor.

- 1 Make sure that the rotor and lid are clean and show no signs of corrosion or cracking.
- 2 Make sure the rotor is equipped with the correct overspeed disk (refer to [CHAPTER 1](#)). If the disk is missing or damaged, replace it as described in [CHAPTER 7](#).

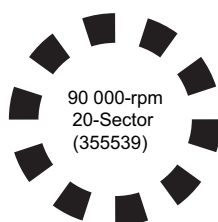

- 3 Check the chemical compatibilities of all materials used. (Refer to IN-175.)
- 4 Verify that tubes, bottles, and accessories being used are listed in the appropriate rotor manual.

Table 4.2 Maximum Run Speeds and Tube Volumes for Uncapped Tubes in Fixed-Angle Rotors

Nominal Dimensions (mm)	Part Number		Maximum Volume (mL)	Maximum Capless Speed ^a (rpm)		Rotor Type ^b
	Polycarbonate	Polypropylene		Polycarbonate	Polypropylene	
7 × 20	343775	343621	230 mL	42,000	42,000	42.2 Ti
7 × 20	342303 ^c	—	230 mL	42,000		42.2 Ti
8 × 51	355657		1	45,000	—	50.4 Ti
				25,000	—	25
13 × 64	355645	355644	4	50,000	30,000	50.4 Ti
16 × 76	355630	355640	10 mL	40,000	40,000	70 Ti, 50.2 Ti
25 × 89	355631	355642	16.5	45,000	20,000	70 Ti, 50.2 Ti
38 × 102	355628	355643	47	15,000	15,000	45 Ti

a. Maximum speeds are those for capless tubes, tested at 25°C for 24 hours.


b. Rotors are not listed for tubes used with adapters.

c. Cellulose propionate

Rotor Preparation and Loading


- 1 Be sure that metal threads in the rotor are clean and lightly but evenly lubricated with Spinkote lubricant (306812).
 - a. Also ensure that O-rings are lightly but evenly coated with silicone vacuum grease (335148).

- 2 Dry the exterior of the tubes.
 - (Moisture between the tube and the rotor cavity may lead to tube collapse and increase the force required to extract the tube.)
 - a. Slide the filled and capped or sealed tubes into the tube cavities.
 - Tubes must be arranged symmetrically in the rotor (see [Figure 1.5](#)).
 - Opposing tubes must be filled to the same level with liquid of the same density.
 - Refer to [Rotor Balance](#) in [CHAPTER 1](#).

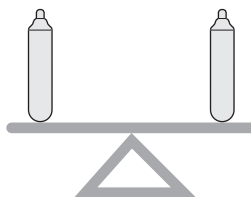


NOTE Place filled tubes in at least two opposing cavities. Make sure that cavities in use also have the proper spacers inserted before installing the rotor lid. *Do not put spacers in cavities that do not contain tubes.*


- 3 Use the required spacers and/or floating spacers, if necessary, to complete the loading operation.
 - a. If OptiSeal tubes are being used, install a spacer over each plugged tube (refer to the applicable rotor manual).
 - b. Leave cavities without tubes completely empty.

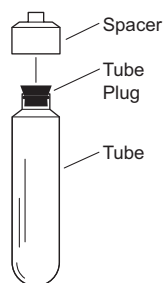

- c. If Quick-Seal tubes are being used, install spacers and/or floating spacers over sealed tubes (refer to the applicable rotor manual).
 - The particular type of tube support for Quick-Seal tubes in fixed-angle rotors depends on the length of the tube, but the top of the tube must be supported.
- d. Leave cavities without tubes completely empty.

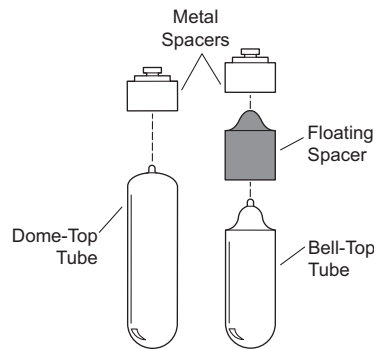

- 4 Place the lid on the rotor and tighten it, as firmly as possible, with the handle.
 - a. Screw the handle down clockwise to fully compress the O-rings

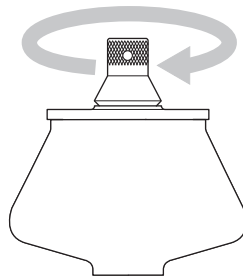

CAUTION


The lid should not touch the tube caps. If the lid touches the caps, the caps are not seated properly on the tubes. Remove the tubes from the rotor and recap them (refer to [CHAPTER 3](#)). Check the tube cavity for foreign matter.

Operation

For runs at other than room temperature, refrigerate or warm the rotor beforehand for fast equilibration.

Installing the Rotor

- 1 Carefully lower the rotor straight down onto the drive hub.
 - a. If the rotor has drive pins, install it so that the pins are at a 90-degree angle to the pins in the drive hub.
 - Careful installation will prevent disturbing the sample or tripping the imbalance detector.
 - b. Refer to the centrifuge instruction manual for detailed operating information.

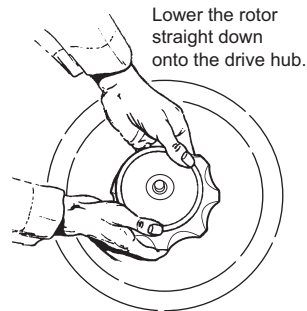


Removal and Sample Recovery

CAUTION

If disassembly reveals evidence of leakage, you should assume that some fluid escaped the rotor. Apply appropriate decontamination procedures to the centrifuge and accessories.

- 1 Remove the rotor from the centrifuge by lifting it straight up and off the drive hub.
- 2 Unscrew the handle counterclockwise and remove the lid.
 - Some rotor handles have holes so that a screwdriver or metal rod can be used to loosen the lid.

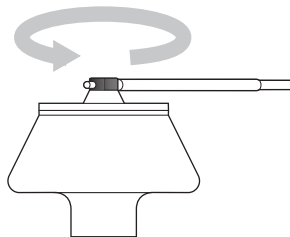


3 Remove spacers and/or floating spacers with a removal tool (338765) or hemostat.

4 Remove tubes or bottles from the rotor using one of the following procedures.

- Refer to [Figure 4.2](#) for removal tools.

NOTE When removing a tube cap, do not remove the cap nut, or the stem may drop into the tube contents and disturb the separation. Instead, loosen the nut just enough to remove the cap assembly as a unit.


- Extract capped tubes using the appropriate removal tool.
 - 1) Insert the threaded end of the tool into the cap and screw at least one turn.
 - If necessary, turn the tube slightly to break any vacuum seal created between the tube and the cavity, and pull the tube out.
 - 2) Use the hex-key end of the removal tool to remove the cap setscrew, but try not to squeeze the tube.
 - With the setscrew removed, supernatant liquid can be withdrawn from the tube, or the tube bottom can be punctured for fraction collection.
- Extract capless tubes using forceps or a hemostat, and OptiSeal or Quick-Seal tubes with the removal tool (361668).
- To remove polycarbonate bottles with black Noryl caps, insert the crossbar end of the removal tool (335381) into the cap slot and turn until the crossbar is past the slot.
 - 1) Pull the bottle out.
- For bottles with red aluminum caps, depress the button of the removal tool (878133) and insert the end of the tool into the cap hole.
 - 1) Release the button and pull the bottle out.

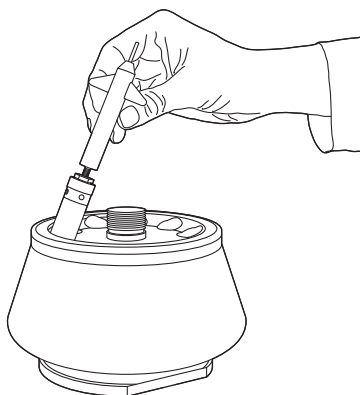

5 Remove adapters using the appropriate removal tool.

6 Refer to [CHAPTER 3](#), for sample recovery methods.

Figure 4.2 Removal Tools Used in Fixed-Angle Rotors

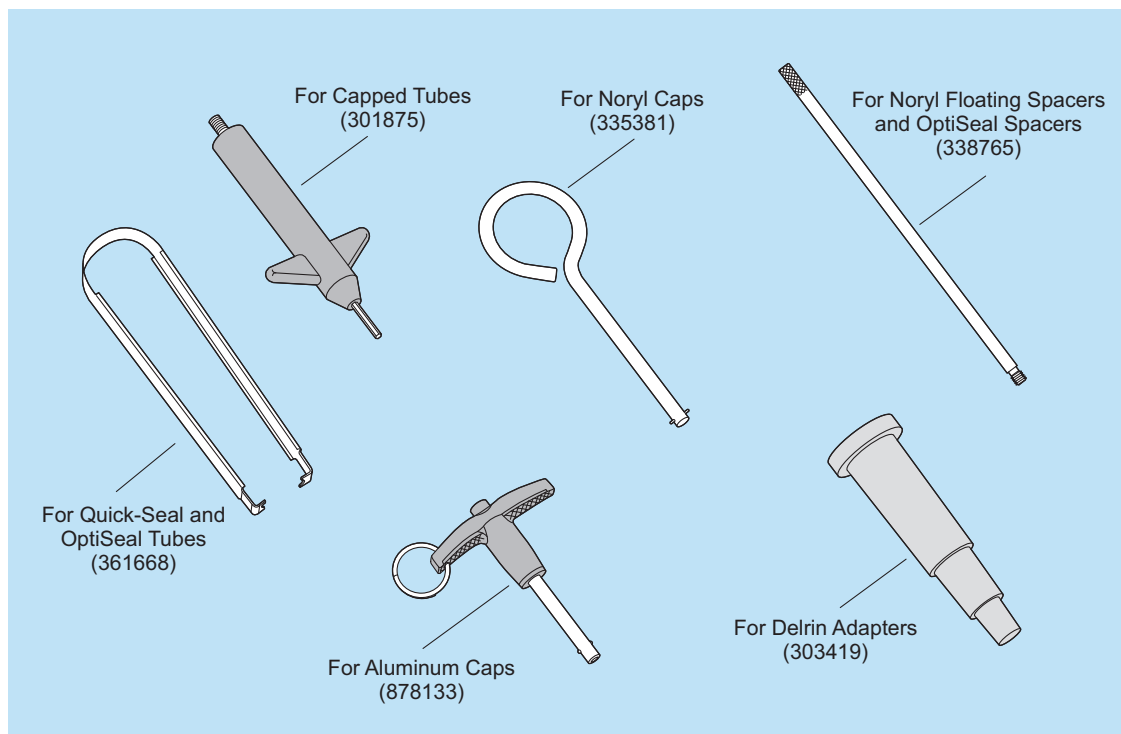


Using Swinging-Bucket Rotors

Introduction

This chapter contains instructions for using swinging-bucket rotors in preparative ultracentrifuges. In addition to these instructions, observe procedures and precautions provided in the applicable rotor and ultracentrifuge manuals.

Refer to [CHAPTER 2](#) for labware selection information, and [CHAPTER 3](#) for recommended filling and sealing or capping requirements and for sample recovery procedures. Refer to [CHAPTER 7](#) for information on the care of rotors and accessories.

Description

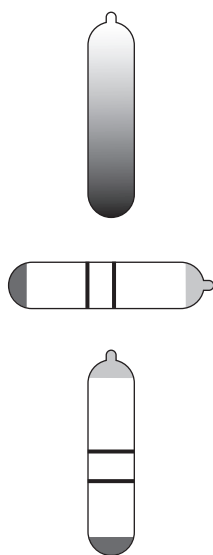


Swinging-bucket rotors (see [Figure 5.1](#)) are most frequently used for density gradient separations, either isopycnic or rate zonal. Refer to [Table 5.1](#) for general rotor specifications.


Tubes in swinging-bucket rotors are held in the rotor buckets. Buckets are attached to the rotor body by hinge pins or a crossbar. The buckets swing out to a horizontal position as the rotor accelerates, then seat against the rotor body for support. Bucket and rotor body positions are numbered for operator convenience.

Each bucket is sealed by an O-ring or gasket between the bucket and the bucket cap. Caps are either a small, flat cap, tightened with a screwdriver, or a cap that is integral with the hanger mechanism, screwed into the bucket by hand.


Figure 5.1 Swinging-Bucket Rotors

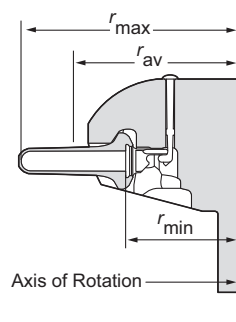


SW 60 Ti

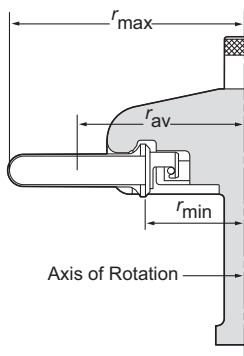

SW 40 Ti

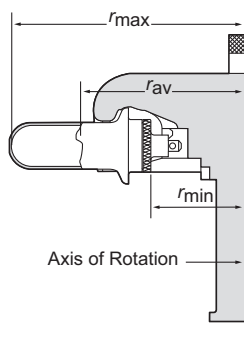

SW 28

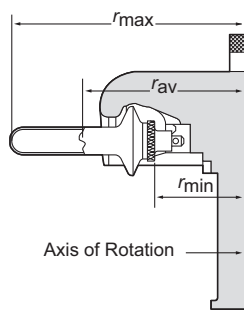

SW 28.1



Table 5.1 General Specifications for Beckman Coulter Preparative Swinging-Bucket Rotors

Rotor	Maximum Speed ^a (rpm)	Relative Centrifugal Field ($\times g$) at $r_{\max}$	Radial Distances (mm)			k Factor	$k\varsigma$ Factors ^b			Number of Tubes $\times$ Tube Capacity (mL)
			$r_{\max}$	r_{av}	$r_{\min}$		(g/mL) $\rho = 1.3$	(g/mL) $\rho = 1.5$	(g/mL) $\rho = 1.7$	
SW 60 Ti	60,000	485,000	120.3	91.7	63.1	45	126	115	111	6×4
SW 55 Ti	55,000	368,000	108.5	84.6	60.8	48	135	123	118	6×5
SW 41 Ti	41,000	288,000	153.1	110.2	67.4	124	335	307	295	6×13.2
SW 40 Ti	40,000	285,000	158.8	112.7	66.7	137	368	338	325	6×14
SW 32 Ti	32,000	175,000	152.5	109.7	66.8	204	468	428	412	6×38.5
SW 32.1 Ti	32,000	187,000	162.8	113.6	64.4	228	613	560	536	6×17
SW 28.1	28,000	150,000	171.3	122.1	72.9	276	757	694	668	6×17
SW 28	28,000	141,000	161.0	118.2	75.3	246	680	622	600	6×38.5

a. Maximum speeds are based on a solution density of 1.2 g/mL in all swinging-bucket rotors.

b. Calculated for 5 to 20% (wt/wt) sucrose at 5°C.

NOTE Although rotor components and accessories made by other manufacturers may fit in the Beckman Coulter rotor you are using, their safety in the rotor cannot be ascertained by Beckman Coulter. Use of other manufacturers' components or accessories in a Beckman Coulter rotor may void the rotor warranty, and should be prohibited by your laboratory safety officer. Only the components and accessories listed in the applicable rotor manual should be used.

Tubes and Bottles

Swinging-bucket rotors can accommodate a variety of tube types, listed in the applicable rotor manual. Refer to [CHAPTER 3](#) for tube filling and sealing or capping requirements. Observe the maximum rotor speeds and fill volumes listed in the rotor instruction manual.

Rotor Preparation and Loading

For runs at other than room temperature, refrigerate or warm the rotor beforehand for fast equilibration.

NOTE All buckets, loaded or empty, must be positioned on the rotor body for every run.


Prerun Safety Checks



Read all safety information in the rotor manual before using the rotor.

- 1 Make sure that the rotor and lid are clean and show no signs of corrosion or cracking.

-
- 2** Make sure the rotor is equipped with the correct overspeed disk (refer to [CHAPTER 1](#)). If the disk is missing or damaged, replace it as described in [CHAPTER 7](#).

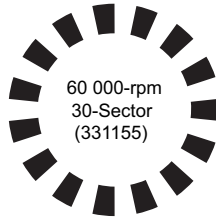


-
- 3** Check the chemical compatibilities of all materials used. (Refer to IN-175.)
-
- 4** Verify that tubes, bottles, and accessories being used are listed in the appropriate rotor manual.
-

Rotor Preparation and Loading

-
- 1** If the rotor has hinge pins, replace any pin that has stripped threads.
-
- 2** Be sure that bucket threads are clean and lightly but evenly lubricated with Spinkote lubricant (306812), as required.

- 3 Remove the bucket gaskets or O-rings and coat them lightly but evenly with silicone vacuum grease (335148).
 - a. Install gaskets or O-rings in the buckets (refer to [Figure 5.2](#)).

Figure 5.2 Install Gaskets Inside the Buckets

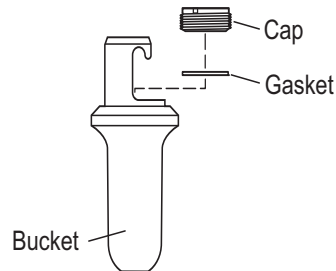


CAUTION

Never run a filled bucket without a gasket or O-ring, as the bucket contents may be lost, leading to rotor imbalance and possible failure.


- 4 Dry the exterior of the tubes.
 - (Moisture between the tube and the bucket may lead to tube collapse and increase the force required to extract the tube.)
 - a. Slide the filled and sealed tubes into the buckets.
 - Loaded buckets can be supported in the bucket holder rack available for each rotor.
- 5 Use the required spacers and/or floating spacers, if necessary, to complete the loading operation.
 - a. If OptiSeal tubes are being used, install a spacer over each plugged tube (refer to the applicable rotor manual and [Figure 5.3](#)).
 - 1) Leave buckets without tubes completely empty.

Figure 5.3 Install a Spacer Over the Plugged Tube

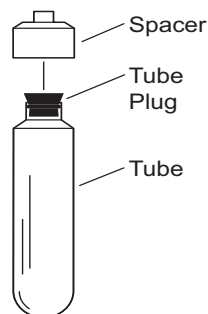

- b. If Quick-Seal tubes are being used, install spacers and/or floating spacers over sealed tubes (refer to the applicable rotor manual and [Figure 5.4](#)).
 - The particular type of tube support for Quick-Seal tubes in swinging-bucket rotors depends on the length of the tube, but the top of the tube must be supported.
- 1) Leave buckets without tubes completely empty.

Figure 5.4 Install Spacers Over Sealed Tubes

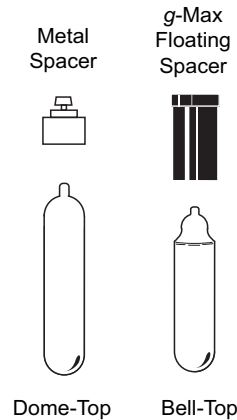

6 Match numbered caps with numbered buckets.

- a. Screw the caps into the bucket until there is metal-to-metal contact.
- b. Tighten flat caps with a screwdriver.


NOTE For SW 32 Ti and SW 32.1 Ti rotors—use a lint-free cotton swab to apply Spinkote lubricant (396812) to cap grooves in the bucket tops. Match bucket caps with numbered buckets. Align the pins on each side of the cap with the guide slots in the bucket. Twist the cap clockwise until it stops (one-quarter turn).

7 Attach all buckets, loaded or empty, to the rotor.

- Loaded buckets must be arranged symmetrically on the rotor (see [Figure 1.5](#)).
- Opposing tubes must be filled to the same level with liquid of the same density.
 - Refer to [Rotor Balance](#) in [CHAPTER 1](#).
- a. If the rotor has hook-on buckets, make certain that both hooks are on the crossbar and that buckets are placed in their proper labeled positions.
- b. If the rotor has hinge pins, lightly lubricate the pin threads with Spinkote.
 - 1) Attach each bucket using the hinge pin tool (330069 and 330070).

NOTE Place filled tubes in at least two opposing buckets. Do not put spacers in buckets that do not contain tubes (refer to [Figure 5.5](#)).

Figure 5.5 Place Filled Tubes in at Least Two Opposing Buckets

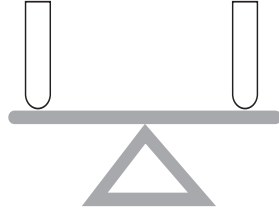


Operation

For runs at other than room temperature, refrigerate or warm the rotor beforehand for fast equilibration.


- 1** Note the location of the two small indentations on the rotor adapter (or the mechanical overspeed devices on older rotors).
 - Their position indicates the location of the drive pins (refer to [Figure 5.6](#)).

Figure 5.6 Rotor Adapter Indentations Indicate Drive Pin Locations

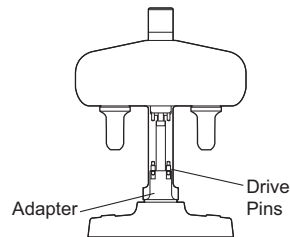

- 2** Carefully lift the rotor with both hands (do not carry a rotor with hook-on buckets by the rotor adapter; the buckets may be dislocated, resulting in an unbalanced rotor, spilled sample, and failed or collapsed tubes) and lower it straight down onto the drive hub (refer to [Figure 5.7](#)).
 - a.** Make sure that the rotor pins are at a 90-degree angle to the drive hub pins.
 - Careful installation will prevent disturbing the sample or tripping the imbalance detector.

Figure 5.7 Lift Rotor With Two Hands and Lower Onto the Drive Hub

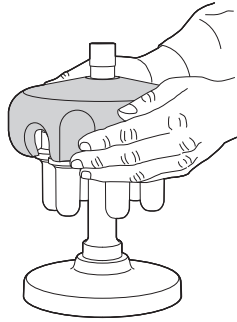

CAUTION

If hook-on buckets have been jarred during installation, check them with a mirror for proper vertical positioning (see [Figure 5.8](#)). Remove the rotor to correct any unhooked buckets.

- 3 Refer to the centrifuge instruction manual for detailed operating information.

Figure 5.8 Checking Hook-on Bucket Positions After the Rotor is Installed*



* Note the partially unhooked bucket on the right.


Removal and Sample Recovery

CAUTION

If disassembly reveals evidence of leakage, you should assume that some fluid escaped the rotor. Apply appropriate decontamination procedures to the centrifuge and accessories.

- 1 Remove the rotor from the centrifuge by lifting it straight up and off the drive hub (refer to [Figure 5.9](#)).


Figure 5.9 Lift Rotor Straight Up

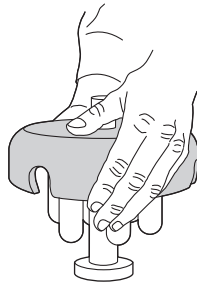

- 2 Set the rotor on the rotor stand and carefully remove the buckets—lift buckets off crossbars or unscrew the hinge pins.
- 3 Remove the bucket caps and use the appropriate removal tool to remove the spacers and tubes.
- 4 Remove adapters using the appropriate removal tool (refer to [Figure 5.10](#)).

NOTE If conical-shaped adapters that support konical tubes are difficult to remove after centrifugation, an extractor tool (354468) is available to facilitate removal.

Figure 5.10 Remove Adapters With the Removal Tool

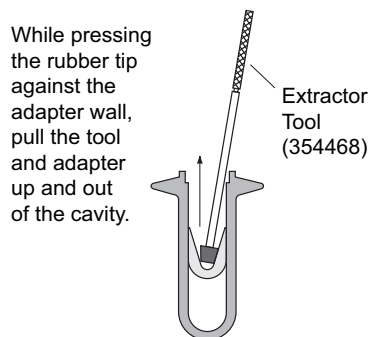

5 Refer to [CHAPTER 3](#) for sample recovery methods.

CHAPTER 6

Using Vertical-Tube and Near-Vertical Tube Rotors

Introduction


This chapter contains instructions for using vertical-tube and near-vertical tube rotors in preparative ultracentrifuges. In addition to these instructions, observe procedures and precautions provided in the applicable rotor and ultracentrifuge manuals.

Refer to [CHAPTER 2](#) for labware selection information, and [CHAPTER 3](#) for recommended filling and sealing or capping requirements and for sample recovery procedures. Refer to [CHAPTER 7](#) for information on the care of rotors and accessories.

Description

Vertical-tube and near-vertical tube rotors are especially useful for isopycnic banding and rate zonal experiments. Some rotors have fluted bodies, designed to eliminate unnecessary weight and minimize stresses. Refer to [Table 6.1](#) for general rotor specifications.

Vertical-Tube Rotors

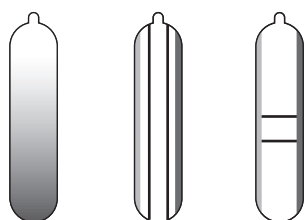

Tubes in vertical-tube rotors (see [Figure 6.1](#)) are held parallel to the axis of rotation in numbered tube cavities. These rotors have plugs that are screwed into the rotor cavities over sealed OptiSeal or Quick-Seal tubes. The plugs (with spacers, when required) restrain the tubes in the cavities and provide support against the hydrostatic force generated by centrifugation.


Table 6.1 General Specifications for Beckman Coulter Preparative Vertical-Tube and Near-Vertical Tube Rotors

Rotor Type	Maximum Speed ^a (rpm)	Relative Centrifugal Field (× g) at r_{max}	Tube Angle (degrees)	Radial Distances (mm)			k Factor	Number of Tubes × Tube Capacity (mL)
				r_{max}	r_{av}	r_{min}		
Vertical Tube								
VTi 90	90,000	645,000	0	71.1	64.5	57.9	6	8 × 5.1
VTi 65.2	65,000	416,000	0	87.9	81.3	74.7	10	16 × 5.1
VTi 65.1	65,000	401,700	0	84.9	76.7	68.5	13	8 × 13.5
VTi 50.1	50,000	251,000	0	89.7	76.8	63.9	34	12 × 39
VTi 50	50,000	242,000	0	86.6	73.7	60.8	36	8 × 39
Near Vertical Tube								
NVT 100	100,000	750,000	8	67.0	57.6	48.3	8	8 × 5.1
NVT 90	90,000	645,000	8	71.1	61.8	52.4	10	8 × 5.1
NVT 65.2	65,000	416,000	8.5	87.9	78.4	68.8	15	16 × 5.1
NVT 65	65,000	402,000	7.5	84.9	72.2	59.5	21	8 × 13.5

a. Maximum speeds are based on a solution density of 1.7 g/mL in all vertical tube and near vertical tube rotors.


NOTE Although rotor components and accessories made by other manufacturers may fit in the Beckman Coulter rotor you are using, their safety in the rotor cannot be ascertained by Beckman Coulter. Use of other manufacturers' components or accessories in a Beckman Coulter rotor may void the rotor warranty, and should be prohibited by your laboratory safety officer. Only the components and accessories listed in the applicable rotor manual should be used.

Near-Vertical Tube Rotors

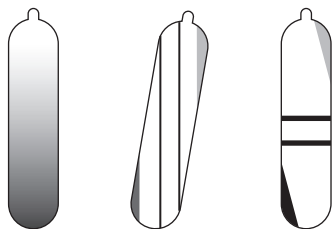


Tubes in near-vertical tube rotors (see [Figure 6.2](#)) are held in numbered tube cavities at an angle to the axis of rotation (typically 7 to 10 degrees). The slight angle of the rotor significantly reduces run times from fixed angle rotors (with tube angles of 20 to 35 degrees) while allowing components that do not band under separation conditions to either pellet to the bottom or float to the top of the tube. Like the vertical-tube rotors, these rotors have plugs to restrain and support sealed OptiSeal or Quick-Seal tubes.

Figure 6.1 Vertical-Tube Rotors

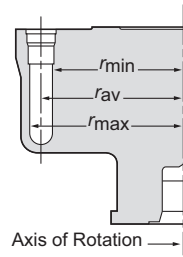


VTi 65.2

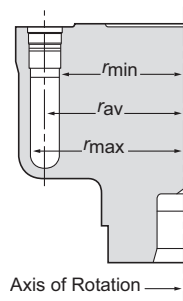


VTi 90



Figure 6.2 Near-Vertical Tube Rotors

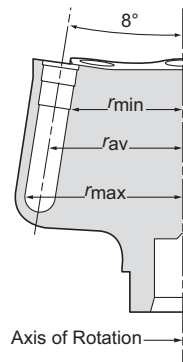

NVT 100

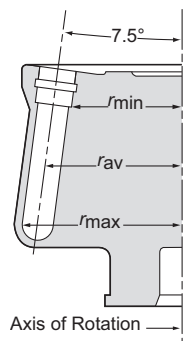

NVT 65



Tubes and Bottles

Only OptiSeal or Quick-Seal tubes are used in these rotors. Refer to [CHAPTER 3](#) for tube filling and sealing or plugging requirements. Observe the maximum rotor speeds and fill volumes listed in the applicable rotor instruction manual.

Rotor Preparation and Loading


For runs at other than room temperature, refrigerate or warm the rotor beforehand for fast equilibration.

Prerun Safety Checks



Read all safety information in the rotor manual before using the rotor.


- 1 Make sure that the rotor, plugs, gaskets, and spacers are clean and show no signs of corrosion or cracking.
 - The high forces generated in these rotors can cause damaged components to fail.
- 2 Make sure the rotor is equipped with the correct overspeed disk (refer to [CHAPTER 1](#)). If the disk is missing or damaged, replace it as described in [CHAPTER 7](#).

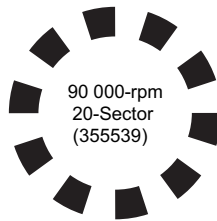


- 3 Check the chemical compatibilities of all materials used. (Refer to IN-175.)
- 4 Verify that tubes, bottles, and accessories being used are listed in the appropriate rotor manual.

Rotor Preparation and Loading

- 1 Be sure that plug threads are clean and lightly but evenly lubricated with Spinkote lubricant (306812).

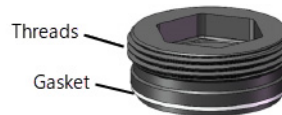


- 2 If using a rotor vise, set the rotor into the vise, which should be bolted or clamped to a rigid surface.
- 3 Dry the exterior of the plugged (OptiSeal) or sealed (Quick-Seal) tubes.
 - (Moisture between the tube and the rotor cavity may lead to tube collapse and increase the force required to extract the tube.)
 - a. Slide the tubes into the tube cavities.
 - Tubes must be arranged symmetrically in the rotor (see [Figure 1.5](#)).
 - Opposing tubes must be filled to the same level with liquid of the same density.
 - Refer to [Rotor Balance](#) in [CHAPTER 1](#).
 - b. Place filled tubes in at least two opposing cavities.

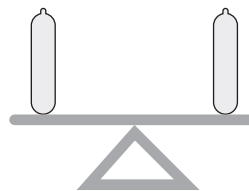


- 4 It is important that each cavity being used is completely filled.
 - Use the required spacers and/or floating spacers, if necessary, to complete the loading operation.
 - a. If OptiSeal tubes are being used, install a spacer over each plugged tube (refer to the applicable rotor manual).

- 1) Leave cavities without tubes completely empty.

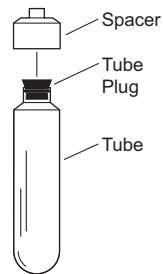

- b. If Quick-Seal tubes are being used, install spacers and/or floating spacers over sealed tubes (refer to the applicable rotor manual).

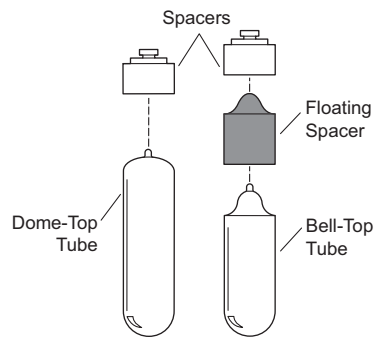

- The particular type of tube support for Quick-Seal tubes depends on the length of the tube, but the top of the tube must be supported.

- 1) Leave cavities without tubes completely empty.

CAUTION

To prevent plug damage, do not put spacers or plugs in cavities that do not contain tubes. Leave unused tube cavities completely empty.


-
- 5 Insert a rotor plug, with the white gasket-end down, over each spacer; screw in the plug.

6 Tighten each rotor plug as shown in Figure 6.3.

- Refer to Table 6.2 for the correct tightening tools and torque values.
- a. To avoid stripping the plugs, apply downward pressure to the plug adapter while tightening the plugs.
 - The top surface of each rotor plug should be flush with the top surface of the rotor.
 - Plugs are not flush on the NVT 65.2 rotor; when properly torqued the plugs should protrude not more than 1 mm above the rotor top surface.
 - Plugs are not flush on the VTi 50.1 rotor; when properly torqued, the plugs will be recessed about 2 mm below the rotor top surface.
- b. Make sure that all plugs are level with each other.

CAUTION

The VTi 50 and VTi 50.1 rotors and rotor plugs must be cooled or warmed to the operating temperature prior to torquing, or leakage may occur.

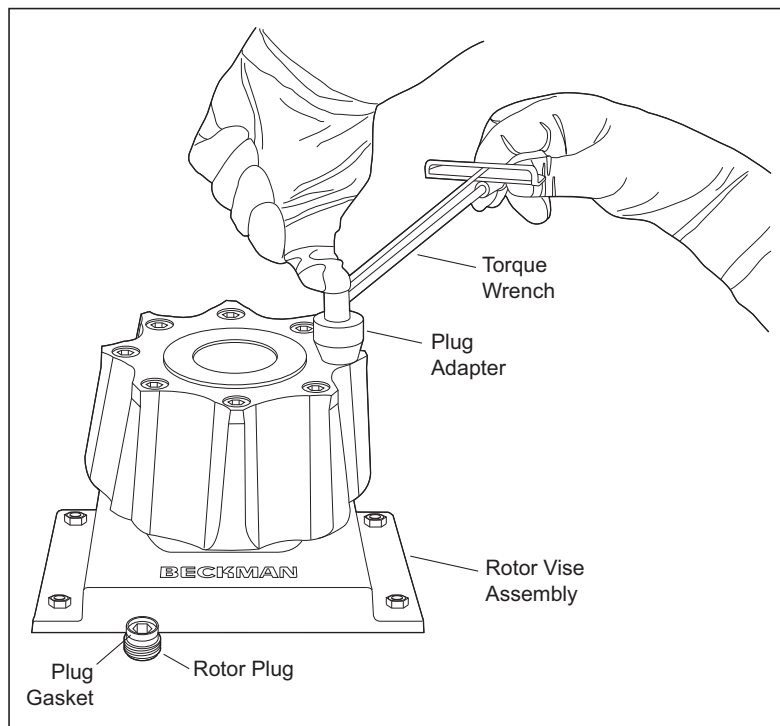
7 Remove the rotor from the vise.**Figure 6.3** Preparing a Vertical-Tube or Near-Vertical Tube Rotor*

* See Table 6.2 for the correct tightening tools and torque values.

Operation

For runs at other than room temperature, refrigerate or warm the rotor beforehand for fast equilibration.

- 1 Carefully lower the rotor straight down onto the drive hub.
 - Careful installation will prevent disturbing the sample or tripping the imbalance detector.

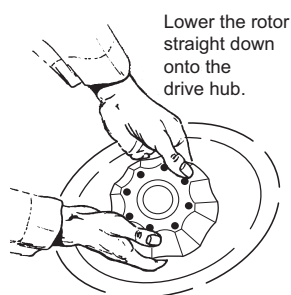

- 2 Refer to the centrifuge instruction manual for detailed operating information.

Table 6.2 Rotor Plugs and Tools Used for Vertical-Tube and Near-Vertical Tube Rotors^a

Rotor	Rotor Plug Part No.	Torque Wrench ^b Part No.	Plug Adapter Part No.	Torque Value
NVT 100	368546	858121	976959	11.3 N•m (100 in.-lb)
NVT 90	368546	858121	976959	13.6 N•m (120 in.-lb)
VTi 90	368546	858121	976959	13.6 N•m (120 in.-lb)
NVT 65.2	368546	858121	976959	13.6 N•m (120 in.-lb)
NVT 65	392084	858121	976959	13.6 N•m (120 in.-lb)
VTi 65.2	368546	858121	976959	13.6 N•m (120 in.-lb)
VTi 65.1	392084	858121	976959	13.6 N•m (120 in.-lb)
VTi 50	355587	889096	355588	17.5 N•m (155 in.-lb) ^c
VTi 50.1	355587	889096	355588	17.5 N•m (155 in.-lb) ^c

a. Rotors listed in parentheses are no longer manufactured.

b. Part number 858121 is a 1/4-in. drive torque wrench; part number 889096 is a 3/8-in. drive torque wrench.


c. The VTi 50 and VTi 50.1 rotors and rotor plugs must be cooled or warmed to operating temperature before torquing or leakage may result.

Removal and Sample Recovery


CAUTION

If disassembly reveals evidence of leakage, you should assume that some fluid escaped the rotor. Apply appropriate decontamination procedures to the centrifuge and accessories.

- 1 Remove the rotor from the centrifuge by lifting it straight up and off the drive hub.
- 2 If a rotor vise is required, set the rotor in the rotor vise.
- 3 Remove the rotor plugs, taking care to apply downward pressure on the plug adapter to avoid stripping the plugs.

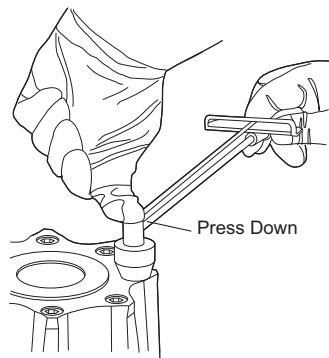

- 4 Remove spacers with the appropriate removal tool or a hemostat.
 - a. Use removal tool (338765) to remove floating spacers.

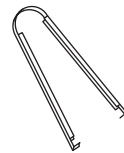

Tube
Removal Tool
(361668)

- 5 Remove tubes with the extraction tool (361668).
- 6 Refer to [CHAPTER 3](#) for sample recovery methods.

CHAPTER 7

Care and Maintenance

Introduction

This chapter provides information on the care of rotors and accessories. Included is a list of some common operating problems with suggestions for their solutions. Rotors and accessories should be kept in optimal condition to minimize the chance of rotor or labware failure. In addition to these instructions, observe procedures and precautions provided in individual rotor manuals. IN-175 provides the chemical resistances of rotor and accessory materials to various acids, bases, salts, and solvents.

Rotor Care

Rotor care involves not only careful operating procedures but also careful attention to:

- Regular cleaning, decontamination, and/or sterilization as required,
- Frequent inspection,
- Corrosion prevention, and
- Regular and proper lubrication.

Do not use sharp tools on a rotor, as the surface can get scratched. Corrosion begins in scratches and may open fissures in the rotor with continued use. The corrosion process accelerates with speed-induced stresses. The potential for damage from corrosion is greatest in aluminum rotors and components.



Wash rotors and rotor components immediately if salts or other corrosive materials are used or if spillage has occurred. **DO NOT** allow corrosive materials to dry on the rotor.

NOTE Do not wash rotor components or accessories in a dishwasher. Do not soak in detergent solution for long periods, such as overnight.


With normal usage, wash rotors frequently to prevent corrosion that can begin in scratches.

CAUTION

Do not immerse or spray a swinging-bucket rotor body with water because liquid can become trapped in the hanger mechanism and lead to corrosion.

- 1 Use plastic or wooden tools to remove O-rings or gaskets for cleaning—do not use metal tools that could scratch anodized surfaces.

- a. Use a mild detergent such as Beckman Solution 555 (339555), diluted 10 to 1 with water, and a soft brush to wash rotors and rotor components and accessories.
 - (Most laboratory detergents are too harsh for aluminum rotors and components.)
 - The Rotor Cleaning Kit (339558) contains two quarts of Solution 555 and brushes that will not scratch rotor surfaces.

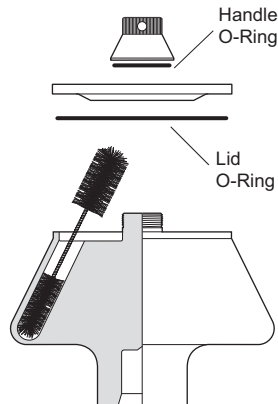

- 2 Rinse thoroughly with water.
- 3 Air-dry the body or buckets upside down.
 - a. Do not use acetone to dry rotors.
- 4 Wipe clean the O-rings or gaskets regularly (lubricate after cleaning).
 - a. Replace them about twice a year or as required.
- 5 Frequently clean all surfaces that contact O-rings.
 - a. Regularly clean the threads of the rotor (lid, handle, buckets, cavities, and so on) with a nonmetal brush and a small amount of concentrated detergent, then rinse, and dry thoroughly.
 - b. Lubricate the threads as directed under [Lubrication](#), below.

Decontamination



Rotors contaminated with radioactive or pathogenic materials must be decontaminated, following appropriate laboratory safety guidelines and/or other regulations.

NOTE Strong bases and/or high-pH solutions can damage aluminum rotors and components.

If a rotor (and/or accessories) becomes contaminated with radioactive material, it should be decontaminated using a solution that will not damage the anodized surfaces. Beckman Coulter has tested Radiacwash.*


While Beckman Coulter has tested this method and found that it does not damage components, no guarantee of decontamination is expressed or implied. Consult your laboratory safety officer regarding the proper decontamination methods to use.

If the rotor or other components are contaminated with toxic or pathogenic materials, follow appropriate decontamination procedures as outlined by appropriate laboratory safety guidelines and/or other regulations. Consult IN-175 to select an agent that will not damage the rotor.

Sterilization and Disinfection

When sterilization or disinfection is a concern, consult your laboratory safety officer regarding proper methods to use. While Beckman Coulter has tested the following methods and found that they do not damage the rotor or components, no guarantee of sterility or disinfection is expressed or implied.

- Rotors and most rotor components (except those made of Noryl) can be autoclaved at 121°C for up to an hour. Remove the lid, bucket caps, or rotor plugs and place the rotor (and/or buckets) in the autoclave upside-down. (O-rings and gaskets can be left in place on the rotor.
- Ethanol (70%)[†] may be used on all rotor components, including those made of plastic. Bleach (sodium hypochlorite) may be used, but may cause discoloration of anodized surfaces. Use the minimum immersion time for each solution, per laboratory standards.

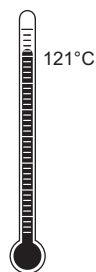

Inspection

Frequent and thorough inspection is crucial to maintaining a rotor in good operating condition.


- 1 Periodically (at least monthly, depending on use) inspect the rotor, especially inside cavities and buckets, for rough spots, cracks, pitting, white powder deposits on aluminum rotors (frequently aluminum oxide), or heavy discoloration.
 - a. If any of these signs are evident, do not run the rotor.

* In U.S., contact Biorad Medical Systems (Shirley, New York); internationally, contact the U.S. office to find the dealer closest to you.

† Flammability hazard. Do not use in or near operating ultracentrifuges.


- b.** Contact your Beckman Coulter representative for information about the Field Rotor Inspection Program and the Rotor Repair Program.

- 2** Regularly check the condition of O-rings or gaskets and replace any that are worn or damaged.

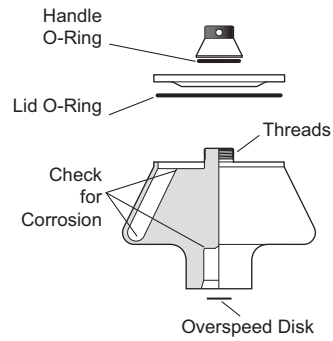

- 3** Regularly check that all sealing surfaces are smooth and undamaged to ensure proper sealing.

- 4 Regularly check the condition of rotor plugs (a component of vertical-tube and near-vertical tube rotors) and rotor plug gaskets.
 - a. Replace worn or damaged gaskets

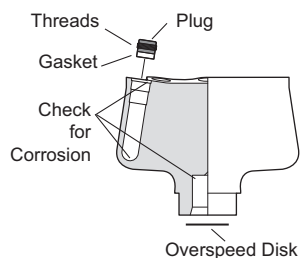

- 5 Regularly inspect the overspeed disk.
 - a. If it is scratched, damaged, or missing, replace it.

Field Rotor Inspection Plan

The Field Rotor Inspection Program (FRIP) has two purposes:

- to prevent premature rotor failures by detecting conditions such as stress, corrosion, metal fatigue, damage, or wear in the anodized coatings; and
- to instruct laboratory personnel in the proper care of rotors.

Beckman Coulter has trained a group of experienced service engineers in the techniques of nondestructive evaluation. For more information about the program, contact your Beckman Coulter representative.

Lubrication


Proper lubrication is essential to obtain specified torque values, where required, and to minimize thread wear.

- Many rotors use O-rings as seals to maintain atmospheric pressure in the rotor during a run. These O-rings and the surfaces they bear against must be kept clean and evenly lubricated. After removing and cleaning rotor O-rings or gaskets, lightly but evenly coat them with silicone vacuum grease (335148) and reposition them in the rotor.
- After cleaning metal threads, lubricate them with Spinkote lubricant (306812). *Failure to keep threads properly lubricated can result in stripped or galled threads and stuck rotor components.*
- Rotor plug gaskets (a component of vertical-tube and near-vertical tube rotors) do NOT require lubrication, but should be checked, cleaned, and or replaced as required.

Overspeed Disk Replacement

The overspeed disk on the rotor bottom is part of the photoelectric overspeed detection system. Replace this disk if it is scratched, damaged, or missing. Start with a dry rotor at room temperature—the disk will not adhere to a damp surface.

- 1 Pry up the edges of the old disk with a scalpel, taking care not to scratch the rotor, then peel the disk off.
- 2 Clean the area around the drive hole with acetone to remove any of the old adhesive.
- 3 Insert the centering tool (331325) into the hole.

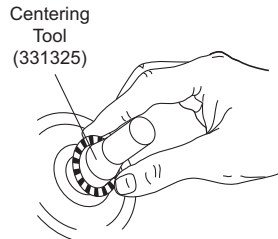

- 4 Peel the paper backing off the new disk, but do not touch the adhesive.
 - a. Fit it, adhesive-side down, around the centering tool.
 - b. Press the disk firmly to the rotor bottom.
- 5 Remove the tool.
 - a. Allow the disk to set for a minimum of 2 hours.

Tube, Bottle, and Accessory Care

Proper care of tubes and bottles involves observing temperature, fill volume, and run speed limitations as well as careful cleaning and sterilization procedures.

Cleaning

Do not wash tubes and bottles in a commercial dishwasher — detergents and temperatures are too harsh.

- Wash tubes, bottles, adapters, and other accessories by hand, using a mild detergent, such as Solution 555 (339555) diluted 10 to 1 with water, and a soft brush.




- Polycarbonate bottles and tubes are vulnerable to attack by alkaline solutions and detergents, so use a detergent with pH less than 9, such as Solution 555. Do not use a brush with exposed metal; scratches in polycarbonate will cause early failure.
- Alcohol and acetone react unsatisfactorily with many tube and accessory materials. If a solvent must be used to rinse, dry, or decontaminate these materials, consult IN-175 to select an appropriate solvent.
- Do not dry tubes, bottles, or accessories in an oven. Labware should be air-dried.
- OptiSeal, Quick-Seal, Ultra-Clear, and thinwall polypropylene tubes are intended for one-time use and should be discarded after use.

Decontamination



Labware contaminated with radioactive or pathogenic solutions should be decontaminated or disposed of following appropriate safety guidelines and/or regulations. Consult IN-175 to select an agent that will not damage the tube or bottle material.

Sterilization and Disinfection

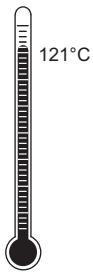

Where sterilization is critical in your application, consider using Beckman Coulter Certified Free & Sterilized Tubes. Beckman Coulter Certified Free & Sterile Tubes are sterilized via Ethylene Oxide in compliance with ISO 11135. Cartons include several peel packages, each containing a typical run quantity of tubes. Packaging meets requirements of ISO 11607. For tubes not available in the sterilized option, refer to [Table 7.1](#), which provides sterilization methods recommended for each container type.

Most tubes and accessories, except those made of Ultra-Clear, polyethylene, Noryl, or cellulose propionate, can be autoclaved at 121°C for about 20 minutes. Note that autoclaving reduces the lifetime of polycarbonate tubes. Also, polypropylene tubes may be permanently deformed if they are autoclaved many times or if they are handled or compressed before they cool. Tubes and bottles should be placed open-end down or supported in a rack if autoclaved. Do not autoclave plastic adapters or spacers.

CAUTION

Do not autoclave tubes or bottles with caps on. Pressure in a sealed container can cause an explosion. Pressures within the autoclave can cause partially sealed containers to collapse when the autoclave vents.

A cold sterilization method, such as immersion in 10% hydrogen peroxide for 30 minutes, may be used on Ultra-Clear tubes. Refer to [Table 7.1](#) to select cold sterilization materials that will not damage tubes and accessories.

While Beckman Coulter has tested these methods and found that they do not damage the components, no guarantee of sterility or disinfection is expressed or implied. When sterilization or disinfection is a concern, consult your laboratory safety officer regarding proper methods to use.

Table 7.1 Tube and Bottle Sterilization and Disinfection^a

Tube/Bottle Material	Autoclave ^b (121°C)	UV Irradiation	Ethylene Oxide	Formaldehyde	Ethanol (70%) ^c	Sodium Hypochlorite (10%)	Hydrogen Peroxide (10%)	Glutaraldehyde (2%)	Phenolic Derivatives
polypropylene	yes	no	yes	yes	yes	yes	yes	yes	no
Ultra-Clear	no	no	yes	yes ^d	yes	yes	yes	yes	no
polycarbonate	yes ^e	no	yes	yes ^d	no	yes ^f	yes	yes	no
polyethylene	no	no	yes	yes	yes ^g	yes	yes	yes	yes
cellulose propionate	no	no	no	no	no	yes	yes	yes	no
stainless steel	yes	yes	yes	yes	yes ^h	no	yes	yes	no

a. This information is provided as a guide to the use of sterilization and disinfection techniques for tube materials. Cold sterilization results shown are for short-duration (10-minute) soak periods; reactions may differ with extended contact. Refer to IN-175 for information about specific solutions.

b. To avoid deformation, autoclave tubes or bottles open-end down in a tube rack at 15 psig for no more than 20 minutes (allow to cool before removing from tube rack). DO NOT autoclave capped or sealed tubes or bottles.

c. Flammable; do not use in or near operating ultracentrifuges.


d. Do not use if there is methanol in the formula.

e. Tube life will be reduced by autoclaving.

f. Discoloration may occur.

g. Below 21°C only.

h. Marginal.

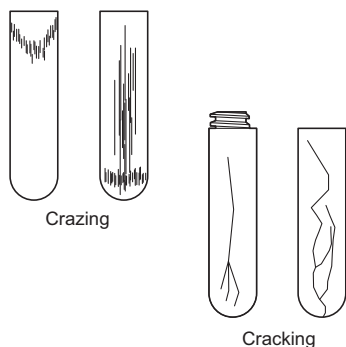

Inspect containers and accessories before use.

- Inspect tubes and bottles for cracks or any major deformities before using them.
- Do not use a tube that has become yellowed or brittle with age or excess exposure to ultraviolet light.
- Crazing—the appearance of fine cracks on tubes and bottles—is the result of stress relaxation. If a crack approaches the outer wall of the tube or bottle, discard it.
- Discard any deformed or cracked adapters.

Tube and Bottle Storage

Tubes and bottles have an indefinite shelf life if properly stored. Store in a dark, cool, dry place away from ozone, chemical fumes, and ultraviolet light sources.

Removing Jammed or Collapsed Tubes

Centrifugal force may collapse improperly sealed or capped thinwall tubes. Observe careful filling and capping procedures to prevent tube collapse.

NOTE Centrifugation often causes a slight vacuum to build up in the tube cavity, occasionally resulting in a suction effect when removing the tubes from the rotor. This effect is especially pronounced in a rotor that has been centrifuged at low temperature. A brief delay (approximately 5 minutes) after the rotor comes to rest before removing the tubes can make tube removal easier. If tubes are difficult to remove from the rotor, use a gentle twisting or rocking motion, and remove the tube slowly to avoid sample mixing.

If a tube is jammed or collapsed in the rotor, try one of the following techniques, but **DO NOT** force the tube. Contact your Beckman Coulter Service representative if you are unsuccessful.



Do not use a hemostat or any metal tool to pry a jammed or collapsed tube out of the rotor. The rotor can be scratched and damaged.

- If an uncapped polycarbonate tube is stuck, remove tube contents and place the rotor or bucket upside-down in an autoclave for about 30 to 60 minutes. When the rotor is cool enough to handle, try to remove the jammed or collapsed tube. Do not autoclave sealed or capped tubes or bottles.
- Pour a solvent in the tube to make the tube material more flexible. Several changes of solvent may be necessary to weaken the tube for easy removal. Refer to the chemical resistances list in IN-175 to select a solvent that will not damage the rotor.

Tube Cap Care

It is very important to keep tube-cap assemblies together as a unit. Do NOT interchange cap components; caps are designed as a unit for a particular tube being centrifuged in a particular rotor. If cap components are separated for cleaning, be sure components are classified according to the tube and rotor for which they are designed. Do not store O-rings or gaskets under compression.

Cleaning

- 1 Disassemble tube caps and wash them in a mild detergent solution, such as Beckman Solution 555 (339555), diluted 10 to 1 with water.
 - a. If necessary, scrub the inside of caps using a cotton-tipped swab or a brush that will not scratch the surface.

NOTE Do not soak aluminum cap parts in a strong detergent solution, as the anodizing may be attacked.

- 2 Clean the nut and stem threads regularly with concentrated Solution 555 and a brush.



- 3 Rinse all parts in distilled water and dry them.
 - 4 Apply a thin, even coat of Spinkote lubricant (306812) to the stem threads.
 - 5 Wipe O-rings and gaskets clean with a tissue.
 - a. Do not lubricate O-rings or gaskets.
-

Decontamination




If the tube caps become contaminated with radioactive material, they should be decontaminated using a solution that will not damage the anodized surfaces. Beckman Coulter has tested Radiacwash.*

While Beckman Coulter has tested this method and found that it does not damage components, no guarantee of decontamination is expressed or implied. Consult your laboratory safety officer regarding the proper decontamination methods to use.

If tube caps are contaminated with toxic or pathogenic solutions, decontaminate or dispose of them as directed by your laboratory safety officer, following appropriate safety guidelines. Check the chemical resistances list in IN-175 to be sure the decontamination method will not damage any part of the rotor.

Sterilization and Disinfection

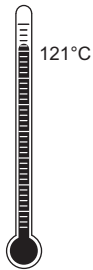

All cap components (except those made of Noryl) can be autoclaved at 121°C for up to 30 minutes. Disassemble caps for autoclaving.

Ethanol (70%)[†] or hydrogen peroxide (6%) may be used on cap components, including those made of plastic. Bleach (sodium hypochlorite) may be used, but may cause discoloration of anodized surfaces. Use the minimum immersion time for each solution, per laboratory standards.

While Beckman Coulter has tested these methods and found that they do not damage components, no guarantee of sterility or disinfection is expressed or implied. When sterilization or disinfection is a concern, consult your laboratory safety officer regarding proper methods to use.

* In U.S., contact Biorad Medical Systems (Shirley, New York); internationally, contact the U.S. office to find the dealer closest to you.

† Flammability hazard. Do not use in or near operating ultracentrifuges.

Lubrication


Keep the stem threads lightly lubricated with Spinkote lubricant (306812). Clean, lubricated threads can be fully tightened without galling or seizing.

The O-ring or gasket must be used dry and without lubrication. A wet or greased O-ring or gasket may allow the stem to rotate when the cap nut is tightened, preventing proper sealing of the cap.


Inspection

Inspect tube-cap components before each use. Refer to [Table 3.3](#) in [CHAPTER 3](#) of this manual for a list of cap components.

- Carefully inspect the crown for deformed or roughened edges. Run your finger around the bottom edge of the crown; surfaces should be flat, squared-off, and not rounded or jagged. Check the top of the crown and the base of the O-ring groove for fine, circular lines or stress cracks. (Do not use damaged wrenches or hex drivers or tools that have burrs. A burred tool can score the crown.) Discard a damaged crown, as it may fail and damage the rotor.

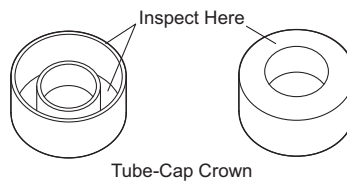

- Inspect the cap stem for evidence of stress cracking. Also, make sure that threads are in good condition and properly lubricated before use. Look at the underside of the stem; the white nylon insert should not protrude below the filling hole. If it does, replace the nylon insert (see replacement procedures below)

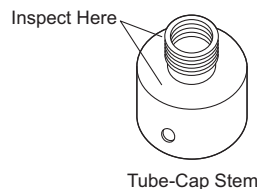


- Check the O-ring or gasket for cuts, excessive abrasions, or flattened areas. It is good practice to replace the O-ring or gasket frequently.
- On caps with filling holes, inspect the filling hole setscrew and threads. If the hex cavity in the setscrew shows signs of wear, replace the setscrew.

Nylon Insert Replacement

The nylon insert fits into the cap stem and can become worn or loose with continued cap use. If the insert is worn, the setscrew will no longer seal the cap and it should be replaced.

- 1 Remove the cap setscrew and fit the nylon insert tool (302460) firmly into the stem.

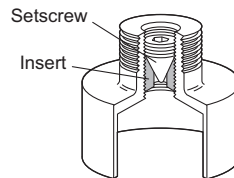

- 2 Unscrew the insert.
- 3 Fit a new insert (302312) on the end of the tool and screw it into the stem until it bottoms firmly against the stem threads.

Returning a Rotor or Accessory to the Factory

Before returning a rotor or accessory for any reason, prior permission must be obtained from Beckman Coulter, Inc. This form may be obtained from your local Beckman Coulter sales office. The form, entitled Returned Material Authorization (RMA) for United States returns or *Returned Goods Authorization* (RGA) for international returns, should contain the following information:

- rotor type and serial number,
- history of use (approximate frequency of use),
- reason for the return,
- original purchase order number, billing number, and shipping number, if possible,
- name and email address of the person to be notified upon receipt of the rotor or accessory at the factory, and
- name and email address of the person to be notified about repair costs, etc.

To protect our personnel, it is the customer's responsibility to ensure that the parts are free from pathogens, chemical hazards, and/or radioactivity. Sterilization and decontamination **MUST** be done before returning the parts. Smaller items (such as tubes, bottles, and so on) should be enclosed in a sealed plastic bag.

All parts must be accompanied by a note, plainly visible on the outside of the box or bag, stating that they are safe to handle and that they are not contaminated with pathogens, chemical hazards, or radioactivity. Failure to attach this notification will result in return or disposal of the items without review of the reported problem.

Use the address label printed on the RMA/RGA form when mailing the rotor and/or accessories.
Customers located outside the United States should contact their local Beckman Coulter office.

Diagnostic Hints

Some of the more common operating problems experienced in centrifugation are listed [Table 7.2](#) with suggestions for their solutions. Contact Beckman Coulter Field Service if a problem cannot be corrected.

NOTE Use only the labware listed in the applicable rotor manual.

Table 7.2 Troubleshooting Chart

Symptom	Possible Cause and Suggested Action
Rotors	
Severe vibration	• Rotor imbalance. To balance the rotor load, fill all opposing tubes to the same level with liquid of the same density. Weight of opposing tubes must be distributed equally. Place tubes in a fixed angle, near vertical tube, or vertical tube rotor symmetrically, as illustrated in CHAPTER 1 (Figure 1.5). • Swinging-bucket rotor — Mishooked bucket, loose bucket cap, wrong type of bucket, mixed bucket types, opposing buckets not filled to the same level with liquids of the same density. Check loading procedures (refer to CHAPTER 5).
Stripped rotor plugs on vertical-tube or near-vertical tube rotors	Rotor vise not used, wrong tool used, incorrect torque, or insufficient pressure on plug adapter, when tightening rotor plugs. Observe careful tightening procedures.
Rotor lid is difficult to remove after centrifugation	Threads contaminated with dirt, dried lubricant, or metal particles, or threads insufficiently lubricated cause rotor components to stick. Do not use excessive force to loosen components. Contact your Beckman Coulter representative. Routinely clean metal threads with concentrated Solution 555 (339555), then lubricate them with Spinkote (306812).
Paint coming off where bucket contacts rotor pocket on swinging-bucket rotor	Not an operational problem.
Tubes	
Tube leakage	
Tubes with cap assemblies	• Caps not properly secured. Caps must be properly seated on tubes and then fully tightened. • Cap components not dry before assembly. Thoroughly dry all components before assembling. • The setscrew may not be sealing the filling hole. The nylon insert may have been driven out by the filling hole setscrew. Check hex cavity. If the threads of the screw are stripped, replace the screw. It may be necessary to replace the stem also. The interface between the setscrew and the nylon insert is critical. Refer to insert replacement procedures in this section. • Insufficient liquid in tube. Observe minimum fill volumes.

Table 7.2 Troubleshooting Chart (*Continued*)

Symptom	Possible Cause and Suggested Action
Tubes with snap-on caps	Tube too full; the meniscus must be kept lower to prevent leakage.
Uncapped tubes	Tube volume exceeds maximum uncapped volume. Refer to the rotor manual for tube volumes and speed reductions.
OptiSeal tubes	Improperly plugged. Make sure that no fluid is trapped in the tube stem, and that the stem is clean and dry before inserting plug. (Refer to publication IN-189 for instructions on filling and plugging OptiSeal tubes.)
Quick-Seal tubes	Improperly sealed. After heat-sealing, squeeze the tube gently (if the tube contents may be disturbed) to test the seal for leaks. If the tube leaks, reseal it.
Tube cracking	• Tubes may crack or become brittle if they are used below their lower temperature limit. Before using tubes at other than stated temperature limits, evaluate them under centrifugation conditions. If sample is frozen in tubes, make sure that tubes are thawed to at least 2°C before centrifugation. • Tubes may become brittle with age and use. Dispose of brittle or cracked tubes.
Tube collapse	• Thinwall tube volume too low to provide tube wall support. Meniscus should be 2 to 3 mm below the tube top. Refer to the rotor manual for tube volumes. • Moisture between the tube and the cavity or bucket can cause the tube to float and collapse. Ensure that tubes and tube cavities or buckets are dry before inserting the tubes. • Reagent used that attacks the tube material. Refer to IN-175 for chemical compatibilities of tube material and chemicals. • Tubes run above their rated speed. Refer to the applicable rotor manual for maximum speeds.
Tube Caps	
Unsure of cap components	For a complete list of cap components, see the Beckman Coulter Ultracentrifuge Rotors, Tubes & Accessories catalog (publication BR-8101), available at www.beckman.com .
Setscrew is difficult to remove	The hex socket or threads of the screw may be stripped. If the screw cannot be removed, replace the cap stem.
Setscrew will not seal the tube cap	Replace the screw and nylon insert if either seems damaged or loose.
Bottles	
Bottle leakage (bottles with cap assemblies)	• Moisture or lubrication on cap or sealing surface. Ensure that the O-ring, plug, and bottle lip are dry and free of lubrication before use. • O-ring or gasket damaged or defective. Replace the O-ring or gasket. • Cap not tightened sufficiently. Tighten cap securely. — Sealing surface of the bottle is not smooth. Replace bottle.

Table 7.2 Troubleshooting Chart (*Continued*)

Symptom	Possible Cause and Suggested Action
Bottle leakage (uncapped bottles)	Bottle too full; the meniscus must be kept lower to prevent leakage. Refer to the rotor manual for fill volumes and speed reductions.
Bottle damage	• Fill volume too low to provide tube wall support. Refer to the rotor manual for fill volumes and speed reduction. • Moisture between the bottle and the cavity or bucket can cause the bottle to float and collapse. Ensure that bottles and cavities or buckets are dry before inserting them. • Reagent used that attacks the bottle material. Refer to IN-175 for chemical compatibilities of bottle material and chemicals. • Bottles may crack or become brittle if they are used below their lower temperature limit. Before using bottles at other than stated temperature limits, evaluate them under centrifugation conditions. If sample is frozen in bottles, make sure that bottles are thawed to at least 2°C before centrifugation. • Bottles may become brittle with age and use. Dispose of brittle or cracked bottles. • Improper cleaning, decontamination, or sterilization procedures used. Refer to Table 7.1 for acceptable procedures and materials.

Use of the ω^2t Integrator

ω^2t Integrator

The centrifugal force applied to a sample in a spinning rotor is shown by ω^2r where r is the radial distance from the axis of rotation and ω is the angular velocity in radians per second ($\omega = 2\pi \text{ rpm}/60$). The sedimentation velocity (dr/dt) is proportional to the centrifugal force; the velocity of a sedimenting particle increases as it moves outward in the tube. Thus, a force-corrected velocity is used to describe the movement of particles under centrifugal force. This is the sedimentation coefficient s , defined as the sedimentation velocity per unit of centrifugal force:

$$s = \frac{dr}{dt} \times \frac{1}{\omega^2 r} \quad (\text{B-1})$$

The integrated form of the equation is:

$$s = \frac{\ln\left(\frac{r_2}{r_1}\right)}{\omega^2 t} \quad (\text{B-2})$$

where r_1 is the initial position of the particle and r_2 is the final position, relative to the axis of rotation. These distances can be readily determined. However, an accurate measure of the centrifugal force applied to the particle necessitates that the value of ω generated during periods of changing speed be calculated, that is from the time when the rotor starts spinning (t_1) until the rotor stops (t_2). The ω^2t integrator automatically computes the total centrifugal effect—during acceleration, constant speed operation, speed changes, and deceleration—and displays this as a continuously updated value of

$$\int_{t_1}^{t_2} \omega^2 dt \quad (\text{B-3})$$

There are two kinds of experiments in which the integrator is particularly useful: duplicating conditions in a series of rate runs, and calculating sedimentation coefficients for rate zonal studies. To duplicate band positions, use the integrator to automatically terminate the run at a preselected value of ω^2t . In this way, even if the set run speed or acceleration is changed for a rotor, band positions will be reproducible. For determining sedimentation coefficients, the value of ω^2t displayed on the integrator at the termination of the run greatly simplifies the arithmetic involved and improves the final result.

Reproducing Band Positions (Refer to Figure B-1)

To achieve the best resolution of particle zones, the centrifugation duration should be set so that the fastest moving zone of particles will move as close as practical to the bottom of the gradient. To determine the centrifugation duration, the following must be known: an estimate of the sedimentation coefficient of the particle of interest, the distance from the axis of rotation it is to travel, its density, and certain properties of the gradient. For example, to position a protein sample characterized by s of 7×10^{-13} seconds (or 7 S) and density of 1.4 g/mL 37 mm down the length of the centrifuge tube in the SW 60 Ti rotor (37 + 63 mm,* or 100 mm from the axis of rotation) through a 10 to 30% gradient at 20°C, the value of $s\omega^2t$ must be 0.92 (from the figure).



$$\omega^2t = \frac{0.92}{7 \times 10^{-13}} \quad (\text{B-4})$$

$$\omega^2t = \frac{1.31 \times 10^{12} \text{ rad}^2}{s}$$

This value can be set into the integrator, and the integrator used to terminate the run when this value is reached. Because of deceleration, however, this value of ω^2t will actually be a little too much. For a more exact approximation, you should make a trial run with an empty rotor and measure the value of ω^2t that accumulates during deceleration from run speed, then subtract that value from the total determined from the charts.

* The radial distance to the tube meniscus in the SW 60 Ti rotor is 63 mm.

Figure A.1 The $s\omega^2t$ Charts for the SW 60 Ti Rotor

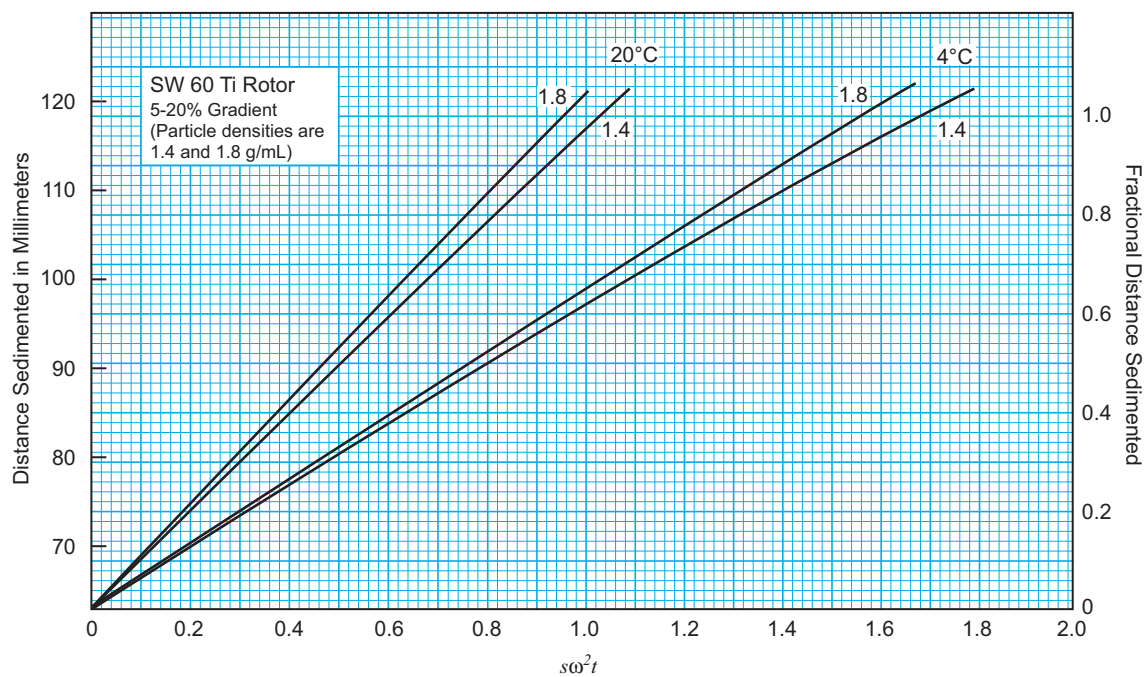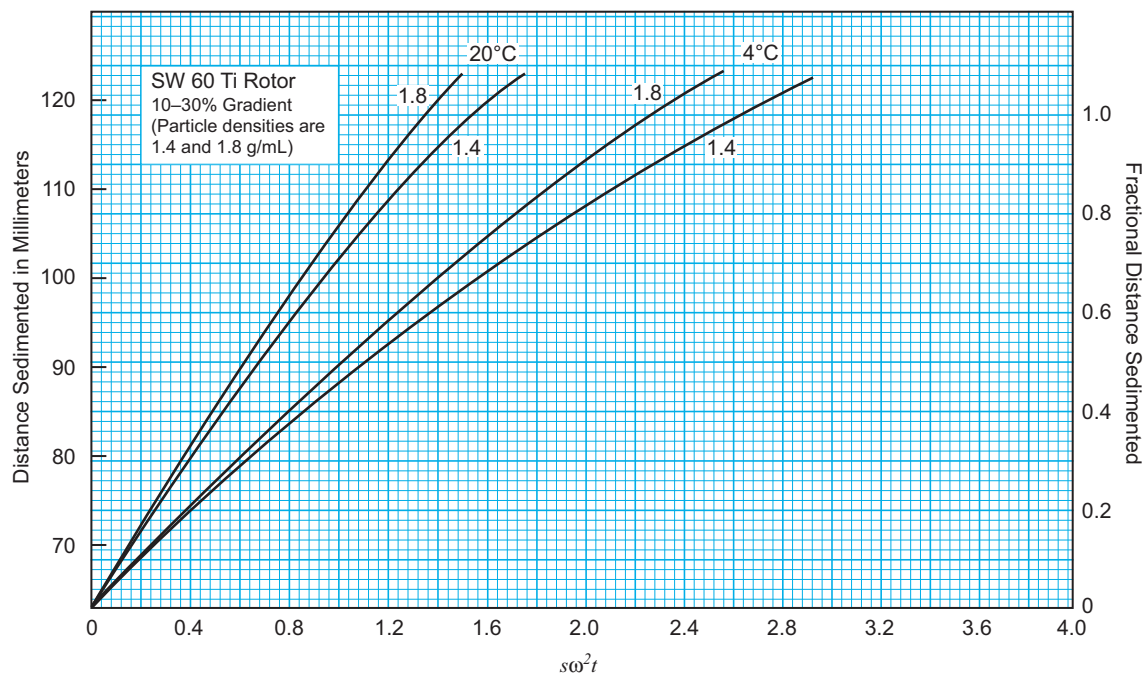

Calculating Sedimentation Coefficients

To calculate sedimentation coefficients, the following must be known: particle density, the distance from the axis of rotation it is to travel, specific properties of the gradient, run speed, and centrifugation time. The value of $\omega^2 t$ is used in place of run speed and time. For example, if a protein of density 1.4 g/mL travels 37 mm down the length of the tube in the SW 60 Ti rotor (37 + 63 mm, or 100 mm from the axis of rotation) through a 10 to 30% gradient at 20°C, the value of $s\omega^2 t$ is 0.92 (from the figure). By dividing the value of $s\omega^2 t$ by the product $\omega^2 t$ (from the integrator), the result is the sedimentation coefficient, in seconds, of the particle:

$$s = \frac{0.92}{\omega^2 t} \quad (\text{B-5})$$

The Use of Cesium Chloride Curves

Cesium Chloride Curves

This Appendix describes how to determine a maximum rotor speed and the final band positions of particles when performing isopycnic separations using cesium chloride gradients. The examples shown here are for the Type 90 Ti rotor only. Similar data and examples for other rotors appear in the applicable rotor manual. Be sure to check the manual for your rotor when calculating run speeds and banding positions.

Rotor speed controls the slope (dr/dr) of a CsCl equilibrium gradient. When planning a separation, gradients should be selected so that the density range from the top to the bottom of the gradient is sufficient to encompass the buoyant densities of particles to be separated. However, speeds must often be limited to avoid precipitation of CsCl at the bottom of the gradient. The density of crystallized CsCl (4 g/mL) produces stresses far in excess of the design limits of most rotors. Also, precipitation will alter the density distribution of the gradient, and the position of sample bands.


The square-root reduction formula—used to determine maximum rotor speeds when centrifuging dense solutions in plastic tubes—does not always guard against CsCl precipitation. The square-root reduction becomes the limiting factor only at relatively high densities and speeds.

Speed and density combinations that intersect on or below the solid curves in [Figure B.1](#) ensure that CsCl will not precipitate in the Type 90 Ti rotor. Curves are provided at two temperatures: 20°C (black lines) and 4°C (gray lines). Note from [Figure B.1](#) that for a given CsCl density, faster rotor speeds can be used as the fill volume in the tube decreases from full to one-quarter filled. Also, for a given rotor speed, the maximum CsCl density that can be safely centrifuged at that speed and temperature increases as the fill volume decreases.

The curves in [Figure B.2](#) show gradient profiles at equilibrium. Each curve was generated for the specific rotor speed shown using the maximum CsCl density (from [Figure B.1](#)) that avoids precipitation at that speed and temperature.* The three-quarter-, one-half-, and one-quarter-filled lines show gradients produced in partially filled tubes. [Figure B.2](#) can be used to approximate banding positions of sample particles. In general, lower speeds generate gradients with shallow slopes; bands will be farther apart. Higher speeds generate gradients with steep slopes where bands will be closer together. Gradient curves not shown can be interpolated.

* Gradients in [Figure B.2](#) result from homogeneous CsCl solutions, but can be more rapidly generated from step or linear gradients, as long as the total amount of CsCl in solution is equal to the amount in the homogeneous solution from the curves in [Figure B.1](#).

Figure B.1 Precipitation Curves for the Type 90 Ti Rotor*

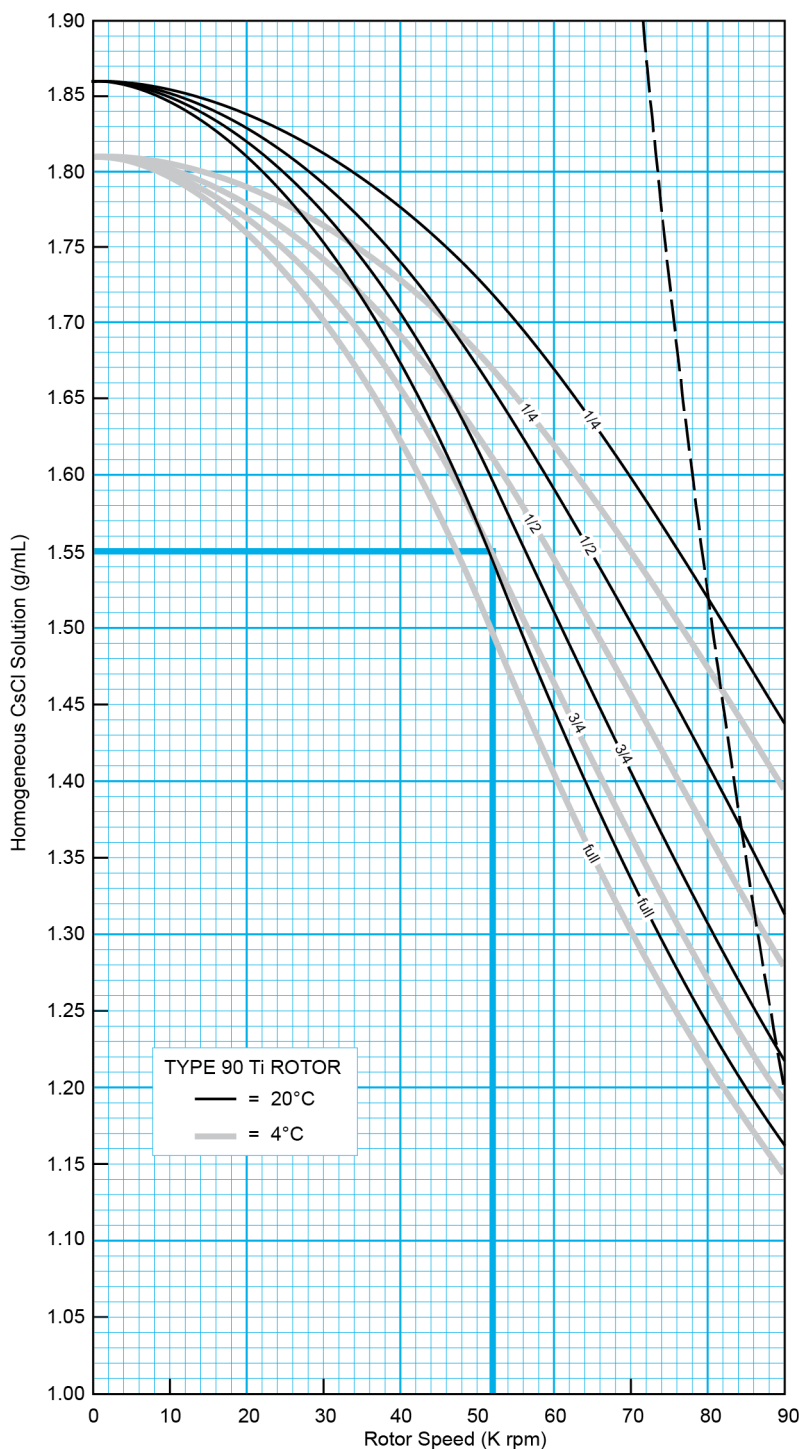

The dashed line is a representation of the equation:


$$\text{RPM} = 90,000 \sqrt{\frac{1.2 \text{ g/mL}}{r}}$$

where r = density of tube contents

and is shown here to illustrate the inability of that equation to guard against CsCl precipitation.

* Using combinations of rotor speeds and homogeneous CsCl solution densities that intersect on or below these curves ensures that CsCl will not precipitate during centrifugation.

Figure B.2 CsCl Gradients at Equilibrium*

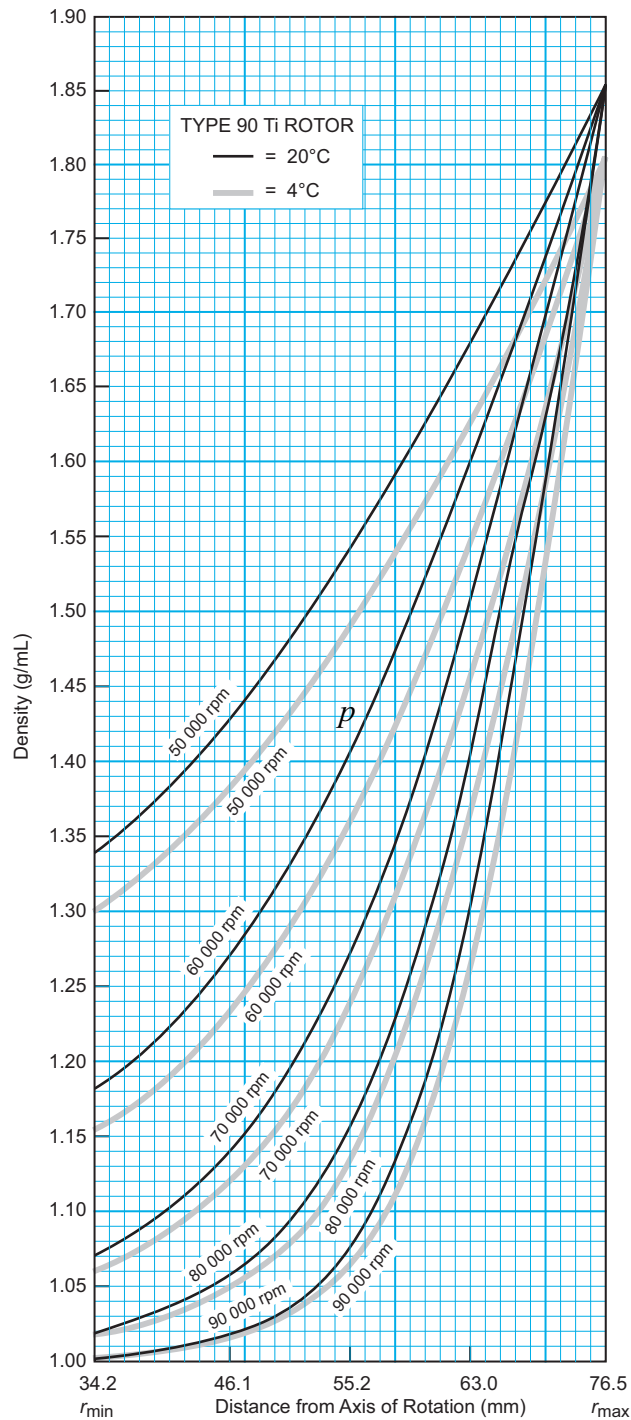

* Centrifugation of homogeneous CsCl solutions at the maximum allowable speeds (from [Figure B.1](#)) results in gradients presented here. Density increases from the top (34.2 mm) to the bottom (76.5 mm) of the tube.

NOTE The curves in [Figure B.1](#) and [Figure B.2](#) are for solutions of CsCl salt only. If other salts are present in significant concentrations, the overall CsCl concentration or the rotor speed must be reduced.

For example, a quarter-filled tube of a 1.52 g/mL homogeneous CsCl solution at 20°C may be centrifuged at 80,000 rpm (see [Figure B.1](#)). The segment of the 80,000 rpm curve ([Figure B.2](#)) from the quarter-filled line to 1.86 g/mL at the tube bottom represents this gradient. The same solution in a half-filled tube ([Figure B.1](#)) may be centrifuged no faster than 68,000 rpm.

Using [Figure B.2](#), interpolate between the 60,000 rpm and 70,000 rpm curves and draw the new 68,000 rpm gradient curve to the half-filled level. The same solution in a three-quarter-filled tube may be centrifuged at 59,000 rpm; [Figure B.2](#) shows the gradient profile (use the three-quarter-filled segment only). A tube full of the 1.52 g/mL CsCl solution may be centrifuged no faster than 53,000 rpm (interpolate and draw in the new gradient profile).

Typical Examples for Determining CsCl Run Parameters

Example A: A separation that is done frequently is the banding of plasmid DNA in cesium chloride with ethidium bromide. The starting density of the CsCl solution is 1.55 g/mL. In this separation the covalently closed, circular plasmid bands at a density of 1.57 g/mL, while the nicked and linear species band at 1.53 g/mL. At 20°C, where will particles band?

- 1 In [Figure B.1](#), find the curve that corresponds to the desired run temperature (20°) and tube fill volume (full).
 - The maximum allowable rotor speed is determined from the point where this curve intersects the homogeneous CsCl density (52,000 rpm).
- 2 In [Figure B.2](#), sketch a horizontal line corresponding to each particle's buoyant density.
- 3 Mark the point where each density intersects the curve corresponding to the maximum speed and selected temperature.
- 4 Particles will band at these points along the tube axis.

In this example, particles will band at about 55.2 and 58.1 mm from the axis of rotation (about 2.9 mm of interband [center-of-band to center-of-band] separation at the 25-degree tube angle). When the tube is held upright, there will be about 3.2 mm of interband separation.

NOTE In swinging bucket rotors, the interband separation after centrifugation is the same as during centrifugation, as there is no gradient reorientation. In fixed angle, near vertical tube, and vertical tube rotors, the gradient must reorient to a horizontal position after centrifugation. Therefore, to determine the interband separation after centrifugation when the tube is held upright (d_{up}) use:

$$d_{up} = \frac{d_{\theta}}{\cos \theta}$$

where d_{θ} is the interband separation achieved during centrifugation, and θ is the tube angle.

Example B: Knowing particle densities (1.50 and 1.52 g/mL), how do you achieve good separation?

- 1 In [Figure B.2](#), sketch in a horizontal line corresponding to each particle's buoyant density.
- 2 Select the curve at the desired temperature (4°C) and tube volume (full) that gives good separation.
- 3 Note the speed indicated along the curve (50,000 rpm).
- 4 From [Figure B.1](#), determine the maximum allowable homogeneous CsCl density that corresponds to the selected temperature, speed, and fill volume from [Figure B.2](#) (in this case 1.51 g/mL).

In this example, particles will band at about 56 and 58 mm from the axis of rotation (about 2 mm of interband separation at the tube angle). When the tube is held upright, there will be about 2.21 mm of interband separation.

To determine the interband volume in millimeters, use:

$$V = \pi r^2 h$$

where r is the tube radius in centimeters and h is the interband separation in centimeters.

The Use of Cesium Chloride Curves

Typical Examples for Determining CsCl Run Parameters

Description

This Appendix contains reference information on commonly used gradient materials. General instructions for filling and sealing tubes, including gradient preparation, are contained in [CHAPTER 3](#).

Gradient material selection depends on a number of factors, including the type of separation to be performed. Sucrose is used for rate zonal and isopycnic separations, and cesium chloride is often used for isopycnic separations. The basic requirement is that the gradient permit the type of separation. Additional considerations in selecting a gradient material include the following.

- Its density range should be sufficient to permit separation of the particles of interest by the chosen density gradient technique, without overstressing the rotor.
- It should not affect the biological activity of the sample.
- It should be neither hyperosmotic or hypoosmotic when the sample is composed of sensitive organelles.
- It should not interfere with the assay technique.
- It should be removable from the purified product.
- It should not absorb in the ultraviolet or visible range.
- It should be inexpensive and readily available; more expensive materials should be recoverable for reuse.
- It should be sterilizable.
- It should not be corrosive to the rotor.
- It should not be flammable or toxic to the extent that its aerosols could be hazardous.

The following charts are provided as a reference for information on commonly used gradient materials.

Table C.1 Commonly Used Gradient Materials with Their Solvents

Materials	Solvent	Maximum Density at 20°C
Sucrose (66%)	H ₂ O	1.32
Sucrose (65%)	D ₂ O	1.37
Silica sols	H ₂ O	1.30
Diodon	H ₂ O	1.37
Glycerol	H ₂ O	1.26
Cesium chloride	H ₂ O	1.91
	D ₂ O	1.98
Cesium formate	H ₂ O	2.10
Cesium acetate	H ₂ O	2.00
Rubidium chloride	H ₂ O	1.49
Rubidium formate	H ₂ O	1.85
Rubidium bromide	H ₂ O	1.63
Potassium acetate	H ₂ O	1.41
Potassium formate	H ₂ O	1.57
	D ₂ O	1.63
Sodium formate	H ₂ O	1.32
	D ₂ O	1.40
Lithium bromide	H ₂ O	1.83
Lithium chloride	D ₂ O	1.33
Albumin	H ₂ O	1.35
Sorbitol	H ₂ O	1.39
Ficoll	H ₂ O	1.17
Metrizamide	H ₂ O	1.46

Table C.2 Density, Refractive Index, and Concentration Data—Cesium Chloride at 25°C,
Molecular Weight = 168.37^a

Density (g/cm ³) ^b	Refractive Index, η_D	% by Weight	mg/mL of Solution ^c	Molarity	Density (g/cm ³) ^b	Refractive Index, η_D	% by Weight	mg/mL of Solution ^c	Molarity
1.0047	1.3333	1	10.0	0.056	1.336	1.3657	34	454.2	2.698
1.0125	1.3340	2	20.2	0.119	1.3496	1.3670	35	472.4	2.806
1.0204	1.3348	3	30.6	0.182	1.363	1.3683	36	490.7	2.914
1.0284	1.3356	4	41.1	0.244	1.377	1.3696	37	509.5	3.026
1.0365	1.3364	5	51.8	0.308	1.391	1.3709	38	528.6	3.140
1.0447	1.3372	6	62.8	0.373	1.406	1.3722	39	548.3	3.257
1.0531	1.3380	7	73.7	0.438	1.4196	1.3735	40	567.8	3.372
1.0615	1.3388	8	84.9	0.504	1.435	1.3750	41	588.4	3.495
1.0700	1.3397	9	96.3	0.572	1.450	1.3764	42	609.0	3.617
1.0788	1.3405	10	107.9	0.641	1.465	1.3778	43	630.0	3.742
1.0877	1.3414	11	119.6	0.710	1.481	1.3792	44	651.6	3.870
1.0967	1.3423	12	131.6	0.782	1.4969	1.3807	45	673.6	4.001
1.1059	1.3432	13	143.8	0.854	1.513	1.3822	46	696.0	4.134
1.1151	1.3441	14	156.1	0.927	1.529	1.3837	47	718.6	4.268
1.1245	1.3450	15	168.7	1.002	1.546	1.3852	48	742.1	4.408
1.1340	1.3459	16	181.4	1.077	1.564	1.3868	49	766.4	4.552
1.1437	1.3468	17	194.4	1.155	1.5825	1.3885	50	791.3	4.700
1.1536	1.3478	18	207.6	1.233	1.601	1.3903	51	816.5	4.849
1.1637	1.3488	19	221.1	1.313	1.619	1.3920	52	841.9	5.000
1.1739	1.3498	20	234.8	1.395	1.638	1.3937	53	868.1	5.156
1.1843	1.3508	21	248.7	1.477	1.658	1.3955	54	859.3	5.317
1.1948	1.3518	22	262.9	1.561	1.6778	1.3973	55	922.8	5.481
1.2055	1.3529	23	277.3	1.647	1.699	1.3992	56	951.4	5.651
1.2164	1.3539	24	291.9	1.734	1.720	1.4012	57	980.4	5.823
1.2275	1.3550	25	306.9	1.823	1.741	1.4032	58	1009.8	5.998
1.2387	1.3561	26	322.1	1.913	1.763	1.4052	59	1040.2	6.178
1.2502	1.3572	27	337.6	2.005	1.7846	1.4072	60	1070.8	6.360
1.2619	1.3584	28	353.3	2.098	1.808	1.4093	61	1102.9	6.550
1.2738	1.3596	29	369.4	2.194	1.831	1.4115	62	1135.8	6.746
1.2858	1.3607	30	385.7	2.291	1.856	1.4137	63	1167.3	6.945
1.298	1.3619	31	402.4	2.390	1.880	1.4160	64	1203.2	7.146
1.311	1.3631	32	419.5	2.492	1.9052	1.4183	65	1238.4	7.355
1.324	1.3644	33	436.9	2.595					

a. Density data are from International Critical Tables.

b. Computed from the relationship $p^{25} = 10.2402 \eta_D^{25} - 12.6483$ for densities between 1.00 and 1.37, and $p^{25} = 10.8601 \eta_D^{25} - 13.4974$ for densities above 1.37 (Bruner and Vinograd, 1965).

c. Divide by 10.0 to obtain % w/v.

Table C.3 Density, Refractive Index, and Concentration Data—Sucrose at 20°C, Molecular Weight = 342.3^a

Density (g/cm ³)	Refractive Index, η_D	% by Weight	mg/mL of Solution ^b	Molarity	Density (g/cm ³)	Refractive Index, η_D	% by Weight	mg/mL of Solution ^b	Molarity
0.9982	1.3330	0			1.1463	1.3883	34	389.7	1.138
1.0021	1.3344	1	10.0	0.029	1.1513	1.3902	35	403.0	1.177
1.0060	1.3359	2	20.1	0.059	1.1562	1.3920	36	416.2	1.216
1.0099	1.3374	3	30.3	0.089	1.1612	1.3939	37	429.6	1.255
1.0139	1.3388	4	40.6	0.119	1.1663	1.3958	38	443.2	1.295
1.0179	1.3403	5	50.9	0.149	1.1713	1.3978	39	456.8	1.334
1.0219	1.3418	6	61.3	0.179	1.1764	1.3997	40	470.6	1.375
1.0259	1.3433	7	71.8	0.210	1.1816	1.4016	41	484.5	1.415
1.0299	1.3448	8	82.4	0.211	1.1868	1.4036	42	498.5	1.456
1.0340	1.3464	9	93.1	0.272	1.1920	1.4056	43	512.6	1.498
1.0381	1.3479	10	103.8	0.303	1.1972	1.4076	44	526.8	1.539
1.0423	1.3494	11	114.7	0.335	1.2025	1.4096	45	541.1	1.581
1.0465	1.3510	12	125.6	0.367	1.2079	1.4117	46	555.6	1.623
1.0507	1.3526	13	136.6	0.399	1.2132	1.4137	47	570.2	1.666
1.0549	1.3541	14	147.7	0.431	1.2186	1.4158	48	584.9	1.709
1.0592	1.3557	15	158.9	0.464	1.2241	1.4179	49	599.8	1.752
1.0635	1.3573	16	170.2	0.497	1.2296	1.4200	50	614.8	1.796
1.0678	1.3590	17	181.5	0.530	1.2351	1.4221	51	629.9	1.840
1.0721	1.3606	18	193.0	0.564	1.2406	1.4242	52	645.1	1.885
1.0765	1.3622	19	204.5	0.597	1.2462	1.4264	53	660.5	1.930
1.0810	1.3639	20	216.2	0.632	1.2519	1.4285	54	676.0	1.975
1.0854	1.3655	21	227.9	0.666	1.2575	1.4307	55	691.6	2.020
1.0899	1.3672	22	239.8	0.701	1.2632	1.4329	56	707.4	2.067
1.0944	1.3689	23	251.7	0.735	1.2690	1.4351	57	723.3	2.113
1.0990	1.3706	24	263.8	0.771	1.2748	1.4373	58	739.4	2.160
1.1036	1.3723	25	275.9	0.806	1.2806	1.4396	59	755.6	2.207
1.1082	1.3740	26	288.1	0.842	1.2865	1.4418	60	771.9	2.255
1.1128	1.3758	27	300.5	0.878	1.2924	1.4441	62	788.3	2.303
1.1175	1.3775	28	312.9	0.914	1.2983	1.4464	62	804.9	2.351
1.1222	1.3793	29	325.4	0.951	1.3043	1.4486	63	821.7	2.401
1.1270	1.3811	30	338.1	0.988	1.3103	1.4509	64	838.6	2.450
1.1318	1.3829	31	350.9	1.025	1.3163	1.4532	65	855.6	2.500
1.1366	1.3847	32	363.7	1.063	1.3224	1.4558	66	872.8	2.550
1.1415	1.3865	33	376.7	1.100	1.3286	1.4581	67	890.2	2.864

a. Density data are from International Critical Tables.

b. Divide by 10.0 to obtain % w/v.

Table C.4 Density Conversion for Cesium and Rubidium Salts at 20°C

% w/w	CsCl	CsBr	CsI	Cs ₂ SO ₄	CsNO ₃	RbCl	RbBr	RbI	Rb ₂ SO ₄	RbNO ₃
1	1.00593	1.00612	1.00608	1.0061	1.00566	1.00561	1.00593	1.00591	1.0066	1.0053
2	1.01374	1.01412	1.01402	1.0144	1.01319	1.01307	1.01372	1.01370	1.0150	1.0125
4	1.02969	1.03048	1.03029	1.0316	1.02859	1.02825	1.02965	1.02963	1.0322	1.0272
6	1.04609	1.04734	1.04707	1.0494	1.04443	1.04379	1.04604	1.04604	1.0499	1.0422
8	1.06297	1.06472	1.06438	1.0676	1.06072	1.05917	1.06291	1.06296	1.0680	1.0575
10	1.08036	1.08265	1.08225	1.0870	1.07745	1.07604	1.08028	1.08041	1.0864	1.0731
12	1.09828	1.10116	1.10071	1.1071	1.09463	1.09281	1.09817	1.09842	1.1052	1.0892
14	1.11676	1.12029	1.11979	1.1275	1.11227	1.11004	1.11661	1.11701	1.1246	1.1057
16	1.13582	1.14007	1.13953	1.1484		1.12775	1.13563	1.13621	1.1446	1.1227
18	1.15549	1.16053	1.15996	1.1696		1.14596	1.15526	1.15605	1.1652	1.1401
20	1.17580	1.18107	1.18112	1.1913		1.16469	1.17554	1.17657	1.1864	1.1580
22	1.19679	1.20362	1.20305	1.2137		1.18396	1.19650	1.19781	1.2083	1.1763
24	1.21849	1.22634	1.22580	1.2375		1.20379	1.21817	1.21980	1.2309	1.1952
26	1.24093	1.24990	1.24942	1.2643		1.22421	1.24059	1.24257	1.2542	1.2146
28	1.26414	1.27435	1.27395			1.24524	1.26380	1.26616	1.2782	1.2346
30	1.28817	1.29973	1.29944			1.26691	1.28784	1.29061	1.3028	1.2552
35	1.35218	1.36764	1.36776			1.32407	1.35191	1.35598	1.3281	1.2764
40	1.42245	1.44275	1.44354			1.38599	1.42233	1.42806		
45	1.49993	1.52626	1.52803			1.45330	1.50010	1.50792		
50	1.58575	1.61970	1.62278			1.52675	1.58639	1.59691		
55	1.68137	1.72492					1.68254	1.69667		
60	1.78859							1.80924		
65	1.90966							1.93722		

List of References

Documents referenced below can be obtained at www.beckman.com or by calling Beckman Coulter at 1-800-742-2345 in the United States, or by contacting your local Beckman Coulter office.*

IN-181	<i>How to Use Quick Seal Tubes with the Beckman Coulter Cordless Tube Topper</i>
IN-189	<i>Using OptiSeal Tubes</i>
IN-192	<i>Use and Care of Centrifuge Tubes and Bottles</i>
IN-197	<i>Rotor Safety (Multi-lingual)</i>
L-ML-5	<i>Master Logbook for Ultracentrifuge Rotors</i>
L-TB-010	<i>Instructions for Using the Beckman Tube Slicer</i>
L5-TB-060	<i>Instructions for Using Aluminum Tube Caps in Fixed Angle Ultracentrifuge Rotors</i>
L5-TB-072	<i>Run Speeds for Stainless Steel Tubes</i>
L5-TB-081	<i>Beckman Fraction Recovery Systems</i>
TLR-IM-9	<i>Rotor and Tubes for Beckman Coulter Tabletop Preparative Ultracentrifuges</i>
BR-8101	<i>Ultracentrifuge Rotors, Tubes & Accessories Catalog</i>
IN-175	<i>Chemical Resistance Table for Beckman Coulter Products</i>

Refer to www.beckman.com for additional information on Beckman Coulter Life Sciences products, resources, service, and support.

* For detailed information on a rotor, see the applicable individual rotor manual.

References

List of References

Glossary

Angular velocity, ω — Rate of rotation, measured in radians per second

$$\omega = \frac{2\pi \text{ rpm}}{60}$$

or

$$\omega = 0.10472 \text{ rpm}$$

Anodized coating — A thin, hard layer of aluminum oxide formed electrochemically on aluminum rotor and/or accessory surfaces as a protective coating for corrosion resistance

Autoclaving — Sterilization by heat (dry or steam)

Buoyant density — The density of a particle in a specified liquid medium

Buna N — Black nitrile rubber used for O-rings and gaskets in rotor assemblies; should be used at temperatures between -34 and 121°C (-30 and 250°F)

Centrifugal effect — Accumulated value of:

$$\int_{t_1}^{t_2} \omega^2 dt$$

where t is time and ω is angular velocity

Centrifugal force — In a centrifugal field, the force which causes a particle to move away from the center of rotation

Clearing factor k — Calculated for all Beckman Coulter high-speed rotors as a measure of the rotor's relative pelleting efficiency:

$$k = \frac{\ln(r_{\max}/r_{\min})}{\omega^2} \times \frac{10^{13}}{3600}$$

or

$$k = \frac{253303 \times \ln(r_{\max}/r_{\min})}{(\text{RPM} / 1000)^2}$$

Clearing time, t — $t = k/s$, where t is time in hours, k is the clearing factor of the rotor, and s is the sedimentation coefficient in Svedberg units (S)

CsCl — Cesium chloride; a high-density salt used in solution in isopycnic separations to separate particles based on their density

Cs₂SO₄ — Cesium sulfate; a salt, similar to CsCl, that will form its own gradient in solution

Delrin — Thermoplastic material (acetal homopolymer) used for most tube adapters (Delrin is a registered trademark of E.I. Du Pont de Nemours & Company.)

Density — Mass per unit volume

Density separation — A centrifugal separation process based on differences in particle densities

Differential separation — A centrifugal separation process based on differences in particle sizes

EPDM — Ethylene propylene rubber used for O-rings and pad adapters; should be used at temperatures between –57 and 120°C (–70 and 250°F)

Ethidium bromide — A fluorescent intercalating orange dye used commonly in the separation of DNA and in gel electrophoresis

Fixed-angle rotor — A rotor in which the tubes are held at an angle (usually 20 to 45 degrees) from the axis of rotation

g-Max — A system of centrifugation using a combination of short Quick-Seal tubes and floating spacers, designed to reduce volumes while maximizing separation efficiency

HDPE — High density polyethylene used for adapters

Isopycnic — A method of particle separation or isolation based on particle buoyant density; particles are centrifuged until they reach a point in the gradient where the density of the particle is the same as the density of the gradient at that point

konical tubes — Thin-walled, polypropylene tubes featuring a conical tip to optimize pelleting separations; the conical tip concentrates the pellet in the narrow base of the tube. Available in both open-top and Quick-Seal bell-top designs.

Maximum volume — The maximum volume at which a tube should be filled for centrifugation (sometimes referred to as maximum fill volume or nominal fill volume)

Mechanical overspeed cartridge — An assembly installed in the bases of some older rotors or swinging-bucket rotor adapters as part of the mechanical overspeed protection system

Meniscus — The curved upper surface of a liquid column that is concave when the container walls are wetted by the liquid and convex when they are not

Near-vertical tube rotor — A rotor in which the tubes are held at a slight angle (usually 7 to 10 degrees)

Neoprene — Black synthetic elastomer used for O-rings in some tube caps and bottle cap assemblies; should be used at temperatures between –54 and 121°C (–65 and 250°F)

Noryl — Modified thermoplastic polyphenylene oxide (PPO) used for floating spacers (part of the g-Max system) and some polycarbonate bottle caps (Noryl is a registered trademark of SHPP GLOBAL TECHNOLOGIES B.V.)

OptiSeal tubes — Capless tubes with sealing plugs inserted in the tube stems; during centrifugation, the combination of g force and hydrostatic pressure seals the tube

Overspeed disk — An adhesive disk, with alternating reflecting and nonreflecting sectors, attached to the bottom of rotors as part of the photoelectric overspeed protection system; the number of sectors on the disk is a function of the rotor's maximum allowable speed

Pelleting — A centrifugal separation process in which particles in a sample sediment to the bottom of the tube (differential separation); differential pelleting separates particles of different sizes by successive centrifugation steps of progressively higher g force and/or longer run duration

PET — polyethylene terephthalate used in some adapters

Polypropylene — Random block copolymer of ethylene and propylene used for certain tubes (Tenite Polypropylene is a registered trademark of Eastman Chemical Co.)

Quick-Seal tubes — bell-top or dome-top thinwall tubes that are heat-sealed and require no caps

Radel — Polyphenylsulfone (PPS) used in plugs, cap closures, cannisters and other accessories

Rate zonal — A method of particle separation, based on differential rate of sedimentation, using a preformed gradient with the sample layered as a zone on top of the gradient

RCF — Relative centrifugal field; the ratio of the centrifugal acceleration at a specified radius and speed ($r\omega^2$) to the standard acceleration of gravity (g) according to the following equation:

$$\text{RCF} = \frac{r\omega^2}{g}$$

where r is the radius in millimeters, ω is the angular velocity in radians per second ($2\pi \text{ RPM}/60$), and g is the standard acceleration of gravity (9807 mm/s^2). Thus the relationship between RCF and RPM is:

$$\text{RCF} = 1.12r \left(\frac{\text{RPM}}{1000} \right)^2$$

r_{max} — (Maximum radius) the position of the liquid in the tube at the maximum distance from the axis of rotation when the rotor is at speed

r_{min} — (Minimum radius) the position of the liquid in the tube at the minimum distance from the axis of rotation when the rotor is at speed

Sedimentation — The settling out of particles from a suspension in the earth's field of gravity; in the centrifuge this process is accelerated and the particles move away from the axis of rotation

Sedimentation coefficient, s — Sedimentation velocity per unit of centrifugal force:

$$s = \frac{dr}{dt} \times \frac{1}{\omega^2 r}$$

Silicone rubber — A large group of silicone elastomers used in various accessories; should be used at temperatures between -59 and 232°C (-75 and 450°F)

Solution 555 — Beckman Coulter concentrated rotor cleaning solution; recommended because it is a mild solution that has been tested and found effective and safe for Beckman Coulter rotors and accessories

Spinkote — Beckman Coulter lubricant for metal-to-metal contacts

Sucrose — A sugar (not a self-forming gradient) used in rate zonal separations; generally used in separating RNA, subcellular organelles, and cell membranes

Supernatant — The liquid above the sedimented material following centrifugation

Svedberg unit, S — A unit of sedimentation velocity:

$$1\text{ S} = 10^{-13}\text{ seconds}$$

Swinging-bucket rotor — A rotor in which the tubes or bottles are carried in buckets, microtiter plate carriers, or racks that swing up to the horizontal position during centrifugation (sometimes referred to as a horizontal or swing-out rotor)

Ultem — Polyetherimide (PEI)—used in adapters, covers, and spacers; should be used at temperatures between –29 and 204°C (–20 and 400°C) (Ultem is a registered trademark of SABIC GLOBAL TECHNOLOGIES B.V.)

Vertical-tube rotor — A rotor in which the tubes or bottles are held parallel to the axis of rotation

Viton — Fluorocarbon elastomer used in high-temperature applications (Viton is a registered trademark of E.I. Du Pont de Nemours & Company.)

Wettable — Tube or bottle material that water or other aqueous solution will adhere to; the more wettable a tube or bottle material is, the more biological material, DNA, protein, cells, and so forth, will adhere to the walls

Ultracentrifuge Rotor Warranty

All Beckman Coulter ultracentrifuge Fixed-Angle, Vertical-Tube, Near-Vertical Tube, Swinging-Bucket, and Airfuge rotors are warranted against defects in materials or workmanship for the time periods indicated below, subject to the Warranty Conditions stated below.

Preparative Ultracentrifuge Rotors	5 years, No Proration
Analytical Ultracentrifuge Rotors	5 years, No Proration
ML and TL Series Ultracentrifuge Rotors	5 years, No Proration
Airfuge Ultracentrifuge Rotors	1 year, No Proration

For Zonal, Continuous Flow, Component Test, and Rock Core ultracentrifuge rotors, see separate warranty.

Warranty Conditions (as applicable)

1. This warranty is valid for the time periods indicated above from the date of shipment to the original Buyer by Beckman Coulter or an authorized Beckman Coulter representative.
2. This warranty extends only to the original Buyer and may not be assigned or extended to a third person without written consent of Beckman Coulter.
3. This warranty covers the Beckman Coulter Centrifuge Systems only (including but not limited to the centrifuge, rotor, and accessories) and Beckman Coulter shall not be liable for damage to or loss of the user's sample, non-Beckman Coulter tubes, adapters, or other rotor contents.
4. This warranty is void if the Beckman Coulter Centrifuge System is determined by Beckman Coulter to have been operated or maintained in a manner contrary to the instructions in the operator's manual(s) for the Beckman Coulter Centrifuge System components in use. This includes, but is not limited to operator misuse, abuse, or negligence regarding indicated maintenance procedures; centrifuge and rotor classification requirements; proper speed reduction for the high density of certain fluids, tubes, and tube caps; speed reduction for precipitating gradient materials; and speed reduction for high-temperature operation.
5. Rotor bucket sets purchased concurrently with or subsequent to the purchase of a Swinging-Bucket Rotor are warranted only for a term co-extensive with that of the rotor for which the bucket sets are purchased.
6. This warranty does not cover the failure of a Beckman Coulter rotor in a centrifuge not of Beckman Coulter manufacture, or if the rotor is used in a Beckman Coulter centrifuge that has been modified without the written permission of Beckman Coulter, or, if the rotor is used with carriers, buckets, belts, or other devices not of Beckman Coulter manufacture.
7. Rotor parts subject to wear, including, but not limited to, rotor O-rings, VTi, NVT, TLV, MLN, and TLN rotor tube cavity plugs and gaskets, tubing, tools, optical overspeed disks, bearings, seals, and lubrication are excluded from this warranty and should be frequently inspected and replaced if they become worn or damaged.
8. Keeping a rotor log is not mandatory, but may be desirable for maintenance of good laboratory practices.

Repair and Replacement Policies

1. If a Beckman Coulter rotor is determined by Beckman Coulter to be defective, Beckman Coulter will repair or replace it, subject to the Warranty Conditions. A replacement rotor will be warranted for the time remaining on the original rotor's warranty.
2. If a Beckman Coulter centrifuge is damaged due to a failure of a rotor covered by this warranty, Beckman Coulter will supply free of charge (i) all centrifuge parts required for repair (except the drive unit, which will be replaced at the then current price less a credit determined by the total number of revolutions or years completed, provided that such a unit was manufactured or rebuilt by Beckman Coulter), and (ii) if

the centrifuge is currently covered by a Beckman Coulter warranty or Full Service Agreement, all labor necessary for repair of the centrifuge.

3. If a Beckman Coulter rotor covered by this warranty is damaged due to a malfunction of a Beckman Coulter ultracentrifuge covered by an Ultracentrifuge System Service Agreement, Beckman Coulter will repair or replace the rotor free of charge.
4. If a Beckman Coulter rotor covered by this warranty is damaged due to a failure of a Beckman Coulter tube, bottle, tube cap, spacer, or adapter, covered under the Conditions of this Warranty, Beckman Coulter will repair or replace the rotor and repair the instrument as per the conditions in policy point (2) above, and the replacement policy.
5. Damage to a Beckman Coulter rotor or instrument due to the failure or malfunction of a non-Beckman Coulter tube, bottle, tube cap, spacer, or adapter is not covered under this warranty, although Beckman Coulter will assist in seeking compensation under the manufacturer's warranty.

Disclaimer

IT IS EXPRESSLY AGREED THAT THE ABOVE WARRANTY SHALL BE IN LIEU OF ALL WARRANTIES OF FITNESS AND OF THE WARRANTY OF MERCHANTABILITY AND BECKMAN COULTER, INC. SHALL HAVE NO LIABILITY FOR SPECIAL OR CONSEQUENTIAL DAMAGES OF ANY KIND WHATSOEVER ARISING OUT OF THE MANUFACTURE, USE, SALE, HANDLING, REPAIR, MAINTENANCE, OR REPLACEMENT OF THE PRODUCT.

Factory Rotor Inspection Service

Beckman Coulter, Inc., will provide free mechanical and metallurgical inspection in Indianapolis, Indiana, USA, of any Beckman Coulter rotor at the request of the user. (Shipping charges to Beckman Coulter are the responsibility of the user.) Rotors will be inspected in the user's laboratory if the centrifuge in which they are used is covered by an appropriate Beckman Coulter Service Agreement. Contact your local Beckman Coulter office for details of service coverage or cost.

Before shipping, contact the nearest Beckman Coulter Sales and Service office and request a Returned Goods Authorization (RGA) form and packaging instructions. Please include the complete rotor assembly, with buckets, lid, handle, tube cavity caps, etc. A SIGNED STATEMENT THAT THE ROTOR AND ACCESSORIES ARE NON-RADIOACTIVE, NON-PATHOGENIC, NON-TOXIC, AND OTHERWISE SAFE TO SHIP AND HANDLE IS REQUIRED.

www.beckman.com

