

ACCESORIES

Empty Cartridges (1 ml)

PRODUCT

Empty Cartridges (1 ml)
Contains: 5 Empty Cartridges (5 Cartridge Bodies and 10 End Plugs)⁽¹⁾



Cartridge Body



End Plug

CAT. No.

EB-Ctg1-5

COLUMN MATERIAL

Polypropylene

FRIT MATERIAL

Polyethylene

FRIT PORE SIZE AVERAGE

12 µm⁽²⁾

INNER COLUMN DIAMETER

6.2 mm

CONNECTIONS

Standard connection compatible to the common chromatography instruments (such as ÄKTA)

CHEMICAL STABILITY

Stable in all commonly used reagents

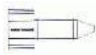
CAPACITY (ml packed resin)

1 ml⁽³⁾

- (1) For storage the packed cartridge, Stop Plugs need be ordered in separate (See accessories Cartridges).
- (2) The customer can use any type of chromatography media (> 20 µm) including ABT range of products for packing their cartridges. We recommend to pack using RAPID RUNTM Fine Resins from ABT (ask for information). After the chosen resin is packed in the cartridge, the cartridge can be readily stored and re-used for many times.
- (3) The packed volume depends on the type of resin. As a guide, the packed volume should be 1.2 – 1.4 ml for agarose based spherical particles. The recommended packed volume is 1.3 ml for Rapid RunTM Fine Resins from ABT (See below).

Note: If longer bed height is needed, 3 cartridges can be inter-connected to give a final bed height of 10 cm. Larger amounts available on request.

ACCESSORY CARTRIDGES:

ACCESSORY		CAT. N°	COMMENTS
Stop Plug (1)	10 / pack	SEB-10	For storage of the packed cartridge, Stop Plugs need be ordered separately.
			
Syringe Connector	1 /pack	SCEB-1	10-32 male / female luer connector. This adapter is used to connect syringes.

Larger amounts available on request.

ABT TS EB-Ctg1 Rev. 2022/A